



# Standard Test Method for Radiochemical Determination of Uranium Isotopes in Urine by Alpha Spectrometry<sup>1</sup>

This standard is issued under the fixed designation C1473; the number immediately following the designation indicates the year of original adoption or, in the case of revision, the year of last revision. A number in parentheses indicates the year of last reapproval. A superscript epsilon ( $\epsilon$ ) indicates an editorial change since the last revision or reapproval.

## 1. Scope

1.1 This test method is applicable to the determination of uranium in urine at levels of detection dependent on sample size, count time, detector background, and tracer yield. It is designed as a screening tool for detection of possible exposure of occupational workers.

1.2 This test method is designed for 50 mL of urine. This test method does not address the sampling protocol or sample preservation methods associated with its use.

1.3 The values stated in SI units are to be regarded as standard. No other units of measurement are included in this standard.

1.4 *This standard does not purport to address all of the safety concerns, if any, associated with its use. It is the responsibility of the user of this standard to establish appropriate safety and health practices and determine the applicability of regulatory limitations prior to use.*

## 2. Referenced Documents

### 2.1 ASTM Standards:<sup>2</sup>

C859 Terminology Relating to Nuclear Materials

C1000 Test Method for Radiochemical Determination of Uranium Isotopes in Soil by Alpha Spectrometry

C1163 Practice for Mounting Actinides for Alpha Spectrometry Using Neodymium Fluoride

C1284 Practice for Electrodeposition of the Actinides for Alpha Spectrometry

D1193 Specification for Reagent Water

D3084 Practice for Alpha-Particle Spectrometry of Water

D3648 Practices for the Measurement of Radioactivity

<sup>1</sup> This test method is under the jurisdiction of ASTM Committee C26 on Nuclear Fuel Cycle and is the direct responsibility of Subcommittee C26.05 on Methods of Test.

Current edition approved Feb. 15, 2011. Published February 2011. Originally approved in 2000. Last previous edition approved in 2005 as C1473 – 05. DOI: 10.1520/C1473-11.

<sup>2</sup> For referenced ASTM standards, visit the ASTM website, [www.astm.org](http://www.astm.org), or contact ASTM Customer Service at [service@astm.org](mailto:service@astm.org). For *Annual Book of ASTM Standards* volume information, refer to the standard's Document Summary page on the ASTM website.

## 3. Summary of Test Method

3.1 A urine sample with <sup>232</sup>U tracer solution added is wet-ashed with nitric acid and hydrogen peroxide to destroy organic material. The uranium-bearing solution is converted to a hydrochloric acid medium. Uranium is absorbed on an anion exchange column from a 9 M hydrochloric acid solution and eluted with 0.1 M hydrochloric acid solution. The separated uranium is prepared for alpha spectrometric measurement either by electrodeposition onto a metal disk or coprecipitation with neodymium fluoride and filtration onto a membrane filter.

## 4. Significance and Use

4.1 This test method is used to detect possible exposures to uranium isotopes from occupational operations.

## 5. Interferences

5.1 The presence of <sup>232</sup>U in the urine sample will be masked by the tracer addition. The alpha energies of <sup>233</sup>U and <sup>234</sup>U cannot be fully resolved by alpha spectrometric measurement. A table of uranium isotope alpha energies is given in [Appendix X1](#). If neptunium is present in the sample in the plus four oxidation state, it will co-elute with the uranium.

## 6. Apparatus

6.1 *Alpha Spectrometry System*—Refer to Test Method C1000, and Practices C1163, D3084, and D3648 for guidance.

6.2 *Electrodeposition Apparatus*—Refer to Practice C1284 for guidance.

6.3 *Neodymium Fluoride Precipitation*—Refer to Practice C1163 for guidance.

6.4 *Borosilicate Beakers or Flasks*, 250 mL.

6.5 *Borosilicate Beakers*, 150 and 250 mL.

6.6 *Borosilicate Graduated Glass Cylinders*, 5, 25, 100, and 1000 mL.

6.7 *Ion Exchange Columns*, disposable polypropylene, with polyethylene frit,  $\geq 5$ -mL capacity.

6.8 *Ion Exchange Column Reservoir*, funnel, polypropylene, 100-mL capacity.

NOTE 1—See [Fig. 1](#) for a typical ion exchange column-reservoir setup.

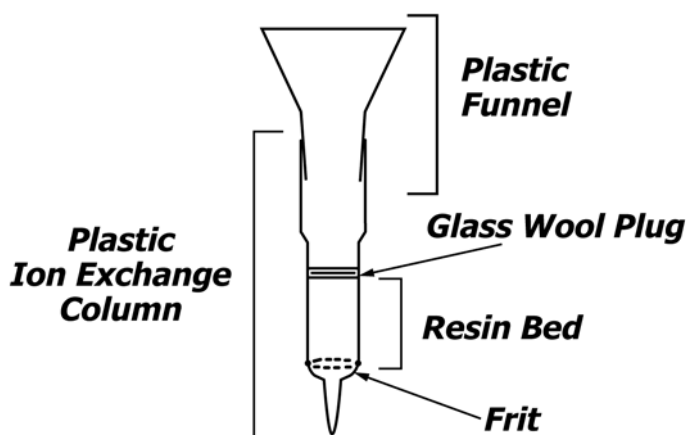


FIG. 1 Typical Ion Exchange Column Arrangement

## 7. Reagents

7.1 *Purity of Reagents*—Reagent grade chemicals shall be used in all tests. Unless otherwise indicated, it is intended that the reagents conform to the specifications of the Committee on Analytical Reagents of the American Chemical Society, where such specifications are available.<sup>3</sup> Other grades may be used, provided it is first ascertained that the reagent is of sufficiently high purity to permit its use without lessening the accuracy of the determination.

7.2 *Purity of Water*—Unless otherwise indicated, references to water shall be understood to mean reagent water as defined by Type III of Specification D1193.

7.3 Reagent purity shall be such that the measured radioactivity of blank samples is minimized.

7.4 *Anion Exchange Resin*—Analytical grade, Type 1, 8 % cross linked, 100-200 mesh, chloride form.

7.5 *Hydrochloric Acid (9 M)*—Mix 750 mL of concentrated HCl with water and dilute to 1 L.

7.6 *Hydrochloric Acid (sp gr 1.19)*—Concentrated hydrochloric acid (HCl).

7.7 *Hydrochloric Acid (0.1 M)*—Mix 8.3 mL of concentrated HCl with water and dilute to 1 L.

7.8 *Hydrogen Peroxide (30 %)*.

7.9 *Nitric Acid (8 M)*—Mix 500 mL of concentrated nitric acid with water and dilute to 1 L.

7.10 *Nitric Acid (sp gr 1.4)*—Concentrated nitric acid (HNO<sub>3</sub>).

7.11 *Uranium-232*, standard solution.<sup>4</sup>

<sup>3</sup> *Reagent Chemicals, American Chemical Society Specifications*, American Chemical Society, Washington, DC. For suggestions on the testing of reagents not listed by the American Chemical Society, see *Analar Standards for Laboratory Chemicals*, BDH Ltd., Poole, Dorset, U.K., and the *United States Pharmacopeia and National Formulary*, U.S. Pharmacopeial Convention, Inc. (USPC), Rockville, MD.

<sup>4</sup> Such tracer radionuclides may be available from respective national standards laboratories as well as commercial vendors that provide radionuclides traceable to their respective national standards laboratories.

## 8. Hazards

8.1 **Warning**—Adequate laboratory facilities, such as fume hoods and controlled ventilation, along with safe techniques must be used in this procedure. Site-specific policies for the handling of biological materials must be adhered to. Extreme care should be exercised in using hydrofluoric and other hot, concentrated acids. The use of proper gloves is recommended.

## 9. Sampling

9.1 Collect the urine sample in accordance with the site-specific protocol.

9.2 Preserve the urine sample in accordance with the site-specific protocol.

## 10. Calibration and Standardization

10.1 A <sup>232</sup>U standard may be used to determine the detection efficiency of the alpha-spectrometry system which, in turn, can be used to calculate the chemical yield of each sample and the lower limit of detection (LLD) of this test method.

## 11. Procedure

### 11.1 Sample Preparation:

11.1.1 Measure 50 mL of urine in a 100-mL graduated cylinder.

11.1.2 Transfer the urine to a 250-mL beaker or flask.

11.1.3 Rinse the cylinder twice with 5 mL of 8 M HNO<sub>3</sub> and add the rinsings to the beaker or flask.

11.1.4 Add an appropriate amount of <sup>232</sup>U tracer solution to the sample (ca. 0.02 Bq or as prescribed in the site-specific protocol) and swirl the vessel to mix.

11.1.5 Evaporate the contents of the vessel to near dryness on a medium hot plate.

11.1.6 Remove the vessel from the hot plate and cool for 1 min.

11.1.7 Add 10 mL of concentrated HNO<sub>3</sub> to the vessel. Swirl to mix.

11.1.8 Return the vessel to the hot plate and evaporate the solution to near dryness.

11.1.9 Remove the vessel from the hot plate and allow to cool for 1 min.

11.1.10 Slowly add 10 to 15 mL of 30 % hydrogen peroxide to the vessel to cover the residue. Swirl the vessel to mix.

11.1.11 Return the vessel to the hot plate and evaporate the solution to near dryness.

11.1.12 Repeat 11.1.9 and 11.1.10.

11.1.13 Return the vessel to the hot plate and evaporate the solution to obtain a white sample residue.

11.1.14 Remove the vessel from the hot plate and allow to cool for 1 min.

11.1.15 Add 10 to 15 mL of concentrated HCl to the vessel and swirl to mix.

11.1.16 Return the vessel to the hot plate and evaporate the solution to near dryness.

11.1.17 Remove the vessel from the hot plate and allow to cool for 1 min.

11.1.18 Repeat 11.1.15 – 11.1.17 twice more. Proceed to 11.1.19.

11.1.19 Add 75 mL of 9 M HCl to the vessel and heat gently to dissolve the residue.

11.1.20 Remove the vessel from the hot plate and cool to room temperature.

#### 11.2 Anion Exchange Separation:

11.2.1 Measure 3 mL of settled, water slurried anion exchange resin.

11.2.2 Place a 250-mL beaker under an empty column. Transfer the resin to the column with water. Allow the water to drain from the column.

11.2.3 Insert a glass wool plug on the top of the resin bed.

11.2.4 Pour 30 mL of 9 M HCl into the ion exchange column reservoir and allow the acid to drain through to condition the resin column.

11.2.5 Transfer the sample solution to the ion exchange column reservoir. Pass the sample solution through the resin column.

11.2.6 Rinse the sample vessel with 10 mL of 9 M HCl and add the rinse to the column reservoir. Repeat the rinse once and allow the solution to drain to the top of the glass wool plug.

11.2.7 Wash the resin column with 30 mL of 9 M HCl and allow the acid to drain from the column.

11.2.8 Remove the beaker containing the combined rinse and wash solutions from beneath the resin column. Discard the combined solutions in accordance with the site-specific disposal requirements for acids.

11.2.9 Place a clean 150-mL beaker under the resin column.

11.2.10 Elute the uranium from the resin column by pouring 50 mL of 0.1 M HCl into the column reservoir and allowing it to drain through.

11.2.11 Place the beaker containing the eluted uranium on a hot plate and evaporate the solution to just dryness. Remove the beaker from the hot plate and cool for 1 min.

11.2.12 Add 2 to 3 mL of concentrated HNO<sub>3</sub> to the beaker.

11.2.13 Repeat 11.2.11 and 11.2.12 until the residue is white.

11.2.14 Remove the beaker from the hot plate and cool for 1 min.

11.2.15 Add 5 mL of concentrated HCl to the beaker and evaporate to dryness. Remove the beaker from the hot plate and cool for 1 min.

11.2.16 Repeat 11.2.15 twice more.

#### 11.3 Electrodeposition or Neodymium Fluoride Precipitation:

11.3.1 For electrodeposition, refer to Practice C1284.

11.3.2 For neodymium fluoride precipitation, refer to Practice C1163.

#### 11.4 Alpha Spectrometry Measurement:

11.4.1 Measure the sample overnight (1000 min) or longer in an alpha spectrometry system.

11.4.2 Measure the background and reagent blank activity for the regions of interest. Correct each uranium isotope activity for background and reagent blank activity.

## 12. Calculation

12.1 The activity of each uranium isotope is calculated as follows:

TABLE 1 Summary of Precision and Bias Data

	Sample Type	
	Blind Quality Control	Reference
Number of analyses	103	140
Concentration <sup>234</sup> U (Bq L <sup>-1</sup> )	0.17	0.025
Precision, %	30	11
Bias, %	4	1.5
Concentration <sup>238</sup> U (Bq L <sup>-1</sup> )	0.042	0.058
Precision, %	20	15
Bias, %	7	6
Reagent Blanks		
Number of analyses	80	
	mean ± SD	
<sup>234</sup> U	0.0007 ± 0.0013	
<sup>238</sup> U	0.0008 ± 0.0013	

$$X_i = (C_i A_i) / (C_i V_s) \quad (1)$$

where:

$X_i$  = concentration of a uranium isotope in Bq/L,

$A_i$  = activity of the yield tracer added in Bq,

$C_i$  = net sample counts in the energy region of the uranium isotope being measured,

$C_t$  = net sample counts in the energy region of the uranium-232 tracer, and

$V_s$  = volume of the sample aliquant in litres.

12.2 The fractional counting efficiency of the alpha spectrometer,  $E$ , must be determined if it is desired to calculate the chemical yield for uranium. Calculate the fractional efficiency as follows:

$$E = R_s A_{std} \quad (2)$$

where:

$R_s$  = net counting rate of the standard source in the region of the alpha emitter(s) of interest in cps, and

$A_{std}$  = activity of the standard source in Bq at the time of the count.

12.3 Calculate the fractional chemical yield for uranium as follows:

$$Y = (C_i) / (T A_i E) \quad (3)$$

where:

$T$  = counting time in s,

$C_i$  = the net sample counts in the energy region of the uranium-232 tracer,

$A_i$  = activity of the yield tracer added in Bq, and

$E$  = fractional counting efficiency of the alpha spectrometer.

12.4 Refer to Practice D3084 for the Minimum Detectable Concentration (MDC), critical level, and uncertainty calculations.

## 13. Precision and Bias

13.1 *Precision*—Replicate reference samples, blind quality control (QC) samples, and reagent blanks were analyzed at the Westinghouse Savannah River Company Radiobioassay Laboratory. These analyses were performed over a 12-month period

by several analysts as blind samples in the routine sample stream. The results of these studies are summarized in **Table 1**.

**13.2 Bias**—The bias is defined as follows:  
 Bias, % = 100\* (observed value-known value)/known value  
 averaged over the entire data set

The uncertainty in the NIST standard of 2 % is included in this estimate.

## 14. Keywords

14.1 alpha spectrometry; ion exchange separation; uranium; urine

## APPENDIX

### (Nonmandatory Information)

#### X1.

**TABLE X1.1 Properties of Uranium Isotopes of Interest in Environmental Samples<sup>A</sup>**

Isotope	Half-Life, years	Alpha Energies, MeV (abundance)
<sup>232</sup> U	68.9	5.320 (68.6 %), 5.263 (31.2 %)
<sup>233</sup> U	1.592 × 10 <sup>5</sup>	4.825 (84.4 %), 4.783 (13.2 %)
<sup>234</sup> U	2.454 × 10 <sup>5</sup>	4.776 (72.5 %), 4.723 (27.5 %)
<sup>235</sup> U	7.037 × 10 <sup>8</sup>	4.395 (85 %), 4.370 (6 %), 4.597, (5 %)
<sup>236</sup> U	2.342 × 10 <sup>7</sup>	4.494 (74 %), 4.445 (26 %)
<sup>238</sup> U	4.468 × 10 <sup>9</sup>	4.196 (77 %), 4.147 (23 %)

<sup>A</sup> Available from Browne, E., and Firestone, R. B., *Table of Radioactive Isotopes*, (V. S. Shirley, Ed.), John Wiley and Sons, Inc., 1986.

*ASTM International takes no position respecting the validity of any patent rights asserted in connection with any item mentioned in this standard. Users of this standard are expressly advised that determination of the validity of any such patent rights, and the risk of infringement of such rights, are entirely their own responsibility.*

*This standard is subject to revision at any time by the responsible technical committee and must be reviewed every five years and if not revised, either reapproved or withdrawn. Your comments are invited either for revision of this standard or for additional standards and should be addressed to ASTM International Headquarters. Your comments will receive careful consideration at a meeting of the responsible technical committee, which you may attend. If you feel that your comments have not received a fair hearing you should make your views known to the ASTM Committee on Standards, at the address shown below.*

*This standard is copyrighted by ASTM International, 100 Barr Harbor Drive, PO Box C700, West Conshohocken, PA 19428-2959, United States. Individual reprints (single or multiple copies) of this standard may be obtained by contacting ASTM at the above address or at 610-832-9585 (phone), 610-832-9555 (fax), or service@astm.org (e-mail); or through the ASTM website (www.astm.org). Permission rights to photocopy the standard may also be secured from the Copyright Clearance Center, 222 Rosewood Drive, Danvers, MA 01923, Tel: (978) 646-2600; http://www.copyright.com/*