



# Standard Test Method for Particle Size Distribution of Alumina or Quartz Powders by Electrical Sensing Zone Technique<sup>1</sup>

This standard is issued under the fixed designation C690; the number immediately following the designation indicates the year of original adoption or, in the case of revision, the year of last revision. A number in parentheses indicates the year of last reapproval. A superscript epsilon ( $\epsilon$ ) indicates an editorial change since the last revision or reapproval.

## 1. Scope\*

1.1 This test method, one of several found valuable for the measurement of particle size, covers the determination of the particle size distribution of alumina or quartz powders (0.6 to 56.0  $\mu\text{m}$ ) using electrical sensing zone particle size analyzers. These instruments use an electric current path of small dimensions which is modulated by individual particle passage through an aperture, and produces individual pulses of amplitude proportional to the particle volume.

1.2 The values stated in SI units are to be regarded as standard. No other units of measurement are included in this standard.

1.3 *This standard does not purport to address all of the safety problems, if any, associated with its use. It is the responsibility of the user of this standard to establish appropriate safety and health practices and determine the applicability of regulatory limitations prior to use.*

## 2. Summary of Test Method

2.1 A carefully dispersed, dilute suspension of the powder in a beaker filled with an electrolyte is placed on the instrument sample stand. The suspension is forced through a restricting aperture. Each particle passing generates an electric pulse that is recorded on an electronic counter.

2.2 The instrument response is essentially related to particle volume (liquid displacement). Equivalent spherical diameter is commonly used to express the particle size. (Comparisons with other techniques have been found to be good for spherical particles; for non-spherical particles results may differ.)

## 3. Significance and Use

3.1 This test method is useful to both sellers and purchasers of alumina and quartz powders for determining particle size distributions for materials specifications, manufacturing control, and development and research.

## 4. Apparatus

4.1 *Electrical Sensing Zone Particle Counter.*

4.2 *Aperture Tubes*, diameter ranging from approximately 30 to 140  $\mu\text{m}$ . The diameter required is dependent upon the particle size distribution of the sample. Generally any given tube will cover a particle size range from 2 to 60 % of its aperture diameter.

NOTE 1—In certain cases, apertures up to 300  $\mu\text{m}$  are usable.

4.3 *Sample Beaker*, capable of maintaining all particles uniformly in suspension (for example, round-bottom).

4.4 *Blender*, capacity 1-L glass container. A means to control speed is required.

4.5 *Beakers*, 100, 500, and 1000-mL.

4.6 *Pipet.*

4.7 *Wash Bottles.*

4.8 *Membrane Filtering Device*, rated at 0.45- $\mu\text{m}$  filters or finer.

## 5. Reagents

5.1 *Purity of Reagents*—Reagent grade chemicals shall be used in all tests. Unless otherwise indicated, it is intended that all reagents shall conform to the specifications of the Committee on Analytical Reagents of the American Chemical Society, where such specifications are available.<sup>2</sup> Other grades may be used, provided it is first ascertained that the reagent is of sufficiently high purity to permit its use without lessening the accuracy of the determination.

5.2 *Dispersing Media*—Ten percent solution of purified or reagent grade sodium hexametaphosphate in distilled water twice filtered through the membrane filtering device.

NOTE 2—Deionized water may be substituted for distilled water.

NOTE 3—This liquid should not be retained longer than 1 month and should not be pH modified or heated.

<sup>1</sup> This test method is under the jurisdiction of ASTM Committee C21 on Ceramic Whitewares and Related Products and is the direct responsibility of Subcommittee C21.04 on Raw Materials.

Current edition approved Dec. 1, 2014. Published December 2014. Originally approved in 1971. Last previous edition approved in 2009 as C690 – 09. DOI: 10.1520/C0690-09R14.

<sup>2</sup> *Reagent Chemicals, American Chemical Society Specifications*, American Chemical Society, Washington, DC. For suggestions on the testing of reagents not listed by the American Chemical Society, see *Analar Standards for Laboratory Chemicals*, VWR International Ltd., U.K., and the *United States Pharmacopoeia*, USPC, Rockville, MD..

5.3 *Electrolyte*—Dissolve 10.0 g of reagent grade sodium chloride (NaCl) in 1000 mL of distilled water and filter twice through the membrane filtering device.

5.4 *Wash Water*—Distilled water twice filtered through the membrane filtering device.

5.5 *Calibration Particles*—NIST or NIST traceable monosized particle standards.

## 6. Procedure

6.1 *Summary*—Disperse the test powder in the electrolyte with a blender. Transfer a representative portion to the sample beaker that contains filtered electrolyte. Place sample beaker in the apparatus and obtain particle size distribution in a chosen size range. Obtain relative weight fraction by assuming constant particle density.

6.2 Precalibrate the aperture and electrolyte combination following the manufacturer's instruction manual.

NOTE 4—Calibration should be performed in accordance with the instruction manual. Monosized NIST or NIST traceable calibration standards should be selected from Fig. A1.1. Mutual agreement on the source and size of calibration standards is necessary for interlaboratory comparisons.

6.3 Check background counts by filling the sample beaker with filtered electrolyte and taking counts without any sample added. Follow 6.6, 6.7, and 6.8.

6.4 Disperse approximately 0.7 g of sample in 200 mL of electrolyte containing 5 drops of dispersing media, by mixing at high speed on the blender or its equivalent for 5 min.

NOTE 5—The proper dispersion conditions for a given mixer or blender should be predetermined by obtaining a time-speed versus median diameter curve (see typical curve in Fig. A1.2) while ensuring that grinding does not occur. The position of the plateau will indicate the proper dispersion conditions for the sample. Experience has shown that full speed on the Waring Blender may cause size reduction. Slightly less than full speed should be used. For some suspensions ultrasonic treatment from 1 to 5 min is effective.

6.5 With a pipet, transfer an appropriate aliquot of dispersed sample into the sample beaker containing electrolyte with dispersing media added in the ratio of 3 drops/200 mL of electrolyte. The aliquot size is dependent on the aperture size used. Wash down the pipet by rinsing with electrolyte several times (see 6.9.3).

NOTE 6—The blender or mixer should be stirring just rapidly enough to maintain a uniform particle suspension while withdrawing the sample. The pipet should deliver all of the withdrawn slurry to ensure a representative transfer of sample in the event of any size classification during the transfer.

6.6 Place the sample beaker in position on the sample stand.

6.7 Adjust the speed of the stirrer to furnish sufficient agitation to maintain a uniform particle suspension, but below air bubble generation speeds.

6.8 Use the apparatus control software to set the measurement parameters. Make three measurements in which each measurement counts and measures at least 5000 particles. Average the particle size distribution from the three measurements and report the statistical parameters from the averaged results.

### 6.9 Precautions:

6.9.1 Before each analysis, using wash bottle and filtered wash water, wash all surfaces coming in contact with sample.

6.9.2 Ensure that the calibration of the instrument is correct by checking the calibration factor at least once a week.

6.9.3 The number of particles per unit volume in the sample beaker should not exceed that which will give a 5 % coincidence correction for the aperture tube being used (see Fig. A1.1).

## 7. Presentation of Data

7.1 Convert data to cumulative weight percent greater than stated particle size according to instrument instruction manual. Coincidence is insignificant if total counts are limited to Fig. A1.1.

NOTE 7—For all electrical sensing zone counters the conversion is actually to volume percent. If all particles in the sample have the same density the volume percent and weight percent are interchangeable.

7.2 Report size distribution graphs, tables, and statistics such as: weight percent, count percent, volume percent, mean, median, mode, quartiles, and standard deviation.

## 8. Precision and Bias

8.1 *Intralaboratory, Same Operator*—Experience of several laboratories indicates that the test method is capable of a precision of  $\pm 1$  % (95 % confidence level) for all size values.

8.2 *Interlaboratory*—Experience of several laboratories indicates that the test method is capable of a precision of  $\pm 3$  % (95 % confidence level) for all size values.

8.3 *Bias*—Instrument calibrations shall be performed using NIST or NIST traceable uniform spheres with relative standard deviation of 5 % or less.

## 9. Keywords

9.1 alumina; particle size; quartz; sensing zone

ANNEX

(Mandatory Information)

A1. APERTURE CHART AND DISPERSION CURVE

Figs. A1.1 and A1.2 are examples of the charts that should be employed in conjunction with this test method.

Nominal Aperture Size, μm	Total Cumulative Count for 5 % Coincidence Correction			Approximate Particle Size Range in Equivalent Spheri- cal Diameter	Suggested Calibra- tion Particles <sup>A</sup>
	Manometer Volume				
	50 μL	500 μL	2000 μL	μm	Diameter, μm
30	74 100	<sup>A</sup>	<sup>A</sup>	0.6 to 12.0	1.5 to 6.0
50	15 950	159 500	<sup>A</sup>	1.0 to 20.0	2.5 to 10.0
70	5 810	58 100	232 400	1.4 to 28.0	3.5 to 14.0
100	2 000	20 000	80 000	2.0 to 40.0	5.0 to 20.0
140	725	7 250	29 000	2.8 to 56.0	7.0 to 28.0
200	<sup>A</sup>	2500	10 000	4.0 to 80.0	10.0 to 40.0
280	<sup>A</sup>	<sup>A</sup>	3 630	5.6 to 112.0	14.0 to 56.0

<sup>A</sup> Aperture size and manometer volume combination not recommended.

FIG. A1.1 Typical Aperture Chart

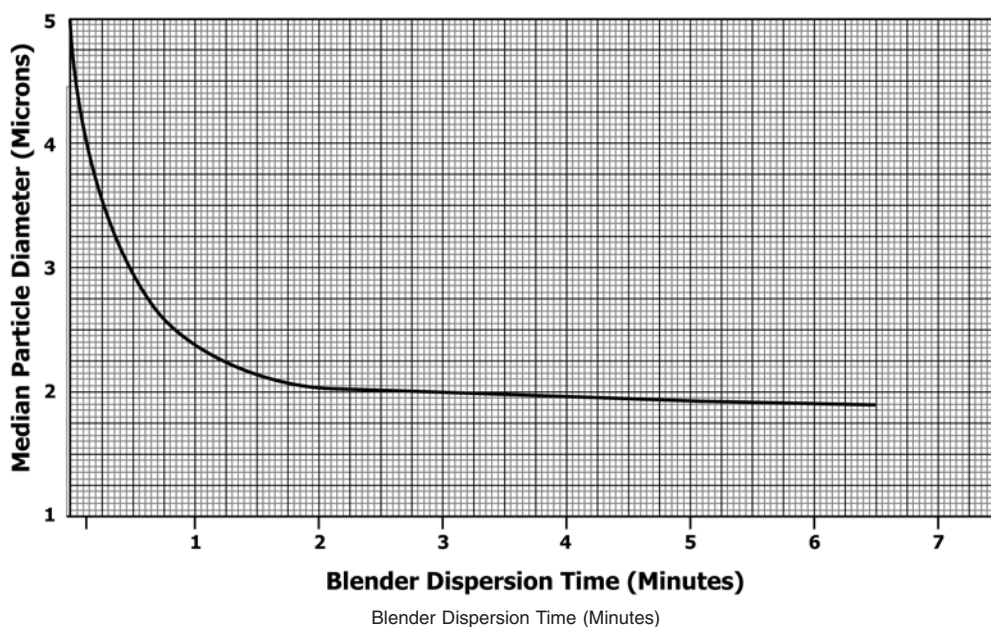


FIG. A1.2 Example of a Sample Dispersion Curve

SUMMARY OF CHANGES

Committee C21 has identified the location of selected changes to this standard since the last issue (C690–03) that may impact the use of this standard.

(1) Instruments are available from multiple suppliers and thus reference to particular instruments in sections 4.1 and 4.8, along with footnotes 2 and 3 must be removed according to ASTM policy. All other footnotes are renumbered accordingly. (2) **Note 1** is revised to indicate aperture openings up to 300µm may be used, instead of 280µm, since 300µm are now available with performance essentially the same as 280µm tubes.

(3) The word “Nominal” has been added to the heading for the first column in **Fig. A1.1** since aperture tubes are available in sizes slightly different from those given in the table.

*ASTM International takes no position respecting the validity of any patent rights asserted in connection with any item mentioned in this standard. Users of this standard are expressly advised that determination of the validity of any such patent rights, and the risk of infringement of such rights, are entirely their own responsibility.*

*This standard is subject to revision at any time by the responsible technical committee and must be reviewed every five years and if not revised, either reapproved or withdrawn. Your comments are invited either for revision of this standard or for additional standards and should be addressed to ASTM International Headquarters. Your comments will receive careful consideration at a meeting of the responsible technical committee, which you may attend. If you feel that your comments have not received a fair hearing you should make your views known to the ASTM Committee on Standards, at the address shown below.*

*This standard is copyrighted by ASTM International, 100 Barr Harbor Drive, PO Box C700, West Conshohocken, PA 19428-2959, United States. Individual reprints (single or multiple copies) of this standard may be obtained by contacting ASTM at the above address or at 610-832-9585 (phone), 610-832-9555 (fax), or service@astm.org (e-mail); or through the ASTM website (www.astm.org). Permission rights to photocopy the standard may also be secured from the Copyright Clearance Center, 222 Rosewood Drive, Danvers, MA 01923, Tel: (978) 646-2600; <http://www.copyright.com/>*