

Standard Test Method for Temperature and Enthalpy Measurement Validation of Differential Scanning Calorimeters¹

This standard is issued under the fixed designation E2253; the number immediately following the designation indicates the year of original adoption or, in the case of revision, the year of last revision. A number in parentheses indicates the year of last reapproval. A superscript epsilon (ε) indicates an editorial change since the last revision or reapproval.

1. Scope*

- 1.1 This test method provides procedures for validating temperature and enthalpy measurements of differential scanning calorimeters (DSC) and analytical methods based upon the measurement of temperature or enthalpy (or heat), or both, by DSC. Performance parameters determined include temperature and calorimetric repeatability (precision), detection limit, quantitation limit, linearity, and bias. This test method is applicable to exothermic events.
- 1.2 Validation of apparatus performance and analytical methods is requested or required for quality initiatives or where results may be used for legal purposes.
- 1.3 The values stated in SI units are to regarded as standard. No other units of measurement are included in this standard.
- 1.4 This standard does not purport to address all of the safety concerns, if any, associated with its use. It is the responsibility of the user of this standard to establish appropriate safety and health practices and determine the applicability of regulatory limitations prior to use.

2. Referenced Documents

2.1 ASTM Standards:²

E177 Practice for Use of the Terms Precision and Bias in ASTM Test Methods

E473 Terminology Relating to Thermal Analysis and Rheology

E967 Test Method for Temperature Calibration of Differential Scanning Calorimeters and Differential Thermal Analyzers

E968 Practice for Heat Flow Calibration of Differential Scanning Calorimeters

E1142 Terminology Relating to Thermophysical Properties

E1860 Test Method for Elapsed Time Calibration of Thermal Analyzers

E1970 Practice for Statistical Treatment of Thermoanalytical Data

E2161 Terminology Relating to Performance Validation in Thermal Analysis and Rheology

2.2 FDA Publications:³

Q2B Validation of Analytical Procedures Methodology, 62 FR 27464, May 19, 1997

3. Terminology

3.1 Technical terms used in this standard are defined in Practice E177 and in Terminologies E473, E1142, and E2161 including analyte, detection limit, differential scanning calorimetry, enthalpy, extrapolated onset value (temperature), linearity, mean, precision, quantitation limit, relative standard deviation, repeatability, standard deviation, thermal curve, and validation.

4. Summary of Test Method

- 4.1 Temperature and time are the primary independent parameters and heat flow is the primary dependent experimental parameter provided by DSC. Integration of heat flow, as a function of time, yields enthalpy (heat).
- 4.1.1 Time, measured by the DSC apparatus, shall conform to better than $0.1\,\%$ verified by Test Method E1860 and reported.
- 4.1.2 Temperature is directly measured by a temperature sensor that is an integral part of the differential scanning calorimetry apparatus.
- 4.1.3 Heat flow, a measured value, is validated by its integration over time to obtain the desired calorimetric (enthalpic) information of interest.
- 4.2 Calorimetric validation of a differential scanning calorimetric apparatus at a single temperature is performed using the indium metal melt as an analyte.
- 4.3 Validation of a DSC method based upon enthalpic measurement may be performed using the test specimen as the analyte.

¹ This test method is under the jurisdiction of ASTM Committee E37 on Thermal Measurements and is the direct responsibility of Subcommittee E37.10 on Fundamental, Statistical and Mechanical Properties.

Current edition approved April 1, 2016. Published April 2016. Originally approved in 2003. Last previous edition approved in 2011 as E2253 – 11. DOI: 10.1520/E2253-16.

² For referenced ASTM standards, visit the ASTM website, www.astm.org, or contact ASTM Customer Service at service@astm.org. For *Annual Book of ASTM Standards* volume information, refer to the standard's Document Summary page on the ASTM website.

³ Available from Food and Drug Administration (FDA), 10903 New Hampshire Ave., Silver Spring, MD 20993-0002, http://www.fda.gov.

- 4.4 The enthalpy of three (or more) specimens (nominally representing the maximum, midpoint and minimum of the range of the test method) are measured in triplicate (or more). A fourth blank specimen, containing no analyte, is also measured in triplicate.
- Note 1—Repeatability is determined by performing a sufficient number of determinations to calculate statistically valid estimates of the standard deviation or relative standard deviation of the measurement.
- 4.4.1 Calorimetric linearity and bias are determined from the best-fit (linear regression) straight-line correlation of the results from measurements of the three (or more) specimens.
- 4.4.2 Calorimetric detection limit and quantitation limit are determined from the standard deviation of the blank determination.
- 4.4.3 Calorimetric repeatability is determined from the repeatability measurement of the three (or more) specimens.
- 4.5 The temperature validation of a differential scanning calorimetric apparatus is performed at three temperatures using indium, bismuth, and zinc metal melts as an analyte (see Section 7).
- 4.5.1 The melting temperature of three (or more) materials representing the maximum, midpoint, and minimum of the temperature range of the test method are measured in triplicate (or more) (see Note 1).
- 4.5.2 Temperature linearity and bias are determined from the best-fit (linear regression) straight-line correlation of the results from the temperature measurements at the three (or more) temperatures.
- 4.5.3 Temperature repeatability, detection limit and quantitation limit are determined from the standard deviation of the replicate temperature measurements.

5. Significance and Use

- 5.1 This test method may be used to determine and validate the performance of a particular DSC apparatus.
- 5.2 This test method may be used to determine and validate the performance of a particular method based upon a DSC temperature and enthalpic measurement, or both.
- 5.3 This test method may be used to determine the repeatability of a particular apparatus, operator, or laboratory.
- 5.4 This test method may be used for specification and regulatory compliance purposes.

6. Apparatus

- 6.1 Differential Scanning Calorimeter (DSC)—The essential instrumentation required to provide the minimum differential scanning calorimetric capability for this test method includes:
 - 6.1.1 DSC Test Chamber, composed of:
- 6.1.1.1 *Furnace(s)*, to provide uniform controlled heating or cooling of a specimen and reference to a constant temperature or at a constant rate within the applicable temperature range of the test method.
- 6.1.1.2 *Temperature Sensor*, to provide an indication of the specimen temperature to readability required.
- 6.1.1.3 *Differential Sensor*, to detect a heat flow difference between the specimen and reference.

- 6.1.1.4 A means of sustaining a test chamber environment of an inert purge gas at a rate of 10 to 50 mL/min \pm 10 %.
- Note 2—Typically, 99+% pure nitrogen (or other inert gas, such as argon or helium) is employed when oxidation in air is a concern. Unless the effects of moisture are to be studied, the use of a dry purge gas is recommended, especially for operation at subambient temperatures.
- 6.1.2 *Temperature Controller*, capable of executing a specific temperature program by operating the furnace(s) between selected temperature limits at a rate of temperature change constant to $\pm 1~\%$ or at an isothermal temperature constant to $\pm 0.5^{\circ}\text{C}$.
- 6.1.3 *Data Collection Device*, to provide a means of acquiring, storing, and displaying measured or calculated signals, or both. The minimum output signals required for differential scanning calorimetry are heat flow, temperature, and time.
- 6.1.4 *Containers*, (pans, crucibles, vials, lids, closures, seals, etc.) that are inert to the specimen and reference materials and that are of suitable structural shape and integrity to contain the specimen and reference in accordance with the specific requirements of the test method.
- 6.2 *Balance*, of 100 mg or greater capacity to weigh specimens and containers to $\pm 1 \mu g$.

Note 3—A balance of this high precision is required so that weighing imprecision is not part of the overall method imprecision.

7. Reagents and Materials

- 7.1 *Indium (In) Metal*, 99.99+ % purity, preferably a certified reference material for which the melting temperature and enthalpy of fusion are known.
- 7.2 *Bismuth (Bi) Metal*, 99.99+ % purity, preferably a certified reference material for which the melting temperature is known.
- 7.3 Zinc (Zn) Metal, 99.99+ % purity, preferably a certified reference material for which the melting temperature is known.

8. Calibration and Standardization

- 8.1 After turning the power on, allow the instrument to equilibrate for at least one hour prior to any measurements.
- 8.2 Perform any cleaning and calibration procedures described by the manufacturer in the apparatus Operator's Manual.
- 8.3 If not previously established, perform temperature and heat flow calibrations according to Test Method E967 and Practice E968 respectively, using the same purge gas, purge gas flow rate and heating rate (here 10°C/min) to be used for validation experiments.
- 8.4 If not previously established, obtain the instrument's elapsed time conformance using Test Method E1860.

9. Procedure for Determining Calorimetric Repeatability, Detection Limit, Quantitation Limit, Linearity, and Bias

9.1 This process involves characterizing, in triplicate, a blank and three (or more) test specimens taken to represent the low, medium and high extremes of the range over which performance is to be validated.

Note 4—The details of this procedure are written using indium as an analyte. For validation of an enthalpic method, test specimens representing the range of that method shall be used, and steps 9.2 to 9.7 replaced with the enthalpic method procedure.

9.2 Prepare three (or more) indium test specimens covering the enthalpy or mass range of the tests. Nominal mass values might be 1, 10, and 20 mg. Measure the mass of each of these specimens to the nearest 1 μ g and record as M_{min} , M_{mid} , and M_{max} . Enclose each test specimen within clean specimen containers and lids. Also prepare a blank specimen that contains no analyte but otherwise is similar to the specimens prepared above.

Note 5—Most thermoanalytical methods cover 1.5 to 2 decades of range. The mass values selected should be approximately equally distributed over the anticipated range. Other masses and mass ranges may be used but shall be reported.

- 9.3 Load the largest specimen into the instrument chamber, purge the chamber with dry nitrogen (or other inert gas) at a flow rate of 10 to 50 mL/min \pm 10 % throughout the experiment. An empty sample specimen container is loaded in the reference position.
- 9.4 Erase any thermal history in the test specimen by heating the specimen to 180°C, then cool at 5°C/min to 120°C. The thermal curve need not be recorded.
 - 9.5 Equilibrate at 120°C for one minute.
- 9.6 Heat the test specimen at 10° C/min through the indium melting transition to 180° C and record the thermal curve.

Note 6—Other heating rates may be used but shall be reported.

- 9.7 Cool the test specimen to 120°C at 5°C/min, then cool to ambient temperature at any convenient rate. The thermal curve need not be recorded.
- 9.8 Construct a baseline for the melting endotherm by selecting a point on the curve immediately before and another immediately after the endotherm. Record the temperatures of

these two points as T_1 and T_2 . Construct a linear baseline between the two points (see Fig. 1).

- 9.9 Integrate, as a function of time, the heat flow described by the constructed baseline and the melting endotherm. Record this value as enthalpy ($\Delta Q_{max}(1)$, in mJ).
- 9.10 Repeat steps 9.3 through 9.9 for the medium mass test specimen from step 9.2. Use the same integration limits (T_1 and T_2) determined in step 9.8. Record this value as enthalpy ($\Delta Q_{mid}(1)$, in mJ).

Note 7—Loading and unloading of the specimen is required to determine analytical repeatability. If only instrumental repeatability is being determined, the specimen may be left in place between determinations

- 9.11 Repeat steps 9.3 through 9.9 for the small mass test specimen from step 9.2. Use the same integration limits (T_1 and T_2) determined in step 9.8. Record this value as enthalpy ($\Delta Q_{min}(1)$, in mJ).
- 9.12 Repeat steps 9.3 through 9.9 for the blank test specimen from step 9.2. Use the same integration limits (T_1 and T_2) determined in step 9.8 (see Fig. 2). Record this value as enthalpy ($\Delta Q_o(1)$, in mJ).

Note 8—Observe and record the sign of the value for Δ Q_o . It may be positive or negative.

- 9.13 Repeat steps 9.5 through 9.9 two more times for the large mass specimen. Remove the specimen from the DSC sample chamber and reload it between each determination. Record these values as enthalpy ($\Delta Q_{max}(2)$ and $\Delta Q_{max}(3)$, in mJ).
- 9.14 Repeat steps 9.5 through 9.9 two more times for the medium mass specimen. Remove the specimen from the DSC sample chamber and reload it between each determination. Record these values as enthalpy $(\Delta Q_{mid}(2))$ and $\Delta Q_{mid}(3)$, in m.l.


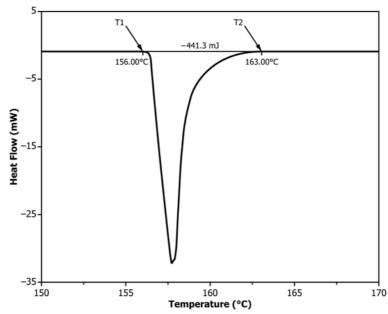


FIG. 1 Integration of Large Indium Melting Endotherm

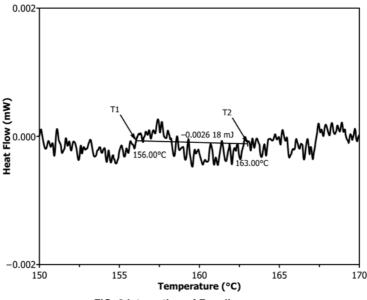

- FIG. 2 Integration of Baseline
- 9.15 Repeat steps 9.5 through 9.9 two more times for the low mass specimen. Remove the specimen from the DSC sample chamber and reload it between each determination. Record these values as enthalpy $(\Delta Q_{min}(2))$ and $\Delta Q_{min}(3)$, in mJ).
- 9.16 Repeat steps 9.5 through 9.9 two more times for the blank specimen. Remove the specimen from the DSC sample chamber and reload it between each determination. Record these values as enthalpy ($\Delta Q_o(2)$ and $\Delta Q_o(3)$, in mJ).
- 9.17 Calculate the mean (Q), standard deviation (s) and relative standard deviation (RSD) for the enthalpy from the replicate determinations made on each of the three mass specimens (see Practice E1970). Record these values as Q_{max} , Q_{mid} , Q_{min} , s_{max} , s_{mid} , s_{min} , and RSD_{max} , RSD_{mid} , and RSD_{min} , respectively. Calculate the mean and standard deviation for the enthalpy determination of the blank. Record these values as Q_o and s_o respectively.
- 9.18 Using the standard deviation for the enthalpy of the blank (s_o) from step 9.17, determine the instrument detection limit (DL) and quantitation limit (QL) as in calculations 11.2 and 11.3. Report DL and QL.
- 9.19 Calculate the pooled relative standard deviation for the enthalpy of fusion from the RSD_{max} , RSD_{mid} , and RSD_{min} , obtained in step 9.17 (see Practice E1970). Report this value as the Calorimetric Repeatability value (r).
- 9.20 Using the three (or more) mass values from step 9.2 as the independent (*X*) values and the three (or more) mean enthalpy values from step 9.17 as the dependent (*Y*) values, determine the least square (linear regression) best fit values for the slope (*m*) in mJ mg⁻¹, and intercept (*b*) in mJ (see Practice E1970).
- Note 9—The units of mJ mg $^{-1}$ for slope are equivalent to the units of J g $^{-1}$.
- 9.21 Calculate the Percent Linearity (L) of the enthalpy of fusion data from the values in step 9.20 and the calculation of Eq 5.

9.22 Compare the slope (mean enthalpy of fusion) from step 9.20 to the accepted reference value for enthalpy of fusion as in calculation of Eq 7 (see Practice E968). Report this percentage value as the Calorimetric Bias (*Bias*).

Procedure for Determining Temperature Repeatability, Detection Limit, Quantitation Limit, Linearity, and Bias

10.1 This process involves characterizing, in triplicate, three (or more) test materials taken to represent the low, medium, and high extremes of the range over which performance is to be validated.

Note 10—The details of this procedure are written using indium, bismuth, and zinc as analytes. For validation of other temperature ranges, test specimens representing the range of that method shall be used and 10.2 to 10.26 replaced with the temperature method procedure. Materials with sufficiently high temperature accuracy are presented in Table 1. The mass of reference materials shall be selected to approximate equal enthalpies. Where practical, analytes with approximately equality temperature spaced transitions shall be used.

10.2 Prepare three (or more) indium test specimens.

10.2.1 Into a tared specimen container, weigh a single 9 to 11 mg of the indium reference material. Seal the specimen container.

TABLE 1 Melting Temperatures^A

Metal	Symbol	Melting Temperature, °C
Indium	In	156.598
Tin	Sn	231.928
Bismuth	Bi	271.442
Lead	Pb	327.502
Zinc	Zn	419.527
Aluminum	Al	660.32 ^B

A Mangum, B.W., Furukawa, G.T., "Guidelines for Realizing the International Practical Temperature Scale of 1990 (ITS-90)", United States National Institute of Technology (NIST) Technical Note 1265, 1990.

^B Should aluminum be used as a reference material, specimen containers other than aluminum shall be used for all test materials.

10.2.2 Repeat 10.2.1 for two additional indium test specimens.

Note 11—The test 10 mg specimen from 9.2 may be used for this purpose.

- 10.3 Prepare three (or more) bismuth test specimens.
- 10.3.1 Into a tared specimen container, weigh a single 4.8 to 5.8 mg of the bismuth reference material. Seal the specimen container.
- 10.3.2 Repeat 10.3 for two additional bismuth test specimens.
 - 10.4 Prepare three (or more) zinc test specimens.
- 10.4.1 Into a tared specimen container, weigh a single 2.4 to 2.8 mg of the zinc reference material. Seal the specimen container.

Note 12—The masses for indium, bismuth, and zinc specimens in steps 10.2 to 10.4 are selected to provide approximately the same melting endothermic peak (valley) heights.

Note 13—The temperatures of the materials selected in steps 10.2 to 10.4 should be approximately equally distributed over the anticipate temperature range. Other materials and masses may be used but shall be reported.

- 10.5 Load an indium (or material with the lowest transition temperature) specimen into the instrument chamber, purge the chamber with dry nitrogen (or other inert gas) at a flow rate of 10 to $50 \, \text{mL/min} \pm 10 \, \%$ throughout the experiment. An empty sample specimen container is loaded in the reference position.
 - 10.6 Equilibrate at 120°C for one minute.
- 10.7 Heat the test specimen at 10°C/min through the indium melting transition to 180°C and record the thermal curve.

Note 14—Other heating rates may be used but shall be reported.

10.8 Cool the test specimen to ambient temperature at any convenient rate. The thermal curve need not be recorded.

- 10.9 Display the thermal curve for the melting endotherm with heat flow on the ordinate (Y-axis) and specimen temperature on the abscissa (X-axis). Determine the extrapolated onset temperature for the indium melting endotherm as T(In1) to two places to the right of the decimal point.
- 10.9.1 Construct a baseline by connecting a straight line from a point on the thermal curve immediately before and another immediately after the melting endotherm.
- 10.9.2 Construct a tangent at the point of maximum rate of change on the leading edge of the melting endotherm. Extrapolate this line to the baseline in 10.9.1.
- 10.9.3 The temperature of intersection of the baseline (in 10.9.1) and the tangent (in 10.9.2) is taken as the extrapolated onset temperature (see Fig. 3).
- 10.10 Load a bismuth specimen (or material with the mid-range transition temperature) from 10.3 into the instrument chamber, purge the chamber with dry nitrogen (or other inert gas) at a flow rate of 10 to $50~\text{mL/min} \pm 10~\%$ throughout the experiment. An empty specimen container is loaded in the reference position.
 - 10.11 Equilibrate at 250°C for one minute.
- 10.12 Heat the test specimen at 10°C/min through the bismuth melting transition to 290°C and record heat flow and specimen temperature (see Note 14).
- 10.13 Cool the test specimen to ambient temperature at any convenient rate. The thermal curve need not be recorded.
- 10.14 Display the thermal curve for the melting endotherm with heat flow on the ordinate (Y-axis) and specimen temperature on the abscissa (X-axis). Determine the extrapolated onset temperature for the bismuth melting endotherm as T(Bil) to two places to the right of the decimal point (see 10.9.1 to 10.9.3).

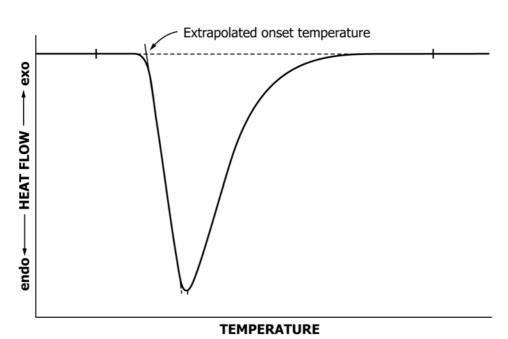


FIG. 3 Determination of the Extrapolated Onset Temperature

10.15 Load a zinc (or material with the highest transition temperature) specimen into the instrument chamber, purge the chamber with dry nitrogen (or other purge gas) at a flow rate of 10 to $50 \,\mathrm{mL/min} \pm 10 \,\%$ throughout the experiment, An empty sample specimen container is loaded in the reference position.

10.16 Equilibrate at 400°C for one minute.

10.17 Heat the test specimen at 10°C/min through the zinc melting temperature to 440°C and record the heat flow, specimen temperature and time (Note 14).

10.18 Cooling the test specimen to ambient temperature at any convenient rate. The thermal curve need not be recorded.

10.19 Display the thermal curve for the melting endotherm with heat flow on the ordinate (Y-axis) and specimen temperature on the abscissa (X-axis). Determine the extrapolated onset temperature for the zinc melting endotherm as T(Zn1) to two places to the right of the decimal point (see 10.9.1 to 10.9.3).

10.20 Repeat steps 10.5 through 10.19 two more times each for the low, medium, and high temperature materials. Remove the specimen from the DSC sample chamber and reload with a fresh specimen for each determination. Record the melting temperature values as T(In2), T(Bi2), T(Zn2), T(In3), T(Bi3) and T(Zn3) in °C.

10.21 Calculate (see Practice E1970) the mean (T) and standard deviation (σ) for the melting temperatures for indium, bismuth, and zinc as T(In), T(Bi), T(Zn), $\sigma(In)$, $\sigma(Bi)$, and $\sigma(Zn)$, respectively.

10.22 Pool (see Practice E1970) the three standard deviation values from 10.21 and record as s.

10.23 Using the pooled standard deviation s from 10.22, determine the detection limit (DL) and quantitation limit (QL)as calculated in Eq 2 and Eq 4. Report Temperature DL and OL.

10.24 Using the three literature values for the melting temperature of indium, bismuth and zinc (for example, from Table 1) as the independent (X) value and the three (or more) mean values for the observed melting temperatures values from 10.9, 10.10, 10.14, and 10.19 as the dependent (Y) values, determine the least-squares (linear regression) best-fit value for the slope (m), and intercept (b) (see Practice E1970).

Note 15—The value for m will be dimensionless, and the value for bwill be in °C.

10.25 Calculate the Percent Linearity (L) of the temperature determination from the values in 10.24 and the calculation of Eq 5.

10.26 Compare the slope (m) from 10.24 to the anticipated value of 1.00000). Report this value as the Temperature Bias (TBias) as a percentage according to Eq 7.

11. Calculations

11.1 When performing these calculations, retain all available decimal places in the measured values and in intermediate calculated values. The final result should be rounded to three significant figures.

11.2 Detection Limit (DL) is given by:

$$DL = 3.3 s_a \tag{1}$$

$$DL = 3.3 s \tag{2}$$

where:

DL = detection limit, mJ or °C, and

= blank enthalpy standard deviation, mJ, and

= temperature standard deviation, °C.

Note 16—To express the Detection Limit in terms of mass of analyte, the value for DL may be divided by the slope value m from 9.20.

11.3 Quantitation Limit (QL) is given by:

$$QL = 10 s_a \tag{3}$$

$$QL = 10 s \tag{4}$$

where:

QL = quantitation limit, mJ or °C.

Note 17-To express the Quantitation Limit in terms of mass of the analyte, the value for QL may be divided by the slope value m from 9.20.

11.4 Percent Linearity (L) is given by:

$$L = \left[\frac{\left(100 \% \times | \text{ Largest } \delta Y | \right)}{m \times M_{max} + b} \right]$$
 (5)

$$L = \left[\frac{\left(100 \% \times | \text{Largest } \delta Y | \right)}{m \times M_{max} + b} \right]$$

$$L = \left[\frac{\left(100 \% \times | \text{Largest } \delta Y | \right)}{\left(m \times \left(T(\text{Zn}) - T(\text{In}) \right) + b \right)} \right]$$
(6)

where:

L= linearity, %,

mg⁻¹, (from = slope, mJ dimensionless, (from 10.23),

= intercept, mJ, (from 9.20), or dimensionless, from 10.23).

 M_{max} = mass of largest test specimen, mg, or highest temperature, °C, (from 10.21),

|Largest δY | = absolute value of the largest deviation from best fit straight line,

T(Zn)= melting temperature of zinc, and T(In)= Melting temperature of indium.

11.5 Calorimetric and temperature Bias (Bias) is given by:

$$Bias = \frac{\left(\Delta H_{ref} - m\right) \times 100\%}{\Delta H_{ref}} \tag{7}$$

$$Bias = [(m - 1.00000) \times 100 \%]$$
 (8)

where:

 ΔH_{ref} = accepted reference value for indium heat of fusion, 28.58 mJ mg⁻¹.4

12. Example Calculations

12.1 For the example calculations described below, the following experimental data is used:

$M_{max} = 15.680 \text{ mg}$	$Q_{max} = 444.33 \text{ mJ}$
$M_{mid} = 8.000 \text{ mg}$	$Q_{mid} = 233.33 \text{ mJ}$
$M_{min} = 0.376 \text{ mg}$	$Q_{min} = 11.133 \text{ mJ}$
$M_o = 0.000 \text{ mg}$	$Q_o = 0.0019093 \text{ mJ}$
$n_{max} = 3$	$s_{max} = 6.4470 \text{ mJ}$
$n_{mid} = 3$	$s_{mid} = 0.14012 \text{ mJ}$
$n_{min} = 3$	$s_{min} = 0.089629 \text{ mJ}$
$n_o = 3$	$s_o = 0.0041285 \text{ mJ}$

⁴ Stølen, S., Grønvold, F., Thermochimica Acta, Vol 327, 1999, p. 1.

where:

n = number of replicate determinations.

12.2 Example calculation for Detection Limit:

$$DL = 3.3 \times 0.0041285 \text{ mJ} = 0.0136 \text{ mJ}$$
 (9)

12.3 Example calculation for Quantitation Limit:

$$QL = 10 \times 0.0041285 \text{ mJ} = 0.0413 \text{ mJ}$$
 (10)

12.4 Example calculation for Calorimetric Repeatability:

$$RSD_{max} = 6.4470 \text{ mJ} \times 100 \% / 444.33 \text{ mJ} = 1.4509 \%$$
 (11)

$$RSD_{mid} = 0.14012 \text{ mJ} \times 100 \% / 233.33 \text{ mJ} = 0.060052 \%$$

$$RSD_{min} = 0.089629 \text{ mJ} \times 100 \% / 11.133 \text{ mJ} = 0.80508 \%$$

$$r = \left[\frac{2 \times 1.4509\%^2 + 2 \times 0.060052\%^2 + 2 \times 0.80508\%^2}{2 + 2 + 2}\right]^{1/2}$$

$$r = \left[\frac{4.2102 + 0.0072 + 1.2963}{6}\right]^{1/2} = \left[\frac{5.5137}{6}\right]^{1/2}$$

$$=0.91895^{1/2}=0.959\%$$

- 12.5 Example calculation of Percent Linearity:
- 12.5.1 From the best fit to the *max*, *mid*, and *min* data of 12.1:

$$m = 28.305 \text{ mJ mg}^{-1} = 28.305 \text{ J g}^{-1}$$
 (12)

$$b = 2.6284 \text{ mJ}$$

12.5.2 From the best fit straight line:

$$\delta Y_{max} = 2.123 \text{ mJ} \tag{13}$$

 $\delta Y_{mid} = -4.261 \text{ mJ} \leftarrow \text{Largest deviation from the best fit line}$

$$\delta Y_{min} = 2.138 \text{ mJ}$$

12.5.3 Example calculation:

$$L = \left[\frac{100\% \times 4.261 \text{ mJ}}{28.305 \text{ mJ} \cdot \text{mg}^{-1} \times 15.68 \text{ mg} + 2.6284 \text{ mJ}} \right]$$
(14)
$$= \left[\frac{426.2 \text{ mJ}}{446.451 \text{ mJ}} \right] = 0.954\%$$

12.6 Example calculation for Calorimetric Bias:

Bias =
$$\frac{(28.58 \text{ mJ mg}^{-1} - 28.305 \text{ mJ mg}^{-1}) \times 100 \%}{28.58 \text{ mJ mg}^{-1}}$$
 (15)
= 0.962 %

13. Report

- 13.1 Manufacturer and model of the DSC apparatus.
- 13.2 Time conformity from Test Method E1860.
- 13.3 Temperature calibration from Test Method E967.
- 13.4 Heat flow calibration from Practice E968.
- 13.5 Some or all of the performance parameters described by this test method useful for the application at hand, including:
 - 13.5.1 Range (Q_{min} to Q_{max} , or M_{min} to M_{max}),
 - 13.5.2 Temperature and Calorimetric Detection Limit (DL),
- 13.5.3 Temperature and Calorimetric Quantitation Limit (QL),
 - 13.5.4 Temperature and Calorimetric Percent Linearity (L),
 - 13.5.5 Temperature and Calorimetric Bias (Bias), and
 - 13.5.6 Calorimetric Repeatability (r).
 - 13.6 The specific dated version of this test method used.

14. Precision and Bias

14.1 Precision and bias are determined by the procedures in this document. No separate precision and bias information is required.

15. Keywords

15.1 detection limit; differential scanning calorimetry; enthalpy; linearity; precision; quantitation limit; range; repeatability; thermal analysis; validation

SUMMARY OF CHANGES

Committee E37 has identified the location of selected changes to this standard since the last issue (E2253 – 11) that may impact the use of this standard. (Approved April 1, 2016.)

(1) Minor editorial changes to 3.1, 4.4.1, 4.5.2, 9.20, 9.21, Section 10 heading, and 10.24.



ASTM International takes no position respecting the validity of any patent rights asserted in connection with any item mentioned in this standard. Users of this standard are expressly advised that determination of the validity of any such patent rights, and the risk of infringement of such rights, are entirely their own responsibility.

This standard is subject to revision at any time by the responsible technical committee and must be reviewed every five years and if not revised, either reapproved or withdrawn. Your comments are invited either for revision of this standard or for additional standards and should be addressed to ASTM International Headquarters. Your comments will receive careful consideration at a meeting of the responsible technical committee, which you may attend. If you feel that your comments have not received a fair hearing you should make your views known to the ASTM Committee on Standards, at the address shown below.

This standard is copyrighted by ASTM International, 100 Barr Harbor Drive, PO Box C700, West Conshohocken, PA 19428-2959, United States. Individual reprints (single or multiple copies) of this standard may be obtained by contacting ASTM at the above address or at 610-832-9585 (phone), 610-832-9555 (fax), or service@astm.org (e-mail); or through the ASTM website (www.astm.org). Permission rights to photocopy the standard may also be secured from the Copyright Clearance Center, 222 Rosewood Drive, Danvers, MA 01923, Tel: (978) 646-2600; http://www.copyright.com/