INTERNATIONAL ISO/IEEE STANDARD

11073-10419

First edition 2016-06-15

Health informatics — Personal health device communication —

Part 10419:

Device specialization — Insulin pump

Informatique de santé — Communication entre dispositifs de santé personnels

Partie 10419: Spécialisation du dispositif — Pompe à insuline





COPYRIGHT PROTECTED DOCUMENT

© ISO 2016 © IEEE 2015

Published in Switzerland

All rights reserved. Unless otherwise specified, no part of this publication may be reproduced or utilized in any form or by any means, electronic or mechanical, including photocopying and microfilm, without permission in writing from either ISO or IEEE at the respective address below.

ISO copyright office
Case postale 56 • CH-1211 Geneva 20
Tel. + 41 22 749 01 11
Fax + 41 22 749 09 47
E-mail copyright@iso.org
Web www.iso.org

Institute of Electrical and Electronics Engineers, Inc. 3 Park Avenue, New York • NY 10016-5997, USA E-mail stds.ipr@ieee.org
Web www.ieee.org

Foreword

ISO (the International Organization for Standardization) is a worldwide federation of national standards bodies (ISO member bodies). The work of preparing International Standards is normally carried out through ISO technical committees. Each member body interested in a subject for which a technical committee has been established has the right to be represented on that committee. International organizations, governmental and non-governmental, in liaison with ISO, also take part in the work. ISO collaborates closely with the International Electrotechnical Commission (IEC) on all matters of electrotechnical standardization.

IEEE Standards documents are developed within the IEEE Societies and the Standards Coordinating Committees of the IEEE Standards Association (IEEE-SA) Standards Board. The IEEE develops its standards through a consensus development process, approved by the American National Standards Institute, which brings together volunteers representing varied viewpoints and interests to achieve the final product. Volunteers are not necessarily members of the Institute and serve without compensation. While the IEEE administers the process and establishes rules to promote fairness in the consensus development process, the IEEE does not independently evaluate, test, or verify the accuracy of any of the information contained in its standards.

The main task of technical committees is to prepare International Standards. Draft International Standards adopted by the technical committees are circulated to the member bodies for voting. Publication as an International Standard requires approval by at least 75 % of the member bodies casting a vote.

Attention is called to the possibility that implementation of this standard may require the use of subject matter covered by patent rights. By publication of this standard, no position is taken with respect to the existence or validity of any patent rights in connection therewith. ISO/IEEE is not responsible for identifying essential patents or patent claims for which a license may be required, for conducting inquiries into the legal validity or scope of patents or patent claims or determining whether any licensing terms or conditions provided in connection with submission of a Letter of Assurance or a Patent Statement and Licensing Declaration Form, if any, or in any licensing agreements are reasonable or non-discriminatory. Users of this standard are expressly advised that determination of the validity of any patent rights, and the risk of infringement of such rights, is entirely their own responsibility. Further information may be obtained from ISO or the IEEE Standards Association.

ISO/IEEE 11073-10419 was prepared by the IEEE 11073 Standards Comittee of the IEEE Engineering in Medicine and Biology Society (as IEEE Std 11073-10419-2015). It was adopted by Technical Committee ISO/TC 215, *Health informatics*, in parallel with its approval by the ISO member bodies, under the "fast-track procedure" defined in the Partner Standards Development Organization cooperation agreement between ISO and IEEE. IEEE is responsible for the maintenance of this document with participation and input from ISO member bodies.

ISO/IEEE 11073 consists of the following parts, under the general title Health informatics — Personal health device communication (text in parentheses gives a variant of subtitle):

- Part 10101: (Point-of-care medical device communication) Nomenclature
- Part 10102: (Point-of-care medical device communication) Nomenclature: Annotated ECG
- Part 10103: (Point-of-care medical device communication) Nomenclature: Implantable device, cardiac
- Part 10201: (Point-of-care medical device communication) Domain information model
- Part 10404: Device specialization Pulse oximeter
- Part 10406: Device specialization Basic electrocardiograph (ECG) (1- to 3-lead ECG)

- Part 10407: Device specialization Blood pressure monitor
- Part 10408: Device specialization Thermometer
- Part 10415: Device specialization Weighing scale
- Part 10417: Device specialization Glucose meter
- Part 10418: Device specialization International Normalized Ratio (INR) monitor
- Part 10420: Device specialization Body composition analyzer
- Part 10421: Device specialization Peak expiratory flow monitor (peak flow)
- Part 10471: Device specialization Independent living activity hub
- Part 10472: Device specialization Medication monitor
- Part 20101: (Point-of-care medical device communication) Application profiles Base standard
- Part 20601: Application profile Optimized exchange protocol
- Part 30200: (Point-of-care medical device communication) Transport profile Cable connected
- Part 30300: (Point-of-care medical device communication) Transport profile Infrared wireless
- Part 30400: (Point-of-care medical device communication) Interface profile Cabled Ethernet
- Part 90101: (Point-of-care medical device communication) Analytical instruments Point-of-care test
- Part 91064: (Standard communication protocol) Computer-assisted electrocardiography
- Part 92001: (Medical waveform format) Encoding rules [Technical Specification]

Abstract: Within the context of the ISO/IEEE 11073 family of standards for device communication, a normative definition of communication between personal telehealth insulin pump devices and compute engines (e.g., cell phones, personal computers, personal health appliances, set top boxes) in a manner that enables plug-and-play interoperability, is established in this standard. Appropriate portions of existing standards including ISO/IEEE 11073 terminology, information models, application profile standards, and transport standards are leveraged. The use of specific term codes, formats, and behaviors in telehealth environments restricting optionality in base frameworks in favor of interoperability are specified. A common core of communication functionality for personal telehealth insulin pump devices is defined.

Keywords: IEEE 11073-10419[™], insulin pump, medical device communication, personal health devices

Copyright © 2015 by The Institute of Electrical and Electronics Engineers, Inc. All rights reserved. Published 10 April 2015. Printed in the United States of America.

IEEE is a registered trademark in the U.S. Patent & Trademark Office, owned by The Institute of Electrical and Electronics Engineers, Incorporated.

PDF: ISBN 978-0-7381-9610-7 STD20158 Print: ISBN 978-0-7381-9611-4 STDPD20158

IEEE prohibits discrimination, harassment, and bullying.

For more information, visit http://www.ieee.org/web/aboutus/whatis/policies/p9-26.html.

No part of this publication may be reproduced in any form, in an electronic retrieval system or otherwise, without the prior written permission of the publisher.

The Institute of Electrical and Electronics Engineers, Inc. 3 Park Avenue, New York, NY 10016-5997, USA

Important Notices and Disclaimers Concerning IEEE Standards Documents

IEEE documents are made available for use subject to important notices and legal disclaimers. These notices and disclaimers, or a reference to this page, appear in all standards and may be found under the heading "Important Notice" or "Important Notices and Disclaimers Concerning IEEE Standards Documents."

Notice and Disclaimer of Liability Concerning the Use of IEEE Standards Documents

IEEE Standards documents (standards, recommended practices, and guides), both full-use and trial-use, are developed within IEEE Societies and the Standards Coordinating Committees of the IEEE Standards Association ("IEEE-SA") Standards Board. IEEE ("the Institute") develops its standards through a consensus development process, approved by the American National Standards Institute ("ANSI"), which brings together volunteers representing varied viewpoints and interests to achieve the final product. Volunteers are not necessarily members of the Institute and participate without compensation from IEEE. While IEEE administers the process and establishes rules to promote fairness in the consensus development process, IEEE does not independently evaluate, test, or verify the accuracy of any of the information or the soundness of any judgments contained in its standards.

IEEE does not warrant or represent the accuracy or content of the material contained in its standards, and expressly disclaims all warranties (express, implied and statutory) not included in this or any other document relating to the standard, including, but not limited to, the warranties of: merchantability; fitness for a particular purpose; non-infringement; and quality, accuracy, effectiveness, currency, or completeness of material. In addition, IEEE disclaims any and all conditions relating to: results; and workmanlike effort. IEEE standards documents are supplied "AS IS" and "WITH ALL FAULTS."

Use of an IEEE standard is wholly voluntary. The existence of an IEEE standard does not imply that there are no other ways to produce, test, measure, purchase, market, or provide other goods and services related to the scope of the IEEE standard. Furthermore, the viewpoint expressed at the time a standard is approved and issued is subject to change brought about through developments in the state of the art and comments received from users of the standard.

In publishing and making its standards available, IEEE is not suggesting or rendering professional or other services for, or on behalf of, any person or entity nor is IEEE undertaking to perform any duty owed by any other person or entity to another. Any person utilizing any IEEE Standards document, should rely upon his or her own independent judgment in the exercise of reasonable care in any given circumstances or, as appropriate, seek the advice of a competent professional in determining the appropriateness of a given IEEE standard.

IN NO EVENT SHALL IEEE BE LIABLE FOR ANY DIRECT, INDIRECT, INCIDENTAL, SPECIAL, EXEMPLARY, OR CONSEQUENTIAL DAMAGES (INCLUDING, BUT NOT LIMITED TO: PROCUREMENT OF SUBSTITUTE GOODS OR SERVICES; LOSS OF USE, DATA, OR PROFITS; OR BUSINESS INTERRUPTION) HOWEVER CAUSED AND ON ANY THEORY OF LIABILITY, WHETHER IN CONTRACT, STRICT LIABILITY, OR TORT (INCLUDING NEGLIGENCE OR OTHERWISE) ARISING IN ANY WAY OUT OF THE PUBLICATION, USE OF, OR RELIANCE UPON ANY STANDARD, EVEN IF ADVISED OF THE POSSIBILITY OF SUCH DAMAGE AND REGARDLESS OF WHETHER SUCH DAMAGE WAS FORESEEABLE.

Translations

The IEEE consensus development process involves the review of documents in English only. In the event that an IEEE standard is translated, only the English version published by IEEE should be considered the approved IEEE standard.

Official statements

A statement, written or oral, that is not processed in accordance with the IEEE-SA Standards Board Operations Manual shall not be considered or inferred to be the official position of IEEE or any of its committees and shall not be considered to be, or be relied upon as, a formal position of IEEE. At lectures, symposia, seminars, or educational courses, an individual presenting information on IEEE standards shall make it clear that his or her views should be considered the personal views of that individual rather than the formal position of IEEE.

Comments on standards

Comments for revision of IEEE Standards documents are welcome from any interested party, regardless of membership affiliation with IEEE. However, IEEE does not provide consulting information or advice pertaining to IEEE Standards documents. Suggestions for changes in documents should be in the form of a proposed change of text, together with appropriate supporting comments. Since IEEE standards represent a consensus of concerned interests, it is important that any responses to comments and questions also receive the concurrence of a balance of interests. For this reason, IEEE and the members of its societies and Standards Coordinating Committees are not able to provide an instant response to comments or questions except in those cases where the matter has previously been addressed. For the same reason, IEEE does not respond to interpretation requests. Any person who would like to participate in revisions to an IEEE standard is welcome to join the relevant IEEE working group.

Comments on standards should be submitted to the following address:

Secretary, IEEE-SA Standards Board 445 Hoes Lane Piscataway, NJ 08854 USA

Laws and regulations

Users of IEEE Standards documents should consult all applicable laws and regulations. Compliance with the provisions of any IEEE Standards document does not imply compliance to any applicable regulatory requirements. Implementers of the standard are responsible for observing or referring to the applicable regulatory requirements. IEEE does not, by the publication of its standards, intend to urge action that is not in compliance with applicable laws, and these documents may not be construed as doing so.

Copyrights

IEEE draft and approved standards are copyrighted by IEEE under U.S. and international copyright laws. They are made available by IEEE and are adopted for a wide variety of both public and private uses. These include both use, by reference, in laws and regulations, and use in private self-regulation, standardization, and the promotion of engineering practices and methods. By making these documents available for use and adoption by public authorities and private users, IEEE does not waive any rights in copyright to the documents.

Photocopies

Subject to payment of the appropriate fee, IEEE will grant users a limited, non-exclusive license to photocopy portions of any individual standard for company or organizational internal use or individual, non-commercial use only. To arrange for payment of licensing fees, please contact Copyright Clearance Center, Customer Service, 222 Rosewood Drive, Danvers, MA 01923 USA; +1 978 750 8400. Permission to photocopy portions of any individual standard for educational classroom use can also be obtained through the Copyright Clearance Center.

Updating of IEEE Standards documents

Users of IEEE Standards documents should be aware that these documents may be superseded at any time by the issuance of new editions or may be amended from time to time through the issuance of amendments, corrigenda, or errata. An official IEEE document at any point in time consists of the current edition of the document together with any amendments, corrigenda, or errata then in effect.

Every IEEE standard is subjected to review at least every ten years. When a document is more than ten years old and has not undergone a revision process, it is reasonable to conclude that its contents, although still of some value, do not wholly reflect the present state of the art. Users are cautioned to check to determine that they have the latest edition of any IEEE standard.

In order to determine whether a given document is the current edition and whether it has been amended through the issuance of amendments, corrigenda, or errata, visit the IEEE-SA Website at http://ieeexplore.ieee.org/xpl/standards.jsp or contact IEEE at the address listed previously. For more information about the IEEE-SA or IEEE's standards development process, visit the IEEE-SA Website at http://standards.jeee.org.

Errata

Errata, if any, for all IEEE standards can be accessed on the IEEE-SA Website at the following URL: http://standards.ieee.org/findstds/errata/index.html. Users are encouraged to check this URL for errata periodically.

Patents

Attention is called to the possibility that implementation of this standard may require use of subject matter covered by patent rights. By publication of this standard, no position is taken by the IEEE with respect to the existence or validity of any patent rights in connection therewith. If a patent holder or patent applicant has filed a statement of assurance via an Accepted Letter of Assurance, then the statement is listed on the IEEE-SA Website at http://standards.ieee.org/about/sasb/patcom/patents.html. Letters of Assurance may indicate whether the Submitter is willing or unwilling to grant licenses under patent rights without compensation or under reasonable rates, with reasonable terms and conditions that are demonstrably free of any unfair discrimination to applicants desiring to obtain such licenses.

Essential Patent Claims may exist for which a Letter of Assurance has not been received. The IEEE is not responsible for identifying Essential Patent Claims for which a license may be required, for conducting inquiries into the legal validity or scope of Patents Claims, or determining whether any licensing terms or conditions provided in connection with submission of a Letter of Assurance, if any, or in any licensing agreements are reasonable or non-discriminatory. Users of this standard are expressly advised that determination of the validity of any patent rights, and the risk of infringement of such rights, is entirely their own responsibility. Further information may be obtained from the IEEE Standards Association.

Participants

At the time this IEEE standard was completed, the Personal Health Devices Working Group had the following membership:

Daidi Zhong, Chair Michael J. Kirwan, Chair Melanie S. Yeung, Vice-Chair

Karsten Aalders Saeed A. Choudhary Charles R. Abbruscato Jinhan Chung Nabil Abujbara Malcolm Clarke Maher Abuzaid John A. Cogan Manfred Aigner John T. Collins Jorge Alberola Cory Condek Murtaza Ali Todd H. Cooper Rolf Ambuehl David Cornejo David Aparisi Douglas Coup Lawrence Arne Nigel Cox Diego B. Arquilo Hans Crommenacker Serafin Arroyo Tomio Crosley Muhammad Asim David Culp Merat Bagha Allen Curtis Doug Baird Ndifor Cyril Fru David Baker Jesús Daniel Trigo Anindya Bakshi Eyal Dassau Ananth Balasubramanian **David Davenport** Sunlee Bang Russell Davis M. Jonathan Barkley Ed Day Gilberto Barrón Sushil K. Deka David Bean Ciro de la Vega

John Bell Pedro de-las-Heras-Ouiros Rudy Belliardi Jim DelloStritto Kathryn M. Bennett Matthew d'Entremont Daniel Bernstein Lane Desborough George A. Bertos Kent Dicks Chris Biernacki Hyoungho Do Ola Björsne Xiaolian Duan Thomas Blackadar Brian Dubreuil Marc Blanchet Jakob Ehrensvard Thomas Bluethner Fredrik Einberg Douglas P. Bogia Roger M. Ellingson Xavier Boniface Michihiro Enokida Shannon Boucousis Javier Escayola Calvo Julius Broma Leonardo Estevez Lyle G. Bullock, Jr. Roger Feeley Bernard Burg Bosco T. Fernandes

Satva Calloji Russell Foster Carole C. Carey Eric Freudenthal Santiago Carot-Nemesio Matthias Frohner Randy W. Carroll Ken Fuchs Simon Carter Jing Gao Marcus Garbe Seungchul Chae Rahul Chauhan John Garguilo James Cheng Rick Geimer Igor Gejdos Peggy Chien

Julian Goldman
Raul Gonzalez Gomez
Chris Gough
Channa Gowda
Charles M. Gropper
Amit Gupta
Jeff Guttmacher
Rasmus Haahr
Christian Habermann
Michael Hagerty
Jerry Hahn
Robert Hall
Nathaniel Hamming

Nicolae Goga

Michael Hagerty
Jerry Hahn
Robert Hall
Nathaniel Hamming
Rickey L. Hampton
Sten Hanke
Jordan Hartmann
Kai Hassing
Marc Daniel Haunschild
Wolfgang Heck

Charles Henderson Jun-Ho Her Takashi Hibino Timothy L. Hirou Allen Hobbs Alex Holland Arto Holopainen Robert Hoy Frank Hsu Anne Huang Sen-Der Huang Zhiqiang Huang Ron Huby David Hughes Robert D. Hughes Jiyoung Huh Hugh Hunter Hitoshi Ikeda Yutaka Ikeda Philip O. Isaacson Atsushi Ito Michael Jaffe Praduman Jain Danny Jochelson Chris Johnson Phaneeth Junga

Akiyoshi Kabe

Tomio Kamioka

Steve Kahle

Kei Kariya

Andy Kaschl

Christoph Fischer

Morten Flintrup

Joseph W. Forler

Ferenc Gerbovics

Chris Burns

Anthony Butt

Jeremy Byford-Rew

Chia-Chin Chong

Junzo Kashihara Tetsu Nishimura Redmond Shouldice Kohichi Kashiwagi Jim Niswander Sternly K. Simon Ralph Kent Hiroaki Niwamoto Marjorie Skubic Laurie M. Kermes Thomas Norgall Robert Smith Ikuo Keshi Anand Noubade Ivan Soh Junhyung Kim Yoshiteru Nozoe Motoki Sone Minho Kim Emily Sopensky Abraham Ofek Min-Joon Kim **Brett Olive** Rajagopalan Srinivasan Taekon Kim Begonya Otal Andreas Staubert Charles Palmer Tetsuya Kimura Nicholas Steblay Alfred Kloos Bud Panjwani Beth Stephen Jeongmee Koh Carl Pantiskas Lars Steubesand Harry P. Pappas Jean-Marc Koller John (Ivo) Stivoric John Koon Mikey Paradis Raymond A. Strickland Patty Krantz Hanna Park Hermanni Suominen Alexander Kraus Jong-Tae Park Lee Surprenant Ravi Swami Ramesh Krishna Myungeun Park Soojun Park Ray Sweidan Geoffrey Kruse Phillip E. Pash Falko Kuester Jin Tan Rafael Lajara TongBi Pei Haruyuyki Tatsumi Pierre Landau Soren Petersen John W. Thomas Brad Tipler Jaechul Lee James Petisce JongMuk Lee Jonas Tirén Peter Piction Kyong Ho Lee Michael Pliskin James Tomcik Rami Lee Jeff Price Janet Traub Gary Tschautscher Sungkee Lee Harald Prinzhorn Woojae Lee John Ouinlan Masato Tsuchid Yonghee Lee Arif Rahman Ken Tubman Joe Lenart Tanzilur Rahman Yoshihiro Uchida Kathryn A. Lesh Steve Ray Sunil Unadkat Qiong Li Phillip Raymond Fabio Urbani Ying Li Tim Reilly Philipp Urbauer Patrick Lichter Barry Reinhold Laura Vanzago Jisoon Lim Brian Reinhold Alpo Värri Melvin I. Reynolds Joon-Ho Lim Dalimar Velez John G. Rhoads John Lin Naveen Verma Jiajia Liu Jeffrey S. Robbins Rudi Voon Moskowitz Robert Wei-Jung Lo Isobel Walker Timothy Robertson David Wang Charles Lowe David Rosales Jerry P. Wang Don Ludolph Christian Luszick Bill Saltzstein Yao Wang Benedikt Salzbrunn **Bob MacWilliams** Yi Wang Srikkanth Madhurbootheswaran Giovanna Sannino Steve Warren Fujio Watanabe Romain Marmot Jose A. Santos-Cadenas Toru Watsuji Sandra Martinez Stefan Sauermann Mike Weng Miguel Martínez de Espronceda John Sawyer Cámara Guillaume Schatz Kathleen Wible Paul Williamson Peter Mayhew Alois Schloegl Jim McCain Paul S. Schluter Jan Wittenber László Meleg Lars Schmitt Jia-Rong Wu Will Wykeham Alexander Mense Mark G. Schnell Richard A. Schrenker Ariton Xhafa Ethan Metsger Antonio Scorpiniti Junjie Yang Yu Miao Jinsei Miyazaki Kwang Seok Seo Ricky Yang Riccardo Serafin Melanie S. Yeung Erik Moll Sid Shaw Done-Sik Yoo Darr Moore Frank Shen Jason Zhang Piotr Murawski Liqun Shen Zhiqiang Zhang Soundharya Nagasubramanian Thomas Zhao Jae-Wook Nah Bozhi Shi Alex Neefus Min Shih Daidi Zhong

Mazen Shihabi

Trong-Nghia Nguyen-Dobinsky

Michael E. Nidd

Miha Zoubek

Szymon Zysko

The following members of the individual balloting committee voted on this standard. Balloters may have voted for approval, disapproval, or abstention.

John Ballingall Noriyuki Ikeuchi Henry Pinto Giberto Barrón Atsushi Ito Melvin I. Reynolds Lyle G. Bullock, Jr. Piotr Karocki Bartien Sayogo Keith Chow Patrick Keith-Hynes Lars Schmitt Joseph El Youssef Patrick Kinney Raymond A. Strickland Randall Groves Robert Kircher Walter Struppler Kai Hassing Michael J. Kirwan Jan Wittenber Werner Hoelzl Nick S. A. Nikjoo Oren Yuen

When the IEEE-SA Standards Board approved this standard on 16 February 2015, it had the following membership:

John Kulick, Chair Jon Walter Rosdahl, Vice-Chair Richard H. Hulett, Past Chair Konstantinos Karachalios, Secretary

Peter Balma Michael Janezic Ron Peterson Farooq Bari Jeffrey Katz Adrian Stephens Ted Burse Joseph L. Koepfinger* Peter Sutherland Clint Chaplain David Law Yatin Trivedi Stephen Dukes Hung Ling Phil Winston Oleg Logvinov Jean-Philippe Faure Don Wright Gary Hoffman T. W. Olsen Yu Yuan Glenn Parsons

*Member Emeritus

Also included are the following non-voting IEEE-SA Standards Board liaisons:

Richard DeBlasio, *DOE Representative* Michael Janezic, *NIST Representative*

Don Messina
IEEE-SA Content Production and Management

Kathryn Bennett
IEEE-SA Technical Program Operations

Introduction

This introduction is not part of IEEE Std 11073-10419-2015, Health informatics—Personal health device communication—Part 10419: Device Specialization—Insulin Pump.

ISO/IEEE 11073 standards enable communication between medical devices and external computer systems. This document uses the optimized framework created in ISO/IEEE 11073-20601:2010¹ and describes a specific, interoperable communication approach for insulin pumps. These standards align with and draw on the existing clinically focused standards to provide support for communication of data from clinical or personal health devices.

_

¹For information on references, see Clause 2.

Contents

1. Overview	1
1.1 Scope	1
1.2 Purpose	2
1.3 Context	
2. Normative references	2
3. Definitions, acronyms, and abbreviations	3
3.1 Definitions	
3.2 Acronyms and abbreviations	4
4. Introduction to ISO/IEEE 11073 personal health devices	5
4.1 General	5
4.2 Introduction to ISO/IEEE 11073-20601 modeling constructs	6
4.3 Compliance with other standards	6
5. Insulin pump device concepts and modalities	7
5.1 General	7
5.2 Device types	8
5.3 Collected data	8
5.4 Stored data	14
5.5 Scheduled data	15
6. Insulin pump domain information model	15
6.1 Overview	15
6.2 Class extensions	
6.3 Object instance diagram	
6.4 Types of configuration	
6.5 Profiles	
6.6 Medical device system object	
6.7 Numeric objects	
6.8 Real-time sample array objects	
6.9 Enumeration objects	
6.10 PM-store objects	
6.11 Schedule-store objects	
6.12 Scanner objects	
6.13 Class extension objects	
6.14 Insulin pump information model extensibility rules	
7. Insulin pump service model	56
7.1 General	
7.2 Object access services	
7.3 Object access event report services	
8. Insulin pump communication model	59
8.1 Overview	
8.2 Communications characteristics	
8.3 Association procedure	
8.4 Configuring procedure	
8.5 Operating procedure	
8.6 Time synchronization	

9. Test associations	64
9.1 Behavior with standard configuration	64
9.2 Behavior with extended configurations	
10. Conformance	64
10.1 Applicability	
10.2 Conformance specification	
10.3 Levels of conformance	
10.4 Implementation conformance statements	
Annex A (informative) Bibliography	71
Annex B (normative) Any additional ASN.1 definitions	72
B.1 Device status and insulin pump status bit mapping	
Annex C (normative) Allocation of identifiers	74
C.1 General	74
C.2 Definitions of terms and codes	74
C.3 Systematic derivations of terms and codes	76
Annex D (informative) Message sequence examples	85
Annex E (normative) Schedule-store class	87
E.1 Schedule-store class	87
E.2 Schedule-segment class	91
Annex F (normative) Schedule class ASN.1 definitions	
F.1 ACTION-method-related data types	
F.2 Data types for new object attributes and object services	
F.3 Data protocol definitions	100
Annex G (informative) The schedule-store concept	102
G.1 General	
G.2 Schedule-store object hierarchy	103
Annex H (informative) Scedule communication model	107
H.1 Operating procedure	107
Annex I (informative) Protocol data unit examples	
I.1 General	
I.2 Association information exchange	
I.3 Configuration information exchange	
I.4 GET MDS attributes service	
I.5 Data reporting	

Health informatics—Personal health device communication

Part 10419: Device Specialization—Insulin Pump

IMPORTANT NOTICE: IEEE Standards documents are not intended to ensure safety, health, or environmental protection, or ensure against interference with or from other devices or networks. Implementers of IEEE Standards documents are responsible for determining and complying with all appropriate safety, security, environmental, health, and interference protection practices and all applicable laws and regulations.

This IEEE document is made available for use subject to important notices and legal disclaimers. These notices and disclaimers appear in all publications containing this document and may be found under the heading "Important Notice" or "Important Notices and Disclaimers Concerning IEEE Documents." They can also be obtained on request from IEEE or viewed at http://standards.ieee.org/IPR/disclaimers.html.

1. Overview

1.1 Scope

The scope of this standard is to establish a normative definition of communication between personal telehealth insulin pump devices (agents) and managers (e.g., cell phones, personal computers, personal health appliances, set top boxes) in a manner that enables plug-and-play interoperability. It leverages work done in other ISO/IEEE 11073 standards including existing terminology, information profiles, application profile standards, and transport standards. It specifies the use of specific term codes, formats, and behaviors in telehealth environments restricting optionality in base frameworks in favor of interoperability. This standard defines a common core functionality of personal telehealth insulin pump devices.

In the context of personal health devices, an insulin pump is a medical device used for the administration of insulin in the treatment of diabetes mellitus, also known as continuous subcutaneous insulin infusion (CSII) therapy.

This standard provides the data modeling according to the ISO/IEEE 11073-20601 standard, and does not specify the measurement method.

Health informatics—Personal health device communication—Part 10419: Device Specialization—Insulin Pump

1.2 Purpose

This standard addresses the need for an openly defined, independent standard that support information exchange to and from personal health devices and compute engines (e.g., cell phones, personal computers, personal health appliances, and set top boxes). Interoperability is key to growing the potential market for these devices and enabling people to be better informed participants in the management of their heath.

1.3 Context

See IEEE Std 11073-20601TM-2014¹ for an overview of the environment within which this standard is written.

This standard defines the device specialization for the insulin pump, being a specific agent type, and provides a description of the device concepts, its capabilities, and its implementation according to this standard.

This standard is based on IEEE Std 11073-20601-2014, which in turn draw information from both ISO/IEEE 11073-10201:2004 [B7]² and ISO/IEEE 11073-20101:2004 [B8]. The medical device encoding rules (MDERs) used within this standard are fully described in IEEE Std 11073-20601-2014.

This standard reproduces relevant portions of the nomenclature found in ISO/IEEE 11073-10101:2004 [B6] and adds new nomenclature codes for the purposes of this standard. Among this standard and IEEE Std 11073-20601-2014, all required nomenclature codes for implementation are documented.

NOTE 1—IEEE Std 11073-20601-2014 is a revision of ISO/IEEE 11073-20601:2010 (and its amendment IEEE Std 11073-20601a). It contains new material and corrections and does not copy the content of ISO/IEEE 11073-20601:2010. Throughout this standard, a reference to IEEE Std 11073-20601-2014 refers to the document that is obtained after applying this new material and corrections to ISO/IEEE 11073-20601:2010.³

NOTE 2—In this standard, ISO/IEEE 11073-104zz is used to refer to the collection of device specialization standards that utilize IEEE Std 11073-20601-2014, where zz can be any number from 01 to 99, inclusive.

2. Normative references

The following referenced documents are indispensable for the application of this document (i.e., they must be understood and used, so that each referenced document is cited in text and its relationship to this document is explained). For dated references, only the edition cited applies. For undated references, the latest edition of the referenced document (including any amendments or corrigenda) applies.

IEEE Std 11073-20601-2014, Health informatics—Personal health device communication—Application Profile—Optimized Exchange Protocol.^{4,5}

See Annex A for all informative material referenced by this standard.

² The numbers in brackets correspond to those of the bibliography in Annex A.

¹ Information on references can be found in Clause 2.

³ Notes in text, tables, and figures are given for information only and do not contain requirements needed to implement the standard.

⁴ The IEEE standards or products referred to in this clause are trademarks of the Institute of Electrical and Electronics Engineers, Inc.

⁵ IEEE publications are available from the Institute of Electrical and Electronics Engineers, 445 Hoes Lane, Piscataway, NJ 08854-4141, USA (http://standards.ieee.org/).

3. Definitions, acronyms, and abbreviations

3.1 Definitions

For the purposes of this document, the following terms and definitions apply. The *IEEE Standards Dictionary Online* should be consulted for terms not defined in this clause.⁶

agent: A node that collects and transmits personal health data to an associated manager.

artificial pancreas: A system combining diabetes devices to provide similar functionality as a pancreas. Examples include linking a continuous glucose monitor to an insulin pump to automatically reduce or increase insulin infusion based upon specified thresholds of measured interstitial glucose.

basal insulin: Insulin required to cover the basic insulin needs of the body.

basal rate: Rate of continuously delivered insulin to cover the basic insulin needs of the body.

basal rate profile: Profile defining the time-dependent course of the basal rate over a periodic time interval (e.g., 24 hours, beginning at midnight).

blood glucose: Glucose concentration in the blood.

bolus: A single dose of a medication typically administered by infusion. See also: bolus insulin.

bolus insulin: Insulin required to cover the intake of food or to correct high blood glucose levels.

class: In object-oriented modeling, it describes the attributes, methods, and events that objects instantiated from the class utilize.

compute engine: See: manager.

correction bolus: A bolus to correct high blood glucose levels. See also: bolus.

daily basal dose: Total amount of basal insulin delivered in 24 hours, beginning at midnight.

device: A term used to refer to a physical apparatus implementing either an agent or a manager role.

extended bolus: A bolus given over an extended period of time with a specified duration. See also: bolus.

fast bolus: A bolus given over a very short period that can be considered as having infinitesimal duration and specified by amount of infusion. *See also:* **bolus**.

glucagon: Naturally occurring hormone released by the pancreas when glucose levels are low.

glucose: Commonly referred to as "sugar," it is the major source of energy used by the body cells.

handle: An unsigned 16-bit number that is locally unique and identifies one of the object instances within an agent.

infusion set: Tubing system with cannula that interfaces the insulin reservoir to the body.

⁶IEEE Standards Dictionary Online subscription is available at: http://www.ieee.org/portal/innovate/products/standard/standards_dictionary.html.

Health informatics—Personal health device communication—Part 10419: Device Specialization—Insulin Pump

insulin: Naturally occurring hormone required for carbohydrate metabolism. Insulin regulates the way sugar, fat, and proteins are moved into the cells and the way they are stored or used for energy.

insulin to carbohydrate ratio (I:CHO): Amount of carbohydrate covered or disposed of by a single unit of insulin.

insulin sensitivity factor (ISF): An estimated amount of blood glucose reduced by a single unit of insulin. Also known as a *correction factor*.

international unit (IU): Unit of measurement used in pharmacology for the measurement of the amount of a substance.

manager: A node receiving data from one or more agent systems. Some examples of managers include a cellular phone, health appliance, set top box, or computer system.

meal bolus: A bolus to compensate for meal intake. See also: bolus.

multi-wave bolus: A bolus combining a fast bolus with an extended bolus; one part of the bolus amount is delivered immediately while the other is delivered over a specific period of time.

obj-handle: See: handle.

object: In object-oriented modeling, a particular instantiation of a class. The instantiation realizes attributes, methods, and events from the class.

occlusion: Total obstruction of the infusion set that prevents administering insulin.

personal health device: A device used in personal health applications.

personal telehealth device: See: personal health device.

piston rod: Mechanical apparatus to propel the plunger of the insulin reservoir.

priming: Filling of the fluidic path with insulin (e.g., after replacement of the reservoir and/or infusion set).

reservoir: The insulin supply of an insulin pump.

temporary basal rate: Temporary adjustment of the basal rate.

total daily dose: The total amount of insulin delivered during one day (i.e., 24 hours, beginning at midnight). This amount contains basal insulin as well as bolus insulin.

3.2 Acronyms and abbreviations

AP artificial pancreas

APDU application protocol data unit
ASN.1 Abstract Syntax Notation One

CSII continuous subcutaneous insulin infusion, commonly known as insulin pump therapy

DIM domain information model

EUI-64 extended unique identifier (64 bits)

IEEE Std 11073-10419-2015

Health informatics—Personal health device communication—Part 10419: Device Specialization—Insulin Pump

HCP health care provider

I:CHO insulin to carbohydrate ratio

ICS implementation conformance statements

ICT intensified conventional therapy

ID identifier

IEEE Institute of Electrical and Electronics Engineers

ISF insulin sensitivity factor

ISO International Organization for Standardization

ITU-T Telecommunication Standardization Sector of International Telecommunication Union

IU international unit

MDC medical device communication

MDER medical device encoding rules

MDS medical device system

MOC managed object class

OID object identifier

PDU protocol data unit
PHD personal health device
PWD person with diabetes
RT-SA real-time sample array
VMO virtual medical object

virtual medical system

4. Introduction to ISO/IEEE 11073 personal health devices

4.1 General

VMS

This standard and the remainder of the series of ISO/IEEE 11073 personal health device (PHD) standards fit in the larger context of the ISO/IEEE 11073 series of standards. The full suite of standards enables agents to interconnect and interoperate with managers and with computerized health-care information systems. See IEEE Std 11073-20601-2014 for a description of the guiding principles for this series of ISO/IEEE 11073 Personal Health Device standards.

IEEE Std 11073-20601-2014 supports the modeling and implementation of an extensive set of personal health devices. This standard defines aspects of the insulin pump device. It describes all aspects necessary to implement the application layer services and data exchange protocol between an ISO/IEEE 11073 PHD insulin pump agent and a manager. This standard defines a subset of the objects and functionality contained in IEEE Std 11073-20601-2014 and extends and adds definitions where appropriate. All new definitions are given in Annex B in Abstract Syntax Notation One (ASN.1 [B9]). Nomenclature codes referenced in this standard, which are not defined in IEEE Std 11073-20601-2014, are normatively defined in Annex C.

Health informatics—Personal health device communication—Part 10419: Device Specialization—Insulin Pump

4.2 Introduction to ISO/IEEE 11073-20601 modeling constructs

4.2.1 General

The ISO/IEEE 11073 series of standards, and in particular IEEE Std 11073-20601-2014, is based on an object-oriented systems management paradigm. The overall system model is divided into three principal components: the domain information model (DIM), the service model, and the communication model. See IEEE Std 11073-20601-2014 for a detailed description of the modeling constructs.

4.2.2 Domain information model

The DIM is a hierarchical model that describes an agent as a set of objects. These objects and their attributes represent the elements that control behavior and report on the status of the agent and data that an agent can communicate to a manager. Communication between the agent and the manager is defined by the application protocol in the IEEE Std 11073-20601-2014.

4.2.3 Service model

The service model defines the conceptual mechanisms for the data exchange services. Such services are mapped to messages that are exchanged between the agent and the manager. Protocol messages within the ISO/IEEE 11073 series of standards are defined in ASN.1. The messages defined in IEEE Std 11073-20601-2014 can coexist with messages defined in other standard application profiles defined in the ISO/IEEE 11073 series of standards.

4.2.4 Communication model

In general, the communication model supports the topology of one or more agents communicating over logical point-to-point connections to a single manager. For each logical point-to-point connection, the dynamic system behavior is defined by a connection state machine as specified in IEEE Std 11073-20601-2014.

4.2.5 Implementing the models

An agent implementing this standard shall implement all mandatory elements of the information, service, and communication models as well as all conditional elements where the condition is met. The agent should implement the recommended elements, and it may implement any combination of the optional elements. A manager implementing this standard shall utilize at least one of the mandatory, conditional, recommended, or optional elements. In this context, *utilize* means to use the element as part of the primary function of the manager device. For example, a manager whose primary function is to display data would need to display a piece of data in the element in order to utilize it.

4.3 Compliance with other standards

Devices that comply with this standard may also be required to comply with other domain-and device-specific standards that supersede the requirements of this standard with respect to issues including safety, reliability, and risk management. A user of this standard is expected to be familiar with all other such standards that apply and to comply with any higher specifications thus imposed. Typically, medical devices should comply with the IEC 60601-1:2005 [B1] base standards with respect to electrical and mechanical

Health informatics—Personal health device communication—Part 10419: Device Specialization—Insulin Pump

safety and any device-specific standard as might be defined in the IEC 60601-2 [B2] series of standards. Software aspects may apply through standards such as IEC 62304:2006/EN 62304:2006 [B3].

Devices that comply with this standard implement higher layers of network software and utilize lower layers as appropriate to the application. The requirements on performance of such applications and conformance are defined elsewhere and are outside the scope of this standard. Moreover, the use of any medical equipment is subject to risk assessment and risk management appropriate to the application. Some relevant examples are ISO 14971:2007 [B5] and IEC 80001-1:2010 [B4]. The requirements of such risk assessment and risk management and conformance are outside the scope of this standard.

5. Insulin pump device concepts and modalities

5.1 General

This clause presents the general concepts of insulin pump devices. In the context of personal health devices in this family of standards, an insulin pump is a device that administers insulin at a defined amount to the human body. Insulin pump devices are primarily used in the continuous subcutaneous insulin infusion (CSII) therapy of type 1 diabetes mellitus. This type of diabetes mellitus is characterized by loss of the insulin-producing beta cells of the islets of Langerhans in the pancreas. Insulin pumps typically inject insulin into the subcutaneous layer of fat tissue under the skin through an infusion set. Preferred sites for the cannula are the abdomen, lumbar region, thighs, buttocks and the upper arms.

Insulin is a hormone that affects the metabolism and other body functions. It regulates, together with glucagon, the glucose level in the blood. Failure to produce sufficient insulin to maintain the blood glucose in the normal range leads to conditions of temporarily or persistently high blood sugar (hyperglycemia).

The objective of CSII therapy, also known as insulin pump therapy, is to mimic the insulin secretion of the healthy body. Insulin pump therapy therefore differentiates between basal insulin and bolus insulin. Basal insulin, sometimes also referred to as background insulin, regulates the blood glucose level without taking into account any food intake. In the context of the CSII therapy, basal insulin is continuously administered at a given rate depending on the fluctuating need throughout the day. By contrast, bolus insulin compensates for an increased need for insulin induced by food intake. This is also known as meal bolus and the insulin dosage can be determined with the knowledge of the user's insulin to carbohydrate ratio (the amount of carbohydrate covered by a unit of insulin). In addition, bolus insulin may also be administered to decrease a temporarily high blood glucose level. This is often referred to as correction bolus. Meal or correction boluses can be delivered as a standard fast bolus, extended bolus, or as a multi-wave bolus. Fast bolus is an infusion of insulin over a very short period similar to a pen injection that can be considered as having infinitesimal duration. An extended bolus is an infusion of insulin over a longer period extending the insulin delivery. A multi-wave bolus combines a fast bolus with an extended bolus; one part of the bolus amount is delivered immediately while the other is delivered over a specific period of time. This type of bolus is suitable for meals that contain both fast and slowly digestible carbohydrates or long meals with several courses.

Although the distinction between basal and bolus insulin has a high relevance for the CSII therapy, the insulin pump typically delivers a single insulin type for both delivery methods. Unlike the Intensified Conventional Therapy (ICT, injection therapy, pen therapy) where typically a combination of injections with fast-acting and slow-acting insulin is utilized, the insulin pump therapy generally uses rapid-or fast-acting insulin. Another term for these types of insulin is insulin analog.

Many insulin pump devices allow the pump user to program basal rate profiles that define the fluctuating delivery of basal insulin of the pump in the course of the day. These profile settings are often referred to as circadian profiles. Well-adjusted basal rate profiles generally cover the insulin pump user's need for basal

Health informatics—Personal health device communication—Part 10419: Device Specialization—Insulin Pump

insulin and there is normally no necessity to adjust these profiles on a daily basis. However, there are situations in which the need for basal insulin deviates at such a rate that a temporary adjustment of the basal delivery is indicated. Phases of stress, intensified physical activity, illness, and menses are counted amongst the reasons for increased or reduced need for basal insulin. The temporary basal rate is a means to make short-term adjustments to the basal insulin delivery. Two different concepts for the temporary basal rate are currently well accepted in the area of insulin pump therapy. The first concept applies a scaling factor to the programmed basal rate for a certain period of time. The second concept consists in overriding the programmed basal rate with a specific value for a certain time period.

The unit of measurement to define the dose of insulin is international unit (IU). One IU of insulin is the biological equivalent of 1/22 mg of pure crystalline insulin. The concentration of insulin in an insulin solution is defined by international units per milliliter (i.e., IU/mL). In pharmacology, it is also very common to indicate the concentration of insulin solutions by the letter U followed by the number of international units per milliliter. For example, a concentration of 100 IU/mL is denominated as U100.

5.2 Device types

Insulin pump devices are generally designed to be portable and permanently connected to the body.

The structural shape of insulin pumps may vary, but insulin pump devices typically include the following components: the pump unit, a disposable insulin reservoir, and a disposable infusion set.

The pump unit typically comprises a computing unit and a dosing mechanism, both powered by a battery. Additionally, pump units may contain several user interface elements such as displays, buttons, or apparatuses to generate visual, haptic, or acoustic signals.

The disposable insulin reservoir typically consists of a glass or plastic ampulla with a volumetric capacity of a few milliliters of insulin solution. Depending on the insulin need of the pump user, this capacity is sufficient for provisioning insulin for several days.

The disposable infusion set interfaces the insulin reservoir to the human body. It generally comprises an adapter that connects to the insulin reservoir, a tubing system and a cannula that is inserted into the subcutaneous fat tissue. It is generally recommended to replace the infusion set after 2 to 3 days of usage for sterility reasons and to help prevent irritation at the insertion site.

5.3 Collected data

5.3.1 General

The insulin pump is a portable device and therefore may not be connected to a manager while collecting data. The two main use cases for an insulin pump agent to connect to a manager and send its data are the following:

— The insulin pump user visits a health care professional (HCP) to examine the adequacy of the therapy. The HCP normally compares the historic data from the insulin pump with the corresponding data from a blood glucose meter or continuous glucose monitor device to derive necessary adjustments to the insulin pump therapy. As the interval between such visits may constitute several months, insulin pumps are typically capable of storing data for such time periods.

Health informatics—Personal health device communication—Part 10419: Device Specialization—Insulin Pump

— The insulin pump user connects the insulin pump agent to a manager at home to examine the adequacy of the therapy and to apply adjustments when indicated. This tends to happen on a more frequent basis (e.g., once per week).

In addition to the two use cases above, an insulin pump agent may also be continuously connected to a manager to report collected data (e.g., artificial pancreas).

5.3.2 Current bolus setting

The current bolus setting represents the amount of insulin requested to be delivered as a bolus and is reported before the pump delivers the bolus. A bolus is typically initiated by the insulin pump user to compensate for an increased need for insulin induced by intake of food containing a relatively high quantity of carbohydrates (e.g., meal bolus). A bolus may also be indicated to decrease a temporarily high blood glucose level (e.g., correction bolus).

A current bolus setting can be delivered in a theoretically infinitesimal period of time (e.g., fast bolus) or can be delivered over a defined period of time (e.g., extended bolus) or as a combination of a fast bolus with an extended bolus (e.g., multi-wave bolus).

A current bolus setting can be set to be delivered at a later time. In this instance, the current bolus setting shall be used with the pending bolus delay together to specify the amount of delay between programming a current bolus setting and the bolus delivered.

A current bolus setting can be cancelled. In this instance, the current bolus setting shall be used with the insulin pump status to specify the bolus canceled.

A bolus amount can be set manually by the user, or recommended by the insulin pump's bolus calculator, or commanded from an artificial pancreas controller. If an insulin pump device has a bolus calculator, the user may accept or change the recommended amount of insulin.

The unit for the bolus setting is international unit (IU) or international unit per hour (IU/hr).

NOTE—This object represents a device setting and not a measurement value.

5.3.3 Pending bolus delay

The pending bolus delay specifies the delay between programming a current bolus setting and the bolus delivered. This setting is helpful for patients with delayed digestion (gastroparesis).

The unit for the pending bolus delay is in minutes (min).

NOTE—This object represents a device setting and not a measurement value.

5.3.4 Bolus delivered

The bolus delivered measurement represents the amount of insulin effectively delivered by the insulin pump as a bolus and is reported after the pump has delivered the bolus.

Health informatics—Personal health device communication—Part 10419: Device Specialization—Insulin Pump

The measurement represents the delivered amount of insulin for a bolus and the modality can be described as a fast bolus, extended bolus, or multi-wave bolus, and the reason for delivering the bolus as meal bolus or correction bolus (Figure 1). The unit for the bolus delivered is international unit (IU).

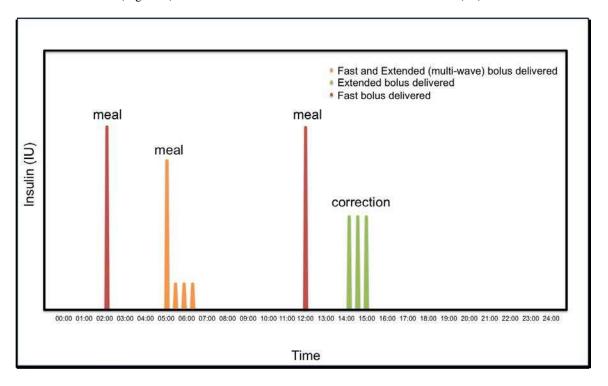


Figure 1—Example of bolus delivered

5.3.5 Basal rate schedule setting

The purpose of a basal rate schedule setting is to define the amount of basal insulin programmed in an insulin pump to be delivered for a specified schedule according to the individually fluctuating physiological needs. One or more basal rate schedule settings form a 24 hour basal profile.

The unit of a basal rate schedule setting is international unit per hour (IU/h).

NOTE—This object represents a device setting and not a measurement value.

5.3.6 Basal profile setting

Insulin pump devices typically offer the option of programming multiple basal profile settings in order to easily meet changing insulin needs (e.g., work day versus weekend).

In the context of continuous subcutaneous insulin infusion (CSII), a basal rate designates the amount of basal insulin programmed to deliver over a specified period of time.

Health informatics—Personal health device communication—Part 10419: Device Specialization—Insulin Pump

The basal profile setting is a schedule-segment that consists of a series of basal rate schedule settings programmed on the insulin pump device. The basal profile setting has a schedule that typically begins at 12 a.m. and ends at 12 a.m. the following day accumulating a total duration of 24 hours (Figure 2).

Insulin pump devices typically offer the option of programming multiple basal profile settings in order to easily meet changing insulin needs (e.g., work day versus weekend). The basal profile setting will also have a unique identifier specified by the manufacturer to identify a basal profile setting.

If an insulin pump device supports programming basal profile settings, an active basal profile setting can be selected.

NOTE—This object represents a device setting and not a measurement value.

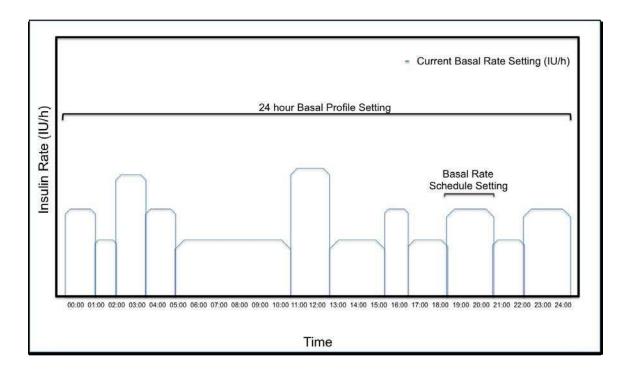


Figure 2—Example of basal profile setting consisting of basal rate schedule settings

5.3.7 Current basal rate setting

The current basal rate setting represents the most recent basal rate to be delivered by the pump and is reported before the pump delivers the basal rate. Current basal rate setting may be initiated by a programmed basal profile or a temporary basal rate adjustment made by the user.

A basal profile is a scheduled sequence of basal rate schedule settings programmed on the insulin device. The active basal profile and current time may define the current basal rate setting.

Alternatively, a user may initiate a temporary basal rate adjustment as a result of a temporarily changed need. This adjustment remains active for a defined period of time and afterwards the current basal rate

IEEE Std 11073-10419-2015

Health informatics—Personal health device communication—Part 10419: Device Specialization—Insulin Pump

setting will revert to the programmed basal profile. The temporary basal rate setting can be initiated by a *relative* or an *absolute* adjustment.

A *relative* temporary basal rate setting applies a dimensionless scaling factor to the current basal rate for the active basal profile. Therefore if the relative temporary basal rate is active, the rate delivery equals the value of the reference factor multiplied by the basal rates given by the active basal profile for a given period of time. A relative temporary basal rate shall be the calculated rate after applying the scaling factor.

An *absolute* temporary basal rate setting overrides the basal rates given by the active basal rate profile for a given period of time. Therefore, if the absolute temporary basal rate is active, the rate delivery equals the value of the absolute temporary basal rate for a given period of time.

This standard assumes that the effective delivery of basal insulin can be represented or approximated by discrete values of constant delivery rate over time. This standard assumes that the temporary basal rate setting overrides the basal rates given by the active basal rate profile. Therefore, if the temporary basal rate (either relative or absolute) is active, the rate of the actual basal insulin delivery equals the value of the temporary basal rate. The unit of the delivery rate is international unit per hour (IU/h).

NOTE—This object represents a device setting and not a measurement value.

5.3.8 Basal delivered

The basal delivered measurement represents the amount of insulin effectively delivered as a basal by an insulin pump and is reported after the pump has delivered the basal. An insulin pump device commonly utilizes semi-continuous infusion and incrementally injects discrete amounts (often called a *micro-bolus*) of insulin to satisfy the current basal rate setting (Figure 3). The insulin pump is responsible for determining the frequency of the basal amount to be reported (e.g., insulin pump may report the amount of insulin delivered at the end of a basal rate or intermittently as micro-boluses).

The measurement represents the delivered amount of basal insulin and can be described as a programmed basal rate, basal rate set by a controller, a temporary absolute basal, temporary relative basal, or undetermined basal rate. The unit for the basal delivered measurement is international unit (IU).

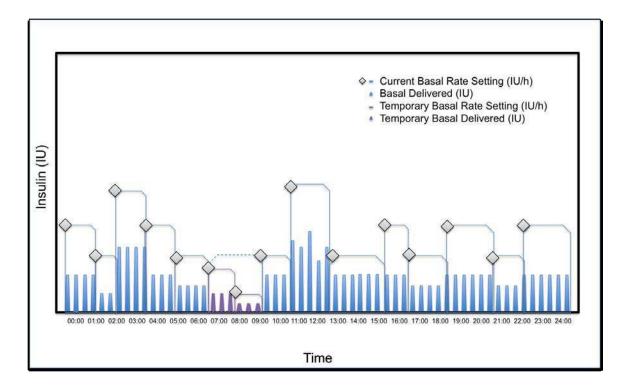


Figure 3—Example of basal rate setting and basal delivered

5.3.9 Insulin reservoir remaining

The insulin reservoir remaining measurement represents the amount of insulin remaining in the reservoir. The reservoir holds insulin inside the insulin pump and can have varying capacity. As the pump delivers insulin, the insulin amount in the reservoir decreases. When the reservoir is empty, users typically refill and insert a new reservoir into the insulin pump. The insulin pump is responsible for determining the frequency of the insulin reservoir remaining amount to be reported (e.g., the insulin reservoir remaining may be reported at the same time as the bolus-delivered or basal-delivered measurements).

The measurement represents the amount of insulin remaining in the reservoir. The unit for the insulin reservoir remaining is international unit (IU).

5.3.10 Insulin concentration

The insulin concentration measurement represents the concentration of insulin in the reservoir. Normally it is specified as the amount of insulin in international units per milliliter of insulin solution. For example, U100 designates an insulin concentration of 100 international units per milliliter of insulin solution.

The measurement represents the insulin concentration. The unit for the insulin concentration is international unit per milliliter (IU/mL).

Health informatics—Personal health device communication—Part 10419: Device Specialization—Insulin Pump

5.3.11 Insulin to carbohydrate ratio profile settings

This represents the patient's insulin to carbohydrate ratio (I:CHO), which is the amount of carbohydrate covered or disposed of by a 1 IU of insulin.

Insulin pump devices may offer the option of storing multiple insulin to carbohydrate ratios as a patient may have varying ratios specific to a time of day or type of meal. The unit for the insulin to carbohydrate ratio profile setting is grams (g).

NOTE—This object represents a device setting and not a measurement value.

5.3.12 Insulin sensitivity factor profile settings

This represents the patient's insulin sensitivity factor (ISF) or sometimes referred to as a *correction factor*, which is an estimated amount of blood glucose reduced by a single unit of insulin.

Insulin pump devices may offer the option of storing multiple insulin sensitivity factors as a patient may have varying factors specific to the time of day. The unit for the insulin sensitivity factor profile setting is millimole/liter (mmol/L) or milligrams/deciliter (mg/dL).

NOTE—This object represents a device setting and not a measurement value.

5.3.13 Operational status

The *operational status* represents the operation and therapy conditions of the insulin pump device. Operational conditions are Off, Stand-By, Preparing, Priming, Waiting, Ready, and Undetermined. Therapy conditions are Pause, Stop, and Run.

5.3.14 PHD DM status

The *PHD DM status* allows generic notification handling for personal health devices. The status of the device is represented in a number of bit flags covering warning, error, service, and undetermined statuses. Changes in status are sent as events notifying the manager.

This represents a generic notification mechanism for a device to report warnings, errors, and other abnormal statuses to a manager that may provide visibility of these status change events to the user.

5.3.15 Insulin pump status

The *insulin pump status* represents the specific notifications given by the insulin pump for handling, warnings, errors, reminders, or confirmations to user-initiated events meant to provide visibility of events to the user and manufacture.

5.4 Stored data

As stated in 5.3.1, an insulin pump may be used over several months of operation without being connected to a manager to send its data. Once an insulin pump is connected to a manager, the manager is able to select

Health informatics—Personal health device communication—Part 10419: Device Specialization—Insulin Pump

which of the agent's stored measurements or observations to retrieve. Depending on the agent's capabilities to organize its data into clusters of chronologically contiguous data, the manager may also select the time ranges of the stored data to retrieve. The agent then transmits the manager's selection in one or several blocks of messages for processing by a manager or other processing apparatus. The manager may also be able to choose a set of data clusters for deletion.

5.5 Scheduled data

An insulin pump may have the ability to program scheduled events such as basal profile settings, insulin to carbohydrate ratio profile settings, and insulin sensitivity factor profile settings. These scheduled data provide information about the actions or events and time points that the actions or events occur. Once an insulin pump is connected to a manager, the manager is able to select which of the agent's scheduled data to retrieve. The agent then transmits the manager's selection in one or several blocks of messages for processing by a manager or other processing apparatus. The manager may also be able to choose a set of data clusters for deletion.

6. Insulin pump domain information model

6.1 Overview

This clause describes the domain information model of the insulin pump.

6.2 Class extensions

In this standard, no class extensions are defined with respect to IEEE Std 11073-20601-2014.

6.3 Object instance diagram

The metric object instance diagram of the insulin pump domain information model, defined for the purposes of this standard, is shown in Figure 4.

The objects of the DIM, as shown in Figure 4, are described in 6.6 through 6.11. This includes the medical device system (MDS) object (see 6.6), the numeric objects (see 6.7), the real-time sample array (RT-SA) objects (see 6.8), the enumeration objects (see 6.9), the persistent metric—store (PM-store) objects (see 6.10), and the scanner objects (see 6.11). See 6.14 for rules for extending the insulin pump information model beyond elements as described in this standard. Each clause that describes an object of the insulin pump contains the following information:

- The nomenclature code used to identify the class of the object. One example where this code is used is the configuration event, where the object class is reported for each object. This allows the manager to determine whether the class of the object being specified is a numeric, real-time sample array, enumeration, scanner, or PM-store class.
- The attributes of the object. Each object has attributes that represent and convey information on the physical device and its data sources. Each object has a handle attribute that identifies the object instance within an agent. Attribute values are accessed and modified using methods such as GET and SET. Attribute types are defined using ASN.1. The ASN.1 definitions for new attribute types

IEEE Std 11073-10419-2015

Health informatics—Personal health device communication—Part 10419: Device Specialization—Insulin Pump

specific to this standard are in Annex B, and the ASN.1 definitions for existing attribute types referenced in this standard are in IEEE Std 11073-20601-2014.

- The methods available on the object.
- The potential events generated by the object. The data are sent to the manager using events.
- The available services such as getting or setting attributes.

The attributes for each class are defined in tables that specify the name of the attribute, its value, and its qualifier. The qualifiers mean:

- M: attribute is mandatory,
- C: attribute is conditional and depends on the condition stated in the Remark or Value column (if IEEE Std 11073-20601-2014 is referenced, then it contains the conditions),
- R: attribute is recommended,
- NR: attribute is not recommended, and
- O: attribute is optional.

Mandatory attributes shall be implemented by an agent. Conditional attributes shall be implemented if the condition applies and may be implemented otherwise. Recommended attributes should be implemented by the agent. Not recommended attributes should not be implemented by the agent. Optional attributes may be implemented by the agent. For attributes with qualifiers set to R or NR, underlying requirements stated in the Remark and Value column in IEEE Std 11073-20601-2014 shall be followed.

The attributes can be either static, meaning that they shall remain unchanged after the configuration is agreed upon, or dynamic, meaning that the attribute may change at some point after configuration.

Health informatics—Personal health device communication—Part 10419: Device Specialization—Insulin Pump

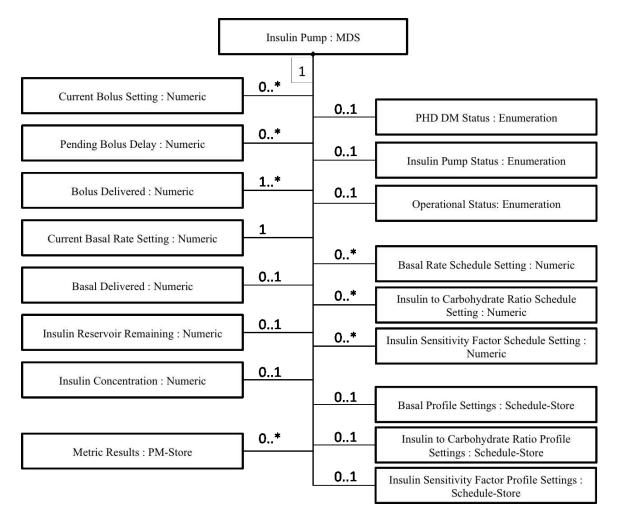


Figure 4—Insulin pump domain information model

6.4 Types of configuration

6.4.1 General

As specified in IEEE Std 11073-20601-2014, there are two styles of configuration available. Subclauses 6.4.2 and 6.4.3 briefly introduce standard and extended configurations.

6.4.2 Standard configuration

Standard configurations are defined in the ISO/IEEE 11073-104zz specializations (such as this standard) and are assigned a well-known identifier (Dev-Configuration-Id). The usage of a standard configuration is negotiated at association time between the agent and the manager. If the manager recognizes and selects to operate using the configuration, then the agent can send measurements immediately. If the manager does not recognize the configuration, the agent provides the configuration prior to transmitting measurement information.

One standard configuration is defined in this standard. The standard configuration 1900 (0x076C) contains the objects bolus delivered (see 6.7.4), and current basal rate setting (see 6.7.5).

17 Copyright © 2015 IEEE. All rights reserved.

Health informatics—Personal health device communication—Part 10419: Device Specialization—Insulin Pump

6.4.3 Extended configuration

In extended configurations, the agent's configuration is not predefined in a standard. The agent determines the objects, attributes, and values that to be used in a configuration and assigns a configuration identifier. When the agent associates with a manager, an acceptable configuration is negotiated. Typically, the manager does not recognize the agent's configuration on the first connection, so the manager responds that the agent needs to send its configuration information as a configuration event report. If, however, the manager recognizes the configuration, either because it was preloaded in some way or the agent had previously associated with the manager, then the manager responds that the configuration is known and no further configuration information needs to be sent.

6.5 Profiles

6.5.1 General

A profile further constrains the objects, services, and communication model of a specialization. By profiling the device specialization, the standard provides more guidance on the specific mandatory objects that shall be implemented, the objects that are optional, and the objects that are not required. This standard does not define profiles for the insulin pump device.

6.6 Medical device system object

6.6.1 MDS object attributes

Table 1 summarizes the attributes of the insulin pump MDS object. The nomenclature code to identify the MDS class is MDC MOC VMS MDS SIMP.

Table 1—MDS object attributes

Walue

Attribute name	Value	Qualifier
Handle	0	M
System-Type	Attribute not present.	С
	See IEEE Std 11073-20601-2014.	
System-Type-Spec-List	{MDC_DEV_SPEC_PROFILE_INSULIN_PUMP, 1}.	M
System-Model	{"Manufacturer", "Model"}	M
System-Id	Extended unique identifier (64-bits) (EUI-64)	M
Dev-Configuration-Id	Standard config: 0x076C (1900)	M
	Extended configs: 0x4000-0x7FFF	
Attribute-Value-Map	See IEEE Std 11073-20601-2014.	C
Production-Specification	See IEEE Std 11073-20601-2014.	O
Mds-Time-Info	See IEEE Std 11073-20601-2014.	C
Date-and-Time	See IEEE Std 11073-20601-2014.	C
Base-Offset-Time	See IEEE Std 11073-20601-2014.	M
Relative-Time	See IEEE Std 11073-20601-2014.	C
HiRes-Relative-Time	See IEEE Std 11073-20601-2014.	C
Date-and-Time-Adjustment	See IEEE Std 11073-20601-2014.	C
Power-Status	onBattery or onMains	R
Battery-Level	See IEEE Std 11073-20601-2014.	R
Remaining-Battery-Time	See IEEE Std 11073-20601-2014.	R
Reg-Cert-Data-List	See IEEE Std 11073-20601-2014.	0
Confirm-Timeout	See IEEE Std 11073-20601-2014.	0

NOTE—See IEEE Std 11073-20601-2014 for information on whether an attribute is static or dynamic.

Health informatics—Personal health device communication—Part 10419: Device Specialization—Insulin Pump

In the response to a Get MDS Object command, only implemented attributes and their corresponding values are returned.

See IEEE Std 11073-20601-2014 for descriptive explanations of the individual attributes as well as for information on attribute ID and attribute type.

The Dev-Configuration-Id attribute holds a locally unique 16-bit identifier that identifies the device configuration. For an insulin pump agent with extended configuration, this identifier is chosen in the range of extended-config-start to extended-config-end (see IEEE Std 11073-20601-2014) as shown in Table 1.

The agent sends the Dev-Configuration-Id during the Associating state (see 8.3) to identify its configuration for the duration of the association. If the manager already holds the configuration information relating to the Dev-Configuration-Id, it recognizes the Dev-Configuration-Id and the Configuring state (see 8.4) is skipped, and the agent and manager then enter the Operating state. If the manager does not recognize the Dev-Configuration-Id, the agent and manager enter the Configuring state.

If an agent implements multiple ISO/IEEE 11073-104zz specializations, System-Type-Spec-List is a list of type/version pairs, each referencing the respective device specialization and version of that specialization.

6.6.2 MDS object methods

Table 2 defines the methods (actions) of the MDS object. These methods are invoked using the Action service. In Table 2, the Subservice type name column defines the name of the method; the Mode column defines whether the method is invoked as an unconfirmed action (i.e., roiv-cmip-action from IEEE Std 11073-20601-2014) or a confirmed action (i.e., roiv-cmip-confirmed-action); the Subservice type (action-type) column defines the nomenclature code to use in the action-type field of an action request and response (see IEEE Std 11073-20601-2014); the Parameters (action-info-args) column defines the associated ASN.1 data structure (see IEEE Std 11073-20601-2014 for ASN.1 definitions) to use in the action message for the action-info-args field of the request; and the Results (action-info-args) column defines the structure to use in the action-info-args of the response.

Table	2-MDS	obiect	methods
-------	-------	--------	---------

Service	Subservice type name	Mode	Subservice type (action-type)	Parameters (action-info-args)	Results (action- info-args)
ACTION	Set-Time	Confirmed	MDC_ACT_SET_TIME	SetTimeInvoke	
ACTION	Set-Base- Offset-Time	Confirmed	MDC_ACT_SET_BO_TIME	SetBOTimeInvoke	_

— Set-Time

This method allows the manager to set a real-time clock in the agent with the absolute time. The agent indicates whether the Set-Time command is valid using the mds-time-capab-set-clock bit in the Mds-Time-Info attribute (see IEEE Std 11073-20601-2014).

If the agent supports the Absolute-Time-Stamp attribute, this method shall be implemented.

— Set-Base-Offset-Time

This method allows the manager to set a real-time clock in the agent with the base time and offset. The agent indicates whether the Set-Base-Offset-Time command is valid using the mds-time-capab-set-clock bit in the Mds-Time-Info attribute (see IEEE Std 11073-20601-2014).

If the agent supports the Base-Offset-Time-Stamp attribute, this method shall be implemented.

Health informatics—Personal health device communication—Part 10419: Device Specialization—Insulin Pump

Agents following only this device specialization and no others shall send event reports using agent initiated measurement data transmission. Agents following this device specialization as well as others shall send event reports in the appropriate fashion. During the association procedure (see 8.3), data-req-mode-capab shall be set to the appropriate value for the event report style. As a result, the manager shall assume the insulin pump agent does not support any of the MDS-Data-Request features (see IEEE Std 11073-20601-2014 for additional information). Thus, implementation of the MDS-Data-Request method/action is not required in this standard and is not shown in Table 2.

6.6.3 MDS object events

Table 3 defines the events that can be sent by the insulin pump MDS object.

Service	Subservice type name	Mode	Subservice type (event-type)	Parameters (event-info)	Results (event-reply- info)
	MDS- Configuration- Event	Confirmed	MDC_NOTI_CONFIG	ConfigReport	ConfigReport Rsp
	MDS- Dynamic- Data-Update- Var	Confirmed	MDC_NOTI_SCAN_ REPORT_VAR	ScanReportInfoVar	_
EVENT REPORT	MDS- Dynamic- Data-Update- Fixed	Confirmed	MDC_NOTI_SCAN_ REPORT_FIXED	ScanReportInfoFixed	
	MDS- Dynamic- Data-Update- MP-Var	Confirmed	MDC_NOTI_SCAN_ REPORT_MP_VAR	ScanReportInfoMP Var	_
	MDS- Dynamic- Data-Update- MP-Fixed	Confirmed	MDC_NOTI_SCAN_ REPORT_MP_FIXED	ScanReportInfoMP Fixed	

Table 3—Insulin pump MDS object events

— MDS-Configuration-Event

This event is sent by the insulin pump agent during the configuring procedure if the manager does not already know the insulin pump agent's configuration from past associations or because the manager has not been implemented to recognize the configuration according to the insulin pump device specialization. The event provides static information about the supported measurement capabilities of the insulin pump agent.

— MDS-Dynamic-Data-Update-Var

This event provides dynamic measurement data from the insulin pump agent for the numeric and enumeration objects. These data are reported using a generic attribute list variable format. The event is sent as an unsolicited message by the agent (i.e., an agent-initiated measurement data transmission). See 8.5.3 for more information on unsolicited event reporting.

— MDS-Dynamic-Data-Update-Fixed

This event provides dynamic measurement data from the insulin pump agent for the numeric and enumeration objects. These data are reported in the fixed format defined by the Attribute-Value-Map attribute of the objects. The event is sent as an unsolicited message by the agent (i.e., an agent-initiated measurement data transmission). See 8.5.3 for more information on unsolicited event reporting.

Health informatics—Personal health device communication—Part 10419: Device Specialization—Insulin Pump

— MDS-Dynamic-Data-Update-MP-Var

This is the same as MDS-Dynamic-Data-Update-Var but allows inclusion of data from multiple people.

— MDS-Dynamic-Data-Update-MP-Fixed

This is the same as MDS-Dynamic-Data-Update-Fixed but allows inclusion of data from multiple people.

NOTE—IEEE Std 11073-20601-2014 requires that managers support all of the MDS Object Events listed above.

6.6.4 Other MDS services

6.6.4.1 GET service

An insulin pump agent shall support the GET service, which is provided by the MDS object to retrieve the values of all implemented MDS object attributes. The GET service can be invoked as soon as the insulin pump agent receives the Association Response and moves to the Associated state, including the Operating and Configuring substates.

The manager may request the MDS object attributes of the insulin pump agent; in which case, the manager shall send the "Remote Operation Invoke | Get" message (see roiv-cmip-get in IEEE Std 11073-20601-2014) with the reserved MDS handle value of 0. The insulin pump agent shall report its MDS object attributes to the manager using the "Remote Operation Response | Get" message (see rors-cmip-get in IEEE Std 11073-20601-2014). See Table 4 for a summary of the GET service including some message fields.

Table 4—Insulin pump MDS object GET service

Service	Subservice type name	Mode	Subservice type	Parameters	Results
GET	<na></na>	<implied< td=""><td><na></na></td><td>GetArgumentSimple</td><td>GetResultSimple</td></implied<>	<na></na>	GetArgumentSimple	GetResultSimple
		confirmed>		= (obj-handle $=$ 0), attribute-	= (obj-handle =
				id-list <optional></optional>	0), attribute-list

See 8.5.2 for details on the procedure for getting the MDS object attributes.

6.6.4.2 SET service

The insulin pump specialization does not require an implementation to support the MDS object SET service.

6.7 Numeric objects

6.7.1 General

The insulin pump DIM for metric objects (see Figure 4) contains numeric objects that represent aspects of insulin delivery.

Sometimes, the interpretation of one attribute value in an object depends on other attribute values in the same object. For example, Unit-Code and Unit-LabelString provide context for the observed values. Whenever a contextual attribute changes, the agent shall report these changes to the manager using an MDS object event (see 6.6.3) prior to reporting any of the dependent values.

Health informatics—Personal health device communication—Part 10419: Device Specialization—Insulin Pump

6.7.2 Current bolus setting

Table 5 summarizes the attributes of the current bolus setting numeric object. The nomenclature code to identify the numeric class is MDC_MOC_VMO_METRIC_NU. The current bolus setting numeric object may be present in the extended configuration.

Table 5—Current bolus setting numeric object attributes

Attribute name	Extended configuration			
Attribute name	Value	Qualifier		
Handle	See IEEE Std 11073-20601-2014.	M		
Туре	{MDC_PART_PHD_DM	M		
	MDC_INS_BOLUS_SETTING}			
Supplemental-Types	Please see following text in Table 6.	R		
Metric-Spec-Small	mss-avail-intermittent mss-avail-stored-data	M		
	mss-upd-aperiodic mss-acc-manager-initiated			
	mss-acc-agent-initiated mss-cat-manual mss-			
	cat-setting mss-cat-calculation			
Metric-Structure-Small	See IEEE Std 11073-20601-2014.	0		
Measurement-Status	See IEEE Std 11073-20601-2014.	0		
Metric-Id	Please see following text in Table 7.	R		
Metric-Id-List	See IEEE Std 11073-20601-2014.	O		
Metric-Id-Partition	See IEEE Std 11073-20601-2014.	NR		
Unit-Code	Please see following text.	M		
Attribute-Value-Map	See IEEE Std 11073-20601-2014.	С		
Source-Handle-Reference	See IEEE Std 11073-20601-2014.	NR		
Label-String	See IEEE Std 11073-20601-2014.	0		
Unit-Label String	See IEEE Std 11073-20601-2014.	0		
Absolute-Time-Stamp	See IEEE Std 11073-20601-2014.	С		
Base-Offset-Time-Stamp	See IEEE Std 11073-20601-2014.	R		
Relative-Time-Stamp	See IEEE Std 11073-20601-2014.	С		
HiRes-Time-Stamp	See IEEE Std 11073-20601-2014.	С		
Measure-Active-Period	Please see following text.	0		
Simple-Nu-Observed-Value	See IEEE Std 11073-20601-2014.	С		
Compound-Simple-Nu-Observed-Value	See IEEE Std 11073-20601-2014.	NR		
Basic-Nu-Observed-Value	See IEEE Std 11073-20601-2014.	R		
Compound-Basic-Nu-Observed-Value	See IEEE Std 11073-20601-2014.	С		
Nu-Observed-Value	See IEEE Std 11073-20601-2014.	С		
Compound-Nu-Observed-Value	See IEEE Std 11073-20601-2014.	С		
Accuracy	See IEEE Std 11073-20601-2014.	NR		

NOTE 1—See IEEE Std 11073-20601-2014 for information on whether an attribute is static or dynamic.

NOTE 2—See 6.3 for a description of the qualifiers.

The value reported in this object is the amount of insulin requested to be delivered as a bolus and is reported before the pump delivers the bolus. Only non-negative numbers shall be used.

An insulin pump agent may contain multiple bolus setting objects if it is designed to transmit multiple modalities of bolus delivery, in which case the modality shall be denoted in the Supplemental-Type attribute.

The Supplemental-Types attribute is used to distinguish the modality and the reason for delivering a bolus of insulin. Table 6 represents the possible values to express a fast or extended bolus delivery mode as well as correction or meal bolus delivery reason.

Certain combinations of Supplemental-Types are allowed. If it is desired to express that a meal bolus was delivered as a fast bolus, the SupplementalTypesList structure of the Supplemental-Types attribute should

Health informatics—Personal health device communication—Part 10419: Device Specialization—Insulin Pump

contain the two values MDC_INS_BOLUS_MEAL and MDC_INS_BOLUS_FAST. Similarly, a correction bolus delivered as an extended bolus should contain MDC_INS_BOLUS_CORR and MDC_INS_BOLUS_EXT in the SupplementalTypesList.

The modality values MDC_INS_BOLUS_FAST and MDC_INS_BOLUS_EXT shall not be combined in the SupplementalTypesList describing the bolus context. However, a multi-wave bolus modality may be represented by reporting separate current bolus setting objects with MDC_INS_BOLUS_FAST and MDC_INS_BOLUS_EXT at the same time.

It is not recommended to combine reason values MDC_INS_BOLUS_MEAL and MDC_INS_BOLUS_CORR in the SupplementalTypesList.

If an acceptable, existing nomenclature term is not available, MDC_INS_BOLUS_UNDETERMINED shall be used.

Туре	Definition
MDC_INS_BOLUS_FAST	Bolus insulin delivered as fast bolus
MDC_INS_BOLUS_EXT	Bolus insulin delivered as an extended bolus
MDC_INS_BOLUS_CORR	Bolus insulin delivered for a correction
MDC_INS_BOLUS_MEAL	Bolus insulin delivered for a meal
MDC INS BOLUS UNDETERMINED	Undetermined bolus

Table 6—Bolus context

The Metric-ID attribute is used to distinguish the bolus setting source and described in Table 7.

If an acceptable, existing nomenclature term is not available, MDC INS BOLUS OTHER shall be used.

Type	Definition
MDC_INS_BOLUS_MANUAL	Manual, user defined bolus
MDC_INS_BOLUS_RECOMMENDED	Recommended bolus (e.g., bolus calculator)
MDC_INS_BOLUS_MANUAL_CHANGE	Recommended bolus changed by a user
MDC_INS_BOLUS_COMMANDED	Commanded bolus (e.g., AP controller)
MDC INS BOLUS OTHER	Other source of holus

Table 7—Bolus source

The Unit-Code attribute shall be MDC_DIM_INTL_UNIT or MDC_DIM_INTL_UNIT_PER_HR, as appropriate (e.g., an MDC_INS_BOLUS_FAST shall use MDC_DIM_INTL_UNIT and an MDC_INS_BOLUS_EXT shall use MDC_DIM_INTL_UNIT_PER_HR).

The Measure-Active-Period attribute shall be used with the MDC_INS_BOLUS_EXT to describe the time duration of the extended bolus setting.

The current bolus setting numeric object does not support any methods, events, or other services.

See IEEE Std 11073-20601-2014 for descriptive explanations on the individual attributes as well as for information on attribute ID and attribute type.

6.7.3 Pending bolus delay

Table 8 summarizes the attributes of the pending bolus delay numeric object. The nomenclature code to identify the numeric class is MDC_MOC_VMO_METRIC_NU. The pending bolus delay numeric object may be present in the extended configuration.

Health informatics—Personal health device communication—Part 10419: Device Specialization—Insulin Pump

Table 8—Pending bolus delay numeric object attributes

A44.9. 4	Extended configuration			
Attribute name	Value	Qualifier		
Handle	See IEEE Std 11073-20601-2014.	M		
Туре	{MDC_PART_PHD_DM	M		
	MDC_INS_BOLUS_PENDING_DELAY}			
Supplemental-Types	See IEEE Std 11073-20601-2014.	NR		
Metric-Spec-Small	mss-avail-intermittent mss-avail-stored-data	M		
	mss-upd-aperiodic mss-acc-manager-initiated			
	mss-acc-agent-initiated mss-cat-manual mss-			
	cat-setting			
Metric-Structure-Small	See IEEE Std 11073-20601-2014.	О		
Measurement-Status	See IEEE Std 11073-20601-2014.	0		
Metric-Id	See IEEE Std 11073-20601-2014.	NR		
Metric-Id-List	See IEEE Std 11073-20601-2014.	NR		
Metric-Id-Partition	See IEEE Std 11073-20601-2014.	NR		
Unit-Code	MDC_DIM_MIN or MDC_DIM_SEC	M		
Attribute-Value-Map	See IEEE Std 11073-20601-2014.	С		
Source-Handle-Reference	Please see following text.	R		
Label-String	See IEEE Std 11073-20601-2014.	O		
Unit-Label String	See IEEE Std 11073-20601-2014.	O		
Absolute-Time-Stamp	See IEEE Std 11073-20601-2014.	С		
Base-Offset-Time-Stamp	See IEEE Std 11073-20601-2014.	R		
Relative-Time-Stamp	See IEEE Std 11073-20601-2014.	С		
HiRes-Time-Stamp	See IEEE Std 11073-20601-2014.	С		
Measure-Active-Period	See IEEE Std 11073-20601-2014.	NR		
Simple-Nu-Observed-Value	See IEEE Std 11073-20601-2014.	С		
Compound-Simple-Nu-Observed-Value	See /IEEE Std 11073-20601-2014.	С		
Basic-Nu-Observed-Value	See IEEE Std 11073-20601-2014.	M		
Compound-Basic-Nu-Observed-Value	See IEEE Std 11073-20601-2014.	С		
Nu-Observed-Value	See IEEE Std 11073-20601-2014.	С		
Compound-Nu-Observed-Value	See IEEE Std 11073-20601-2014.	С		
Accuracy	See IEEE Std 11073-20601-2014.	NR		

NOTE 1—See IEEE Std 11073-20601-2014 for information on whether an attribute is static or dynamic.

NOTE 2—See 6.3 for a description of the qualifiers.

The value reported in this object is the amount of time between programming a bolus and the actual bolus delivered. Only non-negative numbers shall be used.

If the Source-Handle-Reference is defined, it should point to the current bolus setting object.

The pending bolus delay numeric object does not support any methods, events, or other services.

See IEEE Std 11073-20601-2014 for descriptive explanations on the individual attributes as well as for information on attribute ID and attribute type.

6.7.4 Bolus delivered

Table 9 summarizes the attributes of the bolus delivered numeric object. The nomenclature code to identify the numeric class is MDC_MOC_VMO_METRIC_NU. The bolus delivered numeric object shall be supported by an insulin pump agent.

IEEE Std 11073-10419-2015
Health informatics—Personal health device communication—Part 10419: Device Specialization—Insulin Pump

Table 9—Bolus delivered numeric object attributes

Attribute name	Extended configuration	1	Standard configuration (Dev-Configuration-Id = 0x076C)	
Attribute name	Value	Qualifier	Value	Qualifier
Handle	See IEEE Std 11073-20601-2014.	M	1	M
Туре	{MDC_PART_PHD_DM MDC_INS_BOLUS}	M	{MDC_PART_PHD_DM MDC_INS_BOLUS}	M
Supplemental- Types	Please see following text in Table 10.	R	Attribute not initially present. If present follow IEEE Std 11073-20601-2014. Please see following text in Table 10.	О
Metric-Spec- Small	mss-avail-intermittent mss-avail- stored-data mss-upd-aperiodic mss-msmt-aperiodic mss-acc- manager-initiated mss-acc-agent- initiated mss-cat-calculation	М	mss-avail-intermittent mss-avail- stored-data mss-upd-aperiodic mss- msmt-aperiodic mss-acc-manager- initiated mss-acc-agent-initiated mss-cat-calculation	М
Metric- Structure-Small	See IEEE Std 11073-20601-2014.	О	Attribute not initially present. If present follow IEEE Std 11073-20601-2014.	О
Measurement- Status	See IEEE Std 11073-20601-2014.	О	Attribute not initially present. If present follow IEEE Std 11073-20601-2014.	О
Metric-Id	See IEEE Std 11073-20601-2014. Please see following text.	О	Attribute not initially present. If present follow IEEE Std 11073-20601-2014. Please see following text.	О
Metric-Id-List	See IEEE Std 11073-20601-2014.	О	Attribute not initially present. If present follow IEEE Std 11073-20601-2014.	С
Metric-Id- Partition	See IEEE Std 11073-20601-2014.	NR	Attribute not initially present. If present follow IEEE Std 11073-20601-2014.	NR
Unit-Code	MDC_DIM_INTL_UNIT	M	MDC_DIM_INTL_UNIT	M
Attribute-Value- Map	See IEEE Std 11073-20601-2014.	С	MDC_ATTR_NU_VAL_OBS_BAS IC, then MDC_ATTR_TIME_STAMP_BO	M
Source-Handle- Reference	See IEEE Std 11073-20601-2014. Please see following text.	R	Attribute not initially present. If present follow IEEE Std 11073-20601-2014.	О
Label-String	See IEEE Std 11073-20601-2014.	О	Attribute not initially present. If present follow IEEE Std 11073-20601-2014.	О
Unit-Label String	See IEEE Std 11073-20601-2014.	О	Attribute not initially present. If present follow IEEE Std 11073-20601-2014.	О
Absolute-Time- Stamp	See IEEE Std 11073-20601-2014.	С	Attribute not initially present. If present follow IEEE Std 11073-20601-2014.	С
Base-Offset- Time-Stamp	See IEEE Std 11073-20601-2014.	R	If fixed format is used and the standard configuration is not adjusted, this attribute is mandatory; otherwise, the conditions from ISO/IEEE 11073-20601a-2010 apply.	М
Relative-Time- Stamp	See IEEE Std 11073-20601-2014.	С	Attribute not initially present. If present follow IEEE Std 11073-20601-2014.	С

Health informatics—Personal health device communication—Part 10419: Device Specialization—Insulin Pump

Table 9—Bolus delivered numeric object attributes (continued)

Attribute name	Extended configuration		Standard configuration (Dev-Configuration-Id = 0x07	(6C)
	Value	Qualifier	Value	Qualifier
HiRes-Time- Stamp	See IEEE Std 11073-20601-2014.	С	Attribute not initially present. If present follow IEEE Std 11073-20601-2014.	С
Measure- Active-Period	See IEEE Std 11073-20601-2014.	NR	Attribute not initially present. If present follow IEEE Std 11073-20601-2014.	NR
Simple-Nu- Observed-Value	See IEEE Std 11073-20601-2014.	С	Attribute not initially present. If present follow IEEE Std 11073-20601-2014.	С
Compound- Simple-Nu- Observed-Value	See IEEE Std 11073-20601-2014.	С	Attribute not initially present. If present follow IEEE Std 11073-20601-2014	С
Basic-Nu- Observed-Value	See IEEE Std 11073-20601-2014.	M	If fixed format is used and the standard configuration is unchanged, this attribute is mandatory; otherwise, the conditions from IEEE Std 11073-20601-2014 apply.	М
Compound- Basic-Nu- Observed-Value	See IEEE Std 11073-20601-2014.	С	Attribute not initially present. If present follow IEEE Std 11073-20601-2014	С
Nu-Observed- Value	See IEEE Std 11073-20601-2014.	С	Attribute not initially present. If present follow IEEE Std 11073-20601-2014.	С
Compound-Nu- Observed-Value	See IEEE Std 11073-20601-2014.	С	Attribute not initially present. If present follow IEEE Std 11073-20601-2014.	С
Accuracy	See IEEE Std 11073-20601-2014.	NR	Attribute not initially present. If present follow IEEE Std 11073-20601-2014.	NR

NOTE 1—See IEEE Std 11073-20601-2014 for information on whether an attribute is static or dynamic.

NOTE 2—See 6.3 for a description of the qualifiers.

For an insulin pump agent with standard configuration the AttrValMap structure (see IEEE Std 11073-20601-2014) of the Attribute-Value-Map attribute shall contain the attribute ID and attribute length information of the Basic-Nu-Observed-Value and Base-Offset-Time-Stamp attribute in the same order as indicated in Table 9.

The observed value reported in this object is the amount of the insulin effectively delivered as a bolus and is reported after the pump has delivered the bolus. Only non-negative numbers shall be used.

The Supplemental-Types attribute is used to distinguish the modality of a particular bolus and the reason for delivering a bolus of insulin. Table 10 represents the possible values to express a fast or extended bolus delivery mode as well as correction or meal bolus delivery reason.

Certain combinations of Supplemental-Types are allowed. If it is desired to express that a meal bolus was delivered as a fast bolus, the SupplementalTypesList structure of the Supplemental-Types attribute should contain the two values MDC_INS_BOLUS_MEAL and MDC_INS_BOLUS_FAST. Similarly, a correction bolus delivered as an extended bolus should contain MDC_INS_BOLUS_CORR and MDC_INS_BOLUS_EXT in the SupplementalTypesList.

Health informatics—Personal health device communication—Part 10419: Device Specialization—Insulin Pump

The modality values MDC_INS_BOLUS_FAST and MDC_INS_BOLUS_EXT shall not be combined in the SupplementalTypesList describing the bolus context. However, a multi-wave bolus modality may be represented by reporting separate bolus objects with MDC_INS_BOLUS_FAST and MDC_INS_BOLUS_EXT at the same time.

The modality values MDC_INS_BOLUS_MEAL and MDC_INS_BOLUS_CORR shall not be combined in the SupplementalTypesList.

If an acceptable, existing nomenclature term is not available, MDC_INS_BOLUS_UNDETERMINED shall be used.

Table 10 —Bolus context

Type	Definition
MDC_INS_BOLUS_FAST	Bolus insulin delivered as fast bolus
MDC_INS_BOLUS_EXT	Bolus insulin delivered as an extended bolus
MDC_INS_BOLUS_CORR	Bolus insulin delivered for a correction
MDC_INS_BOLUS_MEAL	Bolus insulin delivered for a meal
MDC_INS_BOLUS_UNDETERMINED	Undetermined bolus

If the Source-Handle-Reference is defined, it should point to the current bolus setting object.

The bolus delivered numeric object does not support any methods, events, or other services.

See IEEE Std 11073-20601-2014 for descriptive explanations on the individual attributes as well as for information on attribute ID and attribute type.

6.7.5 Current basal rate setting

Table 11 summarizes the attributes of the current basal rate setting numeric object. The nomenclature code to identify the numeric class is MDC_MOC_VMO_METRIC_NU. The current basal rate setting numeric object shall be supported by an insulin pump agent.

IEEE Std 11073-10419-2015
Health informatics—Personal health device communication—Part 10419: Device Specialization—Insulin Pump

Table 11 —Current basal rate setting numeric object attributes

Attribute	Extended configuration		Standard configuration (Dev-Configuration-Id = 0x0'	76C)
name	Value	Qualifier	Value	Qualifier
Handle	See IEEE Std 11073-20601-2014.	M	2	M
Type	{MDC_PART_PHD_DM MDC_INS_BASAL_RATE_SETTI NG}	М	{MDC_PART_PHD_DM MDC_INS_BASAL_RATE_SETTI NG}	M
Supplemental- Types	Please see following text in Table 12	М	MDC_INS_BASAL_PRGM	M
Metric-Spec- Small	mss-avail-intermittent mss-avail- stored-data mss-upd-aperiodic mss- acc-manager-initiated mss-acc- agent-initiated mss-cat-setting	M	mss-avail-intermittent mss-avail- stored-data mss-upd-aperiodic mss- acc-manager-initiated mss-acc- agent-initiated mss-cat-setting	М
Metric- Structure-Small	See IEEE Std 11073-20601-2014.	NR	Attribute not initially present. If present follow IEEE Std 11073-20601-2014.	NR
Measurement- Status	See IEEE Std 11073-20601-2014.	NR	Attribute not initially present. If present follow IEEE Std 11073-20601-2014.	NR
Metric-Id	Please see following text in Table 13.	0	MDC_INS_BASAL_DEVICE	M
Metric-Id-List	See IEEE Std 11073-20601-2014.	NR	Attribute not initially present. If present follow IEEE Std 11073-20601-2014.	NR
Metric-Id- Partition	See IEEE Std 11073-20601-2014.	NR	Attribute not initially present. If present follow IEEE Std 11073-20601-2014.	NR
Unit-Code	MDC_DIM_INTL_UNIT_PER_HR	M	MDC_DIM_INTL_UNIT_PER_HR	M
Attribute-Value- Map	See IEEE Std 11073-20601-2014.	С	MDC_ATTR_NU_VAL_OBS_BAS IC, then MDC_ATTR_TIME_STAMP_BO.	M
Source-Handle- Reference	See IEEE Std 11073-20601-2014.	NR	Attribute not initially present. If present follow IEEE Std 11073-20601-2014.	NR
Label-String	See IEEE Std 11073-20601-2014.	0	Attribute not initially present. If present follow IEEE Std 11073-20601-2014.	О
Unit- LabelString	See IEEE Std 11073-20601-2014.	0	Attribute not initially present. If present follow IEEE Std 11073-20601-2014.	О
Absolute-Time- Stamp	See IEEE Std 11073-20601-2014.	С	Attribute not initially present. If present follow IEEE Std 11073-20601-2014.	С
Base-Offset- Time-Stamp	See IEEE Std 11073-20601-2014.	R	If fixed format is used and the standard configuration is unchanged, this attribute is mandatory; otherwise, the conditions from IEEE Std 11073-20601-2014 apply.	M
Relative-Time- Stamp	See IEEE Std 11073-20601-2014.	С	Attribute not initially present. If present follow IEEE Std 11073-20601-2014.	С
HiRes-Time- Stamp	See IEEE Std 11073-20601-2014.	С	Attribute not initially present. If present follow IEEE Std 11073-20601-2014.	С
Measure- Active-Period	See IEEE Std 11073-20601-2014.	NR	Attribute not initially present. If present follow IEEE Std 11073-20601-2014.	NR
Simple-Nu- Observed-Value	See IEEE Std 11073-20601-2014.	С	Attribute not initially present. If present follow IEEE Std 11073-20601-2014.	С

Health informatics—Personal health device communication—Part 10419: Device Specialization—Insulin Pump

Table 11—Current basal rate setting numeric object attributes (continued)

Attribute	Extended configuration	l	Standard configuration (Dev-Configuration-Id = 0x076C)	
name	Value	Qualifier	Value	Qualifier
Compound- Simple-Nu- Observed-Value	See /IEEE Std 11073-20601-2014.	С	Attribute not initially present. If present follow IEEE Std 11073-20601-2014	С
Basic-Nu- Observed-Value	See IEEE Std 11073-20601-2014.	R	If fixed format is used and the standard configuration is unchanged, this attribute is mandatory; otherwise, the conditions from IEEE Std 11073-20601-2014 apply.	М
Compound- Basic-Nu- Observed-Value	See IEEE Std 11073-20601-2014.	С	Attribute not initially present. If present follow IEEE Std 11073-20601-2014	С
Nu-Observed- Value	See IEEE Std 11073-20601-2014.	С	Attribute not initially present. If present follow IEEE Std 11073-20601-2014.	С
Compound-Nu- Observed-Value	See IEEE Std 11073-20601-2014.	С	Attribute not initially present. If present follow IEEE Std 11073-20601-2014.	С
Accuracy	See IEEE Std 11073-20601-2014.	NR	Attribute not initially present. If present follow IEEE Std 11073-20601-2014.	NR

NOTE 1—See IEEE Std 11073-20601-2014 for information on whether an attribute is static or dynamic.

NOTE 2—See 6.3 for a description of the qualifiers.

For an insulin pump agent with standard configuration the AttrValMap structure (see IEEE Std 11073-20601-2014) of the Attribute-Value-Map attribute shall contain the attribute ID and attribute length information of the Basic-Nu-Observed-Value and Base-Offset-Time-Stamp attribute in the same order as indicated in Table 11.

The value reported in this object is the rate of insulin requested to be delivered as a basal and is reported prior to the delivery of the basal. Only non-negative numbers shall be used.

The Supplemental-Types attribute is used to distinguish additional context of a particular basal rate. Table 12 represents the possible values to express a programmed, temporary absolute or temporary relative basal rate.

If an acceptable, existing nomenclature term is not available, MDC_INS_BASAL_UNDETERMINED shall be used.

Table 12—Basal rate context

Type	Definition
MDC_INS_BASAL_PRGM	Rate set by the active basal profile schedule setting
MDC_INS_BASAL_TEMP_ABS	Rate set by an absolute temporary change
MDC_INS_BASAL_TEMP_REL	Rate set by a relative temporary change
MDC_INS_BASAL_UNDETERMINED	Undetermined basal rate context

The Metric-Id attribute is used to distinguish the source of a particular basal rate setting. Table 13 represents the possible values to express a basal rate setting by the device, a remote or an artificial pancreas controller.

Health informatics—Personal health device communication—Part 10419: Device Specialization—Insulin Pump

Table 13 — Basal rate source

Туре	Definition
MDC_INS_BASAL_DEVICE	Rate set by the insulin pump
MDC_INS_BASAL_REMOTE	Rate set by a remote controller
MDC_INS_BASAL_AP_CTRL	Rate set by an artificial pancreas controller
MDC_INS_BASAL_OTHER	Rate set by an other source

The current basal rate setting numeric object does not support any methods, events, or other services.

See IEEE Std 11073-20601-2014 for descriptive explanations on the individual attributes as well as for information on attribute ID and attribute type.

6.7.6 Basal delivered

Table 14 defines the attributes for the object that represents the basal delivered. The nomenclature code to identify the numeric class is MDC_MOC_VMO_METRIC_NU. The basal delivered object may be supported by an insulin pump agent with extended configuration.

Table 14—Basal delivered numeric object attributes

Attribute name	Extended configuration	Qualifier
Handle	See IEEE Std 11073-20601-2014.	M
Type	{MDC_PART_PHD_DM	M
	MDC_INS_BASAL}	
Supplemental-Types	See IEEE Std 11073-20601-2014.	O
Metric-Spec-Small	mss-avail-intermittent mss-avail-stored-data mss-	M
	upd-aperiodic mss-msmt-aperiodic mss-acc-	
	manager-initated mss-acc-agent-initiated mss-cat-	
	calculation	
Metric-Structure-Small	See IEEE Std 11073-20601-2014.	NR
Measurement-Status	See IEEE Std 11073-20601-2014.	NR
Metric-Id	See IEEE Std 11073-20601-2014.	0
Metric-Id-List	See IEEE Std 11073-20601-2014.	C
Metric-Id-Partition	See IEEE Std 11073-20601-2014.	О
Unit-Code	MDC_DIM_INTL_UNIT	M
Attribute-Value-Map	See IEEE Std 11073-20601-2014.	C
Source-Handle-Reference	See IEEE Std 11073-20601-2014.	О
Label-String	See IEEE Std 11073-20601-2014.	O
Unit-LabelString	See IEEE Std 11073-20601-2014.	О
Absolute-Time-Stamp	See IEEE Std 11073-20601-2014.	C
Base-Offset-Time-Stamp	See IEEE Std 11073-20601-2014.	R
Relative-Time-Stamp	See IEEE Std 11073-20601-2014.	С
HiRes-Time-Stamp	See IEEE Std 11073-20601-2014.	C
Measure-Active-Period	See IEEE Std 11073-20601-2014.	NR
Simple-Nu-Observed-Value	See IEEE Std 11073-20601-2014.	C
Compound-Simple-Nu-Observed-Value	See IEEE Std 11073-20601-2014.	С
Basic-Nu-Observed-Value	See IEEE Std 11073-20601-2014.	R
Compound-Basic-Nu-Observed-Value	See IEEE Std 11073-20601-2014.	С
Nu-Observed-Value	See IEEE Std 11073-20601-2014.	С
Compound-Nu-Observed-Value	See IEEE Std 11073-20601-2014.	С
Accuracy	See IEEE Std 11073-20601-2014.	NR

NOTE 1—See IEEE Std 11073-20601-2014 for information on whether an attribute is static or dynamic.

NOTE 2—See 6.3 for a description of the qualifiers.

Health informatics—Personal health device communication—Part 10419: Device Specialization—Insulin Pump

The observed value reported in this object is the amount of the insulin effectively delivered as a basal and is reported after the pump has delivered the basal. The insulin pump is responsible for determining the frequency of the basal delivered event to be reported (i.e., insulin pump may report the amount of insulin delivered at the end of a current basal rate setting or intermittently as micro-boluses). Only non-negative numbers shall be used.

The basal delivered numeric object does not support any methods, events, or other services.

See IEEE Std 11073-20601-2014 for descriptive explanations on the individual attributes as well as information on attribute ID and attribute type.

6.7.7 Basal rate schedule setting

Table 15 defines the attributes for the numeric object that represents the basal rate for the schedule.

The basal rate schedule setting object may be supported by an insulin pump agent with extended configuration.

Table 15—Basal rate schedule setting numeric object attributes

Attribute name	Extended configuration	Qualifier
Handle	See IEEE Std 11073-20601-2014.	M
Type	{MDC_PART_PHD_DM	M
	MDC_INS_BASAL_RATE_SCHED}	
Supplemental-Types	See IEEE Std 11073-20601-2014.	NR
Metric-Spec-Small	mss-avail-intermittent mss-avail-stored-data mss-upd-	M
	aperiodic mss-msmt-aperiodic mss-acc-agent-initiated	
	l mss-cat-setting	
Metric-Structure-Small	See IEEE Std 11073-20601-2014.	NR
Measurement-Status	See IEEE Std 11073-20601-2014.	NR
Metric-Id	See IEEE Std 11073-20601-2014.	NR
Metric-Id-List	See IEEE Std 11073-20601-2014.	NR
Metric-Id-Partition	See IEEE Std 11073-20601-2014.	NR
Unit-Code	MDC_DIM_INTL_UNIT_PER_HR.	M
Attribute-Value-Map	See IEEE Std 11073-20601-2014.	С
Source-Handle-Reference	See IEEE Std 11073-20601-2014.	NR
Label-String	See IEEE Std 11073-20601-2014.	О
Unit-LabelString	See IEEE Std 11073-20601-2014.	О
Absolute-Time-Stamp	See IEEE Std 11073-20601-2014.	С
Base-Offset-Time-Stamp	See IEEE Std 11073-20601-2014.	R
Relative-Time-Stamp	See IEEE Std 11073-20601-2014.	С
HiRes-Time-Stamp	See IEEE Std 11073-20601-2014.	С
Measure-Active-Period	See IEEE Std 11073-20601-2014.	NR
Simple-Nu-Observed-Value	See IEEE Std 11073-20601-2014.	M
Compound-Simple-Nu-Observed-Value	See IEEE Std 11073-20601-2014.	NR
Basic-Nu-Observed-Value	See IEEE Std 11073-20601-2014.	С
Compound-Basic-Nu-Observed-Value	See IEEE Std 11073-20601-2014.	С
Nu-Observed-Value	See IEEE Std 11073-20601-2014.	С
Compound-Nu-Observed-Value	See IEEE Std 11073-20601-2014.	С
Accuracy	See IEEE Std 11073-20601-2014.	NR

NOTE 1—See IEEE Std 11073-20601-2014 for information on whether an attribute is static or dynamic.

NOTE 2—See 6.3 for a description of the qualifiers.

Health informatics—Personal health device communication—Part 10419: Device Specialization—Insulin Pump

The value reported in this object is the basal rate schedule setting that represents a basal rate for a specifc schedule setting. The basal rate schedule setting shall be made available only through schedule objects.

The unit of a basal rate schedule setting is international unit per hour (IU/h).

See IEEE Std 11073-20601-2014 for descriptive explanations on the individual attributes as well as information on attribute ID and attribute type.

6.7.8 Insulin to carbohydrate ratio schedule setting

Table 16 defines the attributes for the object that represents the insulin to carbohydrate ratio (I:CHO) schedule values, which include the amount of carbohydrates in grams covered or disposed of by a single unit of insulin. The I:CHO schedule setting may be supported by the agent with extended configuration.

Table 16—Insulin to carbohydrate ratio schedule setting attributes

Attribute name	Extended configuration	Qualifier
Handle	See IEEE Std 11073-20601-2014.	M
Type	{MDC_PART_PHD_DM	M
	MDC_INS_I2CHO_SCHED}	
Supplemental-Types	See IEEE Std 11073-20601-2014.	NR
Metric-Spec-Small	mss-avail-intermittent mss-avail-stored-data mss-upd-	M
	aperiodic mss-msmt-aperiodic mss-acc-agent-initiated	
	mss-cat-setting	
Metric-Structure-Small	See IEEE Std 11073-20601-2014.	NR
Measurement-Status	See IEEE Std 11073-20601-2014.	NR
Metric-Id	See IEEE Std 11073-20601-2014.	NR
Metric-Id-List	See IEEE Std 11073-20601-2014.	NR
Metric-Id-Partition	See IEEE Std 11073-20601-2014.	NR
Unit-Code	MDC_DIM_G.	M
	Please see following text.	
Attribute-Value-Map	See IEEE Std 11073-20601-2014.	C
Source-Handle-Reference	See IEEE Std 11073-20601-2014.	NR
Label-String	See IEEE Std 11073-20601-2014.	О
Unit-LabelString	See IEEE Std 11073-20601-2014.	О
Absolute-Time-Stamp	See IEEE Std 11073-20601-2014.	С
Base-Offset-Time-Stamp	See IEEE Std 11073-20601-2014.	R
Relative-Time-Stamp	See IEEE Std 11073-20601-2014.	С
HiRes-Time-Stamp	See IEEE Std 11073-20601-2014.	С
Measure-Active-Period	See IEEE Std 11073-20601-2014.	NR
Simple-Nu-Observed-Value	See IEEE Std 11073-20601-2014.	M
Compound-Simple-Nu-Observed-Value	See IEEE Std 11073-20601-2014.	NR
Basic-Nu-Observed-Value	See IEEE Std 11073-20601-2014.	С
Compound-Basic-Nu-Observed-Value	See IEEE Std 11073-20601-2014.	С
Nu-Observed-Value	See IEEE Std 11073-20601-2014.	С
Compound-Nu-Observed-Value	See IEEE Std 11073-20601-2014.	С
Accuracy	See IEEE Std 11073-20601-2014.	NR

NOTE 1—See IEEE Std 11073-20601-2014 for information on whether an attribute is static or dynamic.

NOTE 2—See 6.3 for a description of the qualifiers.

The value reported in this object is the I:CHO for a specifc schedule setting. The I:CHO schedule setting shall be made available only through Schedule objects.

Health informatics—Personal health device communication—Part 10419: Device Specialization—Insulin Pump

The unit of a I:CHO schedule setting is grams (g).

See IEEE Std 11073-20601-2014 for descriptive explanations on the individual attributes as well as information on attribute ID and attribute type.

6.7.9 Insulin sensitivity factor schedule setting

Table 17 defines the attributes for the object that represents the insulin sensitivity factor (ISF) schedule values. The ISF schedule setting may be supported by the agent with extended configuration.

Table 17—Insulin sensitivity factor schedule setting attributes

Attribute name	Extended configuration	Qualifier
Handle	See IEEE Std 11073-20601-2014.	M
Type	{MDC_PART_PHD_DM	M
	MDC_INS_ISF_SCHED}	
Supplemental-Types	See IEEE Std 11073-20601-2014.	NR
Metric-Spec-Small	mss-avail-intermittent mss-avail-stored-data mss-	M
	upd-aperiodic mss-msmt-aperiodic mss-acc-agent-	
	initiated mss-cat-setting	
Metric-Structure-Small	See IEEE Std 11073-20601-2014.	NR
Measurement-Status	See IEEE Std 11073-20601-2014.	NR
Metric-Id	See IEEE Std 11073-20601-2014.	NR
Metric-Id-List	See IEEE Std 11073-20601-2014.	NR
Metric-Id-Partition	See IEEE Std 11073-20601-2014.	NR
Unit-Code	MDC_DIM_MILLI_MOLE_PER_L or	M
	MDC_DIM_MILLI_G_PER_DL.	
	Please see following text.	
Attribute-Value-Map	See IEEE Std 11073-20601-2014	С
Source-Handle-Reference	See IEEE Std 11073-20601-2014.	NR
Label-String	See IEEE Std 11073-20601-2014.	О
Unit-LabelString	See IEEE Std 11073-20601-2014.	О
Absolute-Time-Stamp	See IEEE Std 11073-20601-2014.	C
Base-Offset-Time-Stamp	See IEEE Std 11073-20601-2014.	R
Relative-Time-Stamp	See IEEE Std 11073-20601-2014.	С
HiRes-Time-Stamp	See IEEE Std 11073-20601-2014.	C
Measure-Active-Period	See IEEE Std 11073-20601-2014.	NR
Simple-Nu-Observed-Value	See IEEE Std 11073-20601-2014.	M
Compound-Simple-Nu-Observed-Value	See IEEE Std 11073-20601-2014.	NR
Basic-Nu-Observed-Value	See IEEE Std 11073-20601-2014.	C
Compound-Basic-Nu-Observed-Value	See IEEE Std 11073-20601-2014.	С
Nu-Observed-Value	See IEEE Std 11073-20601-2014.	С
Compound-Nu-Observed-Value	See IEEE Std 11073-20601-2014.	С
Accuracy	See IEEE Std 11073-20601-2014.	NR

NOTE 1—See IEEE Std 11073-20601-2014 for information on whether an attribute is static or dynamic.

NOTE 2—See 6.3 for a description of the qualifiers.

The observed value reported in this object is ISF for a specifc schedule setting. The ISF schedule setting shall be made available only through schedule objects.

The unit of ISF schedule setting is millimole per liter (mmol/L) or milligrams per deciliter (mg/dL).

See IEEE Std 11073-20601-2014 for descriptive explanations on the individual attributes as well as information on attribute ID and attribute type.

Health informatics—Personal health device communication—Part 10419: Device Specialization—Insulin Pump

6.7.10 Insulin reservoir remaining

Table 18 defines the attributes for the object that represents the insulin remaining in the reservoir. The insulin reservoir remaining object may be supported by an insulin pump agent with extended configuration.

Table 18—Insulin reservoir remaining numeric object attributes

Attribute name	Extended configuration	Qualifier
Handle	See IEEE Std 11073-20601-2014.	M
Type	{MDC_PART_PHD_DM M	
	MDC_INS_RESERVOIR}	
Supplemental-Types	See IEEE Std 11073-20601-2014.	О
Metric-Spec-Small	mss-avail-intermittent mss-avail-stored-data mss-	M
	upd-aperiodic mss-msmt-aperiodic mss-acc-	
	manager-initiated mss-acc-agent-initiated mss-cat-	
	calculation	
Metric-Structure-Small	See IEEE Std 11073-20601-2014.	NR
Measurement-Status	See IEEE Std 11073-20601-2014.	NR
Metric-Id	See IEEE Std 11073-20601-2014.	0
Metric-Id-List	See IEEE Std 11073-20601-2014.	С
Metric-Id-Partition	See IEEE Std 11073-20601-2014.	О
Unit-Code	MDC_DIM_INTL_UNIT	M
Attribute-Value-Map	See IEEE Std 11073-20601-2014.	С
Source-Handle-Reference	See IEEE Std 11073-20601-2014.	NR
Label-String	See IEEE Std 11073-20601-2014.	О
Unit-LabelString	See IEEE Std 11073-20601-2014.	О
Absolute-Time-Stamp	See IEEE Std 11073-20601-2014.	С
Base-Offset-Time-Stamp	See IEEE Std 11073-20601-2014.	R
Relative-Time-Stamp	See IEEE Std 11073-20601-2014.	C
HiRes-Time-Stamp	See IEEE Std 11073-20601-2014.	C
Measure-Active-Period	See IEEE Std 11073-20601-2014.	NR
Simple-Nu-Observed-Value	See IEEE Std 11073-20601-2014.	C
Compound-Simple-Nu-Observed-Value	See IEEE Std 11073-20601-2014.	С
Basic-Nu-Observed-Value	See IEEE Std 11073-20601-2014.	R
Compound-Basic-Nu-Observed-Value	See IEEE Std 11073-20601-2014.	С
Nu-Observed-Value	See IEEE Std 11073-20601-2014.	С
Compound-Nu-Observed-Value	See IEEE Std 11073-20601-2014.	С
Accuracy	See IEEE Std 11073-20601-2014.	NR

NOTE 1—See IEEE Std 11073-20601-2014 for information on whether an attribute is static or dynamic.

NOTE 2—See 6.3 for a description of the qualifiers.

The observed value reported in this object is the amount of the insulin remaining in the reservoir. The insulin pump is responsible for determining the frequency of the reservoir remaining event to be reported (e.g., the insulin reservoir remaining may be reported at the same time as the bolus delivered or basal delivered measurements). Only non-negative numbers shall be used.

The insulin reservoir remaining numeric object does not support any methods, events, or other services.

See IEEE Std 11073-20601-2014 for descriptive explanations on the individual attributes as well as information on attribute ID and attribute type.

Health informatics—Personal health device communication—Part 10419: Device Specialization—Insulin Pump

6.7.11 Insulin concentration

Table 19 defines the attributes for the object that represents the insulin concentration. The insulin concentration object may be supported by an insulin pump agent with extended configuration.

Table 19—Insulin concentration numeric object attributes

Attribute name	Extended configuration	Qualifier
Handle	See IEEE Std 11073-20601-2014.	M
Type	{MDC_PART_PHD_DM	M
	MDC_INS_CONC}	
Supplemental-Types	See IEEE Std 11073-20601-2014.	NR
Metric-Spec-Small	mss-avail-intermittent mss-avail-stored-data mss-	M
	upd-aperiodic mss-acc-manager-initatied mss-acc-	
	agent-initiated mss-cat-setting	
Metric-Structure-Small	See IEEE Std 11073-20601-2014.	NR
Measurement-Status	See IEEE Std 11073-20601-2014.	NR
Metric-Id	See IEEE Std 11073-20601-2014.	O
Metric-Id-List	See IEEE Std 11073-20601-2014.	C
Metric-Id-Partition	See IEEE Std 11073-20601-2014.	О
Unit-Code	MDC_DIM_INTL_UNIT_PER_ML, or	M
	MDC_DIM_INTL_UNIT_PER_L, or	
	MDC_DIM_INTL_UNIT_PER_M_CUBE or	
	MDC_DIM_INTL_UNIT_PER_CM_CUBE.	
Attribute-Value-Map	See IEEE Std 11073-20601-2014.	C
Source-Handle-Reference	See IEEE Std 11073-20601-2014.	NR
Label-String	See IEEE Std 11073-20601-2014.	O
Unit-LabelString	See IEEE Std 11073-20601-2014.	O
Absolute-Time-Stamp	See IEEE Std 11073-20601-2014.	C
Base-Offset-Time-Stamp	See IEEE Std 11073-20601-2014.	R
Relative-Time-Stamp	See IEEE Std 11073-20601-2014.	C
HiRes-Time-Stamp	See IEEE Std 11073-20601-2014.	С
Measure-Active-Period	See IEEE Std 11073-20601-2014.	NR
Simple-Nu-Observed-Value	See IEEE Std 11073-20601-2014.	С
Compound-Simple-Nu-Observed-Value	See IEEE Std 11073-20601-2014.	С
Basic-Nu-Observed-Value	See IEEE Std 11073-20601-2014.	R
Compound-Basic-Nu-Observed-Value	See IEEE Std 11073-20601-2014.	С
Nu-Observed-Value	See IEEE Std 11073-20601-2014.	С
Compound-Nu-Observed-Value	See IEEE Std 11073-20601-2014.	С
Accuracy	See IEEE Std 11073-20601-2014.	NR

NOTE 1—See IEEE Std 11073-20601-2014 for information on whether an attribute is static or dynamic.

NOTE 2—See 6.3 for a description of the qualifiers.

The observed value reported in this object is the insulin concentration. Only non-negative numbers shall be used.

The insulin concentration numeric object does not support any methods, events, or other services.

See IEEE Std 11073-20601-2014 for descriptive explanations on the individual attributes as well as information on attribute ID and attribute type.

6.8 Real-time sample array objects

Real-time sample array objects are not required by this standard.

Health informatics—Personal health device communication—Part 10419: Device Specialization—Insulin Pump

6.9 Enumeration objects

6.9.1 General

The enumeration objects listed in this clause represent the enumeration observations that can be produced by an insulin pump agent. The nomenclature code to identify an enumeration object class is MDC_MOC_VMO_METRIC_ENUM. Subclauses 6.9.2, 6.9.3, and 6.9.4 define the precise definitions for operational status, general device status, and specific insulin pump status enumeration objects.

Enumeration objects do not support any methods, events, or other services.

6.9.2 Operational status

Table 20 defines the attributes for the object that represents the operational status. The operational status enumeration object may be supported by an insulin pump agent with extended configuration.

Table 20—Operational status enumeration object attributes

Attribute name	Extended configuration	Qualifier
Handle	See IEEE Std 11073-20601-2014.	M
Type	MDC_PART_PHD_DM	M
	MDC_INS_PUMP_OP_STAT	
Supplemental-Types	See IEEE Std 11073-20601-2014.	NR
Metric-Spec-Small	mss-avail-intermittent mss-avail-stored-data mss-	M
	upd-aperiodic mss-acc-manager-initiated mss-acc-	
	agent-initiated	
Metric-Structure-Small	See IEEE Std 11073-20601-2014.	NR
Measurement-Status	See IEEE Std 11073-20601-2014.	NR
Metric-Id	See IEEE Std 11073-20601-2014.	NR
Metric-Id-List	See IEEE Std 11073-20601-2014.	NR
Metric-Id-Partition	See IEEE Std 11073-20601-2014.	NR
Unit-Code	See IEEE Std 11073-20601-2014.	NR
Attribute-Value-Map	See IEEE Std 11073-20601-2014.	C
Source-Handle-Reference	Please see following text.	NR
Label-String	See IEEE Std 11073-20601-2014.	NR
Unit-LabelString	See IEEE Std 11073-20601-2014.	NR
Absolute-Time-Stamp	See IEEE Std 11073-20601-2014.	C
Base-Offset-Time-Stamp	See IEEE Std 11073-20601-2014.	R
Relative-Time-Stamp	See IEEE Std 11073-20601-2014.	C
HiRes-Time-Stamp	See IEEE Std 11073-20601-2014.	C
Measure-Active-Period	See IEEE Std 11073-20601-2014.	NR
Enum-Observed-Value-Simple-OID	See IEEE Std 11073-20601-2014.	NR
Enum-Observed-Value-Simple-Bit-Str	See IEEE Std 11073-20601-2014.	NR
Enum-Observed-Value-Basic-Bit-Str	Please see following text.	R
Enum-Observed-Value-Simple-Str	See IEEE Std 11073-20601-2014.	NR
Enum-Observed-Value	See IEEE Std 11073-20601-2014.	NR
Enum-Observed-Value-Partition	See IEEE Std 11073-20601-2014.	NR

NOTE 1—See IEEE Std 11073-20601-2014 for information on whether an attribute is static or dynamic.

NOTE 2—See 6.3 for a description of the qualifiers.

The observed value reported in this object is the operational status.

Health informatics—Personal health device communication—Part 10419: Device Specialization—Insulin Pump

Because these are essentially flags, the Unit-Code attribute is not appropriate for this object. Similarly, Source-Handle-Reference is inappropriate, as this object monitors the status of the equipment.

The explicit expression of the existence of annunciations is realized by the setting of the appropriate bit in the Enum-Observed-Value-Basic-Bit-Str attribute, as defined in Table 21. If a manager supports the interpretation of this object, it shall be able to interpret the entire set of presented conditions. An agent shall report each time there is a change in status in a condition. An agent may implement any subset of these same conditions. Note that a manager shall interpret these bits only within the context of this attribute and only within this device specialization as other specializations may use corresponding terms for different purposes.

Table 21 defines the nomenclature terms for the operational and therapy conditions. If an acceptable, existing nomenclature term is not available, insulin-device-op-undetermined shall be used.

Table 21 —Mapping of operational and therapy conditions to object Bit-Str attribute

Operational and therapy condition	OpStat mnemonic
Agent is in an undetermined operating condition.	insulin-device-op-undetermined
Agent is powering off (e.g., when the user has requested the insulin	insulin-device-op-off
pump to turn off). No functionality available (e.g., no insulin delivery	
and no configuration of the device is possible).	
Agent is in standby and most functionality can resume faster than from	insulin-device-op-standby
off condition (e.g., the device is in a condition to save power).	
Agent is performing initialization and preparing for infusion therapy	insulin-device-op-preparing
(e.g., the insulin delivery device rewinds the piston rod to enable the	
insertion of a new reservoir or is detecting the position of the plunger	
in the reservoir).	
Agent is preparing the delivery mechanism (e.g., filling the infusion set	insulin-device-op-priming
fluidic path with insulin after replacement of the reservoir and/or	
infusion set).	
Agent is ready to deliver therapy but waiting for an interaction (e.g.,	insulin-device-op-waiting
waiting for the infusion set to be connected to the body after priming or	
waiting for a user confirmation).	
Agent is ready to administer insulin infusion therapy. The current	insulin-device-op-ready
therapy condition is further indicated by bits 9-11.	
Agent has stopped therapy and insulin delivery (e.g., bolus nor basal	insulin-device-therapy-stop
insulin is being administered), and current bolus setting is canceled.	
Agent has paused therapy momentarily and insulin delivery suspended	insulin-device-therapy-pause
(e.g., during preparing operation condition, a therapy pause can be used	
to bridge the time which is typically limited to several minutes during a	
reservoir change. When the agent returns to therapy run the remaining	
insulin amount is delivered).	
Agent is administerating and delivering insulin to patient. The current	insulin-device-therapy-run
operation condition shall indicate "insulin-device-op-ready" (bit 6).	

NOTE 1—The bits in the table above are defined to: 0 = False and 1 = True.

NOTE 2—The specific bit mappings of OpStat are defined in Annex B.

NOTE 3—All bits not defined in Table 21 or Annex B are reserved for future use.

NOTE 4—An agent is not required to implement all the features specified in Table 21.

At least one current therapy condition bit (insulin-device-therapy-stop, insulin-device-therapy-pause, or insulin-device-therapy-run) shall be set for any of the operating conditions, but the insulin-device-therapy-run condition shall only be set when the insulin-device-op-ready condition is set.

Health informatics—Personal health device communication—Part 10419: Device Specialization—Insulin Pump

See IEEE Std 11073-20601-2014 for descriptive explanations on the individual attributes as well as information on attribute ID and attribute type.

6.9.3 PHD DM status

The PHD DM status object allows generic device events to be recorded in order to track important events for the user and troubleshooting information for manufacturers. Table 22 summarizes the attributes of the PHD DM status enumeration object. The PHD DM status enumeration object may be supported by an insulin pump agent.

Table 22—PHD DM status enumeration object attributes

Attribute name	Extended configuration	Qualifier
Handle	See IEEE Std 11073-20601-2014.	M
Type	MDC_PART_PHD_DM	M
	MDC_PHD_DM_DEV_STAT	
Supplemental-Types	See IEEE Std 11073-20601-2014.	NR
Metric-Spec-Small	mss-avail-intermittent, mss-avail-stored-data, mss-upd-	M
	aperiodic, mss-acc-manager-initiated, mss-acc-agent-	
	initiated.	
Metric-Structure-Small	See IEEE Std 11073-20601-2014.	NR
Measurement-Status	See IEEE Std 11073-20601-2014.	NR
Metric-Id	See IEEE Std 11073-20601-2014.	NR
Metric-Id-List	See IEEE Std 11073-20601-2014.	NR
Metric-Id-Partition	See IEEE Std 11073-20601-2014.	NR
Unit-Code	See IEEE Std 11073-20601-2014.	NR
Attribute-Value-Map	See IEEE Std 11073-20601-2014.	C
Source-Handle-Reference	Please see Table 23.	NR
Label-String	See IEEE Std 11073-20601-2014.	NR
Unit-LabelString	See IEEE Std 11073-20601-2014.	NR
Absolute-Time-Stamp	See IEEE Std 11073-20601-2014.	С
Base-Offset-Time-Stamp	See IEEE Std 11073-20601-2014.	R
Relative-Time-Stamp	See IEEE Std 11073-20601-2014.	С
HiRes-Time-Stamp	See IEEE Std 11073-20601-2014.	С
Measure-Active-Period	See IEEE Std 11073-20601-2014.	NR
Enum-Observed-Value-Simple-OID	See IEEE Std 11073-20601-2014.	NR
Enum-Observed-Value-Simple-Bit-Str	Please see Table 23.	M
Enum-Observed-Value-Basic-Bit-Str	See IEEE Std 11073-20601-2014.	С
Enum-Observed-Value-Simple-Str	See IEEE Std 11073-20601-2014.	С
Enum-Observed-Value	See IEEE Std 11073-20601-2014.	С
Enum-Observed-Value-Partition	See IEEE Std 11073-20601-2014.	С

NOTE 1—See IEEE Std 11073-20601-2014 for information on whether an attribute is static or dynamic.

NOTE 2—See 6.3 for a description of the qualifiers.

The observed value reported in this object is the general device status.

Since these are essentially event flags, the Unit-Code attribute is not appropriate for this object. Similarly, Source-Handle-Reference is inappropriate, as this object monitors the status of the equipment.

Health informatics—Personal health device communication—Part 10419: Device Specialization—Insulin Pump

The explicit expression of the existence of annunciations is realized by the setting of the appropriate bit in the Enum-Observed-Value-Simple-Bit-Str attribute, as defined in Table 23. If a manager supports the interpretation of this object, it shall be able to interpret the entire set of presented conditions. Anytime the status changes for any monitored condition, the agent shall report on all monitored conditions.

The detection of the condition change could take some time. As time stamp shall be used the correct start and stop time without the delay of the detection duration.

If an acceptable, existing bit is not available, device-status-undetermined shall be. A manager shall interpret these bits only within the context of this attribute and only within this device specialization, as other specializations may use corresponding terms for different purposes.

Table 23 — Mapping of PHD DM status to object Bit-Str attribute

PHD DM status condition	PHDDMStat mnemonic
Agent reports that an undetermined or not supported	device-status-undetermined
condition occurred.	
Agent reports that a reset has occurred.	device-status-reset
Agent reports that a general fault occurred.	device-status-error
Agent reports that a mechanical fault occurred.	device-status-error-mechanical
Agent reports that an electronic fault occurred.	device-status-error-electronic
Agent reports that a software error occurred.	device-status-error-software
Agent reports that a battery fault occurred.	device-status-error-battery
Agent reports that a general service is required.	device-status-service
Agent reports that a time synchronization is required.	device-status-service-time-sync-required
Agent reports that a calibration is required.	device-status-service-calibration-required
Agent reports that a component replenishment is required.	device-status-service-replenishment-required
Agent reports that battery power is low.	device-status-battery-low
Agent reports that battery is depleted.	device-status-battery-depleted
Agent reports that battery has been replaced.	device-status-battery-replaced
Agent reports that battery is interrupted.	device-status-battery-interrupted

NOTE 1—The bits in the table above are defined to: 0 = False and 1 = True.

NOTE 2—The specific bit mappings of PHDDMStat are defined in Annex B.

NOTE 3—All bits not defined in Table 23 or Annex B are reserved for future use.

NOTE 4—An agent is not required to implement all the features specified in Table 23.

6.9.4 Insulin pump status

The insulin pump status object allows insulin pump specific information to be recorded in order to track important events for the user and troubleshooting information for manufacturers. Insulin pump notifications typically indicate conditions that inhibit proper insulin pump operation until the cause of the condition is eliminated.

Health informatics—Personal health device communication—Part 10419: Device Specialization—Insulin Pump

Table 24 defines the attributes for the object that represents the insulin pump status. The insulin pump status enumeration object may be supported by an insulin pump agent in extended configurations.

If this object is to be implemented, then the nomenclature code to identify the enumeration object class is MDC_MOC_VMO_METRIC_ENUM. Refer to Table 24 for the set of attributes of this object.

Table 24—Insulin pump status enumeration object attributes

Attribute name	Extended configuration	Qualifier
Handle	See IEEE Std 11073-20601-2014.	M
Type	MDC_PART_PHD_DM	M
	MDC_INS_PUMP_DEV_STAT	
Supplemental-Types	See IEEE Std 11073-20601-2014.	NR
Metric-Spec-Small	mss-avail-intermittent mss-avail-stored-data mss-upd-	M
	aperiodic mss-msmt-aperiodic mss-acc-agent-initiated	
Metric-Structure-Small	See IEEE Std 11073-20601-2014.	NR
Measurement-Status	See IEEE Std 11073-20601-2014.	NR
Metric-Id	See IEEE Std 11073-20601-2014.	NR
Metric-Id-List	See IEEE Std 11073-20601-2014.	NR
Metric-Id-Partition	See IEEE Std 11073-20601-2014.	NR
Unit-Code	Please see following text.	NR
Attribute-Value-Map	See IEEE Std 11073-20601-2014.	C
Source-Handle-Reference	Please see following text.	NR
Label-String	See IEEE Std 11073-20601-2014.	O
Unit-LabelString	See IEEE Std 11073-20601-2014.	O
Absolute-Time-Stamp	See IEEE Std 11073-20601-2014.	C
Base-Offset-Time-Stamp	See IEEE Std 11073-20601-2014.	R
Relative-Time-Stamp	See IEEE Std 11073-20601-2014.	C
HiRes-Time-Stamp	See IEEE Std 11073-20601-2014.	С
Measure-Active-Period	See IEEE Std 11073-20601-2014.	NR
Enum-Observed-Value-Simple-OID	See IEEE Std 11073-20601-2014.	C
Enum-Observed-Value-Simple-Bit-Str	Please see following text.	M
Enum-Observed-Value-Basic-Bit-Str	See IEEE Std 11073-20601-2014.	NR
Enum-Observed-Value-Simple-Str	See IEEE Std 11073-20601-2014.	NR
Enum-Observed-Value	Please see following text.	NR
Enum-Observed-Value-Partition	See IEEE Std 11073-20601-2014.	O

NOTE 1—See IEEE Std 11073-20601-2014 for information on whether an attribute is static or dynamic.

NOTE 2—See 6.3 for a description of the qualifiers.

The observed value reported in this object is the insulin pump status.

Because these are essentially flags, the Unit-Code attribute is not appropriate for this object. Similarly, Source-Handle-Reference is inappropriate, as this object monitors the status of the equipment.

The explicit expression of the existence of annunciations is realized by the setting of the appropriate bit in the Enum-Observed-Value-Simple-Bit-Str attribute, as defined in Table 25. The Enum-Observed-Value attribute should not be used, as it unnecessarily complicates the modeling of the object. If a manager supports the interpretation of this object, it shall be able to interpret the entire set of presented conditions, defined in the InsPumpStat entity. An agent shall report each time there is a change in status in a condition. An agent may implement any subset of these same conditions. Note that a manager shall interpret these bits only within the context of this attribute and only within this device specialization as other specializations may use corresponding terms for different purposes.

Health informatics—Personal health device communication—Part 10419: Device Specialization—Insulin Pump

Table 25 defines the nomenclature terms for the insulin pump status.

Table 25 — Mapping of insulin pump status to object Bit-Str attribute

Insulin Pump condition	InsPumpStat mnemonic	
Agent reports the air pressure is outside normal operating range.	air-pressure-out-of-range	
Agent reports that active bolus canceled.	bolus-canceled	
Agent detects delivery has reached a user-defined high threshold for	daliyary may	
delivery (basal, bolus, total).	delivery-max	
Agent reports that infusion set (including tubing and/or cannula) is not	infusion-set-detached	
attached to the body.	illusion-set-detached	
Agent reports that the physical connection between infusion set		
(including tubing and/or cannula) and the insulin delivery device is	infusion-set-incomplete	
incomplete.		
Agent reports an occlusion occurred and clogging of infusion set.	occlusion-detected	
Agent reports an insufficient power to charge the device.	power-insufficient	
Agent reports a priming issue (e.g., after replacement of reservoir	priming-issue	
and/or infusion set).	prinning-issue	
Agent reports that the reservoir is extremely low or empty and cannot	reservoir-empty	
continue insulin delivery.	reservoir-empty	
Agent reports an error related to the replacement or functioning of the		
reservoir occurred (e.g., exceeded time out of the device, not properly	reservoir-issue	
attached or has been dislodged) and cannot continue insulin delivery.		
Agent reports that the reservoir fill level reached a manufacturer	reservoir-low	
defined low threshold.	reservon-tow	
Agent reports that the reservoir was attached.	reservoir-attached	
Agent reports that temporary basal rate canceled (e.g., device changed	temp-basal-canceled	
from operational to standby mode)	temp-basar-canceled	
Agent reports that temporary basal rate expired.	temp-basal-expired	
Agent reports the temperature is outside normal operating range.	temperature-out-of-range	

NOTE 1—The bits in the table above are defined to: 0 = False and 1 = True.

NOTE 2—The specific bit mappings of InsPumpStat are defined in Annex B.

NOTE 3—All bits not defined in Table 25 or Annex B are reserved for future use.

NOTE 4—An agent is not required to implement all the features specified in Table 25.

6.10 PM-store objects

6.10.1 General

In the context of personal health devices, insulin pump agent devices are portable or mobile devices. As stated in 5.1, insulin pump agents may be used to collect measurements or observations at a time when out of the network and agent/manager associations cannot be established. It is also common that a given set of measurements or observations made by insulin pump agents may need to be uploaded to more than one manager, for example, in the home and at a medical facility.

To support dual usage, an agent may provide two or more configurations. One configuration may use a temporary measurement storage model that uploads the most recent data immediately on association (agent initiated) with little user intervention, such as might be used by a typical home user that uploads measurements frequently to a personal computer or a mobile device such as a cell phone. Another

Health informatics—Personal health device communication—Part 10419: Device Specialization—Insulin Pump

configuration may use a long-term measurement storage model that uploads data at the request of the manager, such as might be used by the patient's physician or other HCP.

The long-term storage model is realized using PM-stores. Any configuration that does not include a PM-store object utilizes agent-initiated event reports to transmit the observations. The use of temporarily stored data as defined in IEEE Std 11073-20601-2014 is most useful for small numbers of measurements and is subject to automatic deletion upon upload.

Alternatively, in the case where a large number of measurements may be stored or if automatic deletion is to be avoided, a PM-store configuration should be used. Any configuration with a PM-store for persistent storage shall enable access to the PM-store transmissions. As a result, this standard describes a mechanism using PM-store to hold measurements for longer durations. The data held in PM-store objects may be deleted by user actions via the manager or user interface on the device, and the capacity is limited only by the agent's data storage capabilities.

6.10.2 Persistent store model

The PM-store model defined by this standard utilizes one optional PM-store object, the metric PM-store object for the persistent storage of metric data (see Figure 5).



Figure 5—Insulin pump—persistent store model

Note that the PM-store objects are not part of standard configurations defined in this standard.

Following the guides provided in this standard should enable an implementer to store and retrieve the data within this model, but the specifics for determining the specific nature of the data layout and the subsequent visualization, mining, or other managing of the retrieved data is outside the scope of this standard.

The PM-store model defined by this standard utilizes one or more PM-segments for the data of each object to be persistently stored (see Figure 6 for example).

Health informatics—Personal health device communication—Part 10419: Device Specialization—Insulin Pump

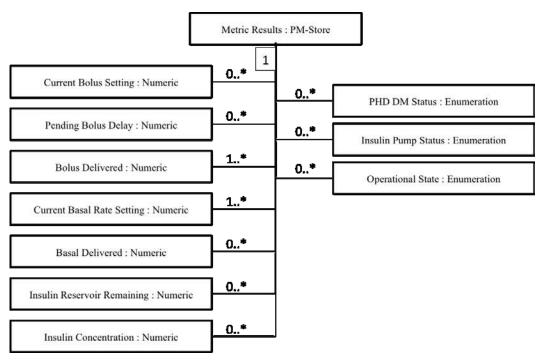


Figure 6—Insulin pump—metric persistent store model

For an insulin pump, the metric PM-store model (see Figure 6) consists of segments holding metric and enumeration readings. The segments holding readings from the objects for bolus delivered and current basal rate setting shall be supported if the metric PM-store is implemented. The other segments are optional and hold observations from the objects that are instantiated.

Note that the metric PM-store object is not part of standard configurations defined in this standard.

Each entry in a segment shall include one of the time formats in the segm-entry-header so a manager can correlate entries across the different segments. If a particular object is not supported, the corresponding segment is not required to exist. Each segment has a cardinality of zero-to-many or one-to-many, as PM-segments are required to contain data from a contiguous period of time (see IEEE Std 11073-20601-2014). Therefore changing time and/or date on the agent typically results in the creation of new segment instances for the supported measurement or observation objects. Furthermore, an insulin pump agent may subdivide data from one contiguous period of time into several segments for further clustering of data (e.g., one segment per day or for an uninterrupted time span of the insulin pump being in operating mode). If a particular segment resulting from such time/date changes or clustering does not contain any entries, it is not required to exist.

Health informatics—Personal health device communication—Part 10419: Device Specialization—Insulin Pump

6.10.3 Metric results PM-store object attributes

Table 26 summarizes the attributes of the metric results PM-store object. The nomenclature code to identify the PM-store object is MDC_MOC_VMO_PMSTORE.

Table 26 — Metric results PM-store object attributes

A 44	Extended configuration		
Attribute name	Value	Qualifier	
Handle	See IEEE Std 11073-20601-2014.	M	
PM-Store-Capab	See IEEE Std 11073-20601-2014.	M	
Store-Sample-Algorithm	See IEEE Std 11073-20601-2014.	M	
Store-Capacity-Count	See IEEE Std 11073-20601-2014.	M	
Store-Usage-Count	See IEEE Std 11073-20601-2014.	M	
Operational-State	See IEEE Std 11073-20601-2014.	M	
PM-Store-Label	See IEEE Std 11073-20601-2014.	О	
Sample-Period	See IEEE Std 11073-20601-2014.	NR	
Number-Of-Segments	See IEEE Std 11073-20601-2014.	M	
Clear-Timeout	See IEEE Std 11073-20601-2014.	M	

The PM-Store-Capab attribute shall set the following bits as indicated:

— pmsc-var-no-of-segm:

If the agent creates new segments either due to storing data of multiple sessions or due to time changes as described in subclause 8.12.2.2, "Comparable time," of IEEE Std 11073-20601-2014, then pmsc-var-no-of-segm shall be set.

— pmsc-epi-seg-entries:

The pmsc-epi-seg-entries bit shall be set.

— pmsc-peri-seg-entries:

The pmsc-peri-seg-entries bit shall not be set.

The remaining bits of the PM-Store-Capab attribute are agent specific and shall be set appropriately.

NOTE 1—See IEEE Std 11073-20601-2014 for information on whether an attribute is static or dynamic.

NOTE 2—See 6.3 for a description of the qualifiers.

6.10.4 PM-store object methods

Table 27 defines the methods of the PM-store objects.

Table 27—PM-store object methods

Service	Subservice type name	Mode	Subservice type (action-type)	Parameters (action-info-args)	Results (action-info-args)
ACTION	Clear- Segments	Confirmed	MDC_ACT_SEG_ CLR	SegmSelection	
	Get-Segment- Info	Confirmed	MDC_ACT_SEG_ GET_INFO	SegmSelection	SegmentInfoList
ACTION	Get-Segment- Id-List	Confirmed	MDC_ACT_SEG_ GET_ID_LIST	(empty)	SegmIdList
	Trig-Segment- Data-Xfer	Confirmed	MDC_ACT_SEG_ TRIG_XFER	TrigSegmDataXferReq	TrigSegmDataXfer Rsp

Health informatics—Personal health device communication—Part 10419: Device Specialization—Insulin Pump

6.10.4.1 Clear-Segments

This method allows the manager to delete all data entries stored in a PM-segment object. The agent shall support the Clear-Segments method by setting the pmsc-clear-segm-by-all-sup bit for the PM-Store-Capab attribute. Deletion of PM-segments is not guaranteed by this method. See IEEE Std 11073-20601-2014 for information on how the agent shall reply in case it decides to protect certain segments from deletion.

— Get-Segment-Info

This method allows the manager to retrieve the PM-segment attributes.

— Get-Segment-Id-List

This method allows the manager to retrieve a list of the instance numbers of all the PM-segments of a PM-store.

— Trig-Segment-Data-Xfer

This method allows the manager to initiate the transfer of the data entries stored in the PM-segment object.

Refer to IEEE Std 11073-20601-2014 for details.

6.10.5 PM-store object events

Table 28 defines the events sent by the PM-store objects.

Table 28 — PM-store object events

Service	Subservice type name	Mode	Subservice type (event-type)	Parameters (event-info)	Results (event-reply-info)
EVENT REPORT	Segment- Data-Event	Confirmed	MDC_NOTI_SEGMENT _DATA	SegmentDataEvent	SegmentDataResult

— Segment-Data-Event

This event allows the agent to send the data entries stored in the PM-segment object. This event is triggered by the manager using the Trig-Segment-Data-Xfer action. Refer to IEEE Std 11073-20601-2014 for details.

6.10.6 PM-store object services

6.10.6.1 GET service

The GET service shall be provided by an agent implementing PM-store objects. This service shall be available only while the agent is in the Operating state. Refer to IEEE Std 11073-20601-2014 for details.

6.10.6.2 SET service

There are currently no SET services defined for PM-store objects in this standard.

Health informatics—Personal health device communication—Part 10419: Device Specialization—Insulin Pump

6.10.7 PM-segment objects

Table 29 defines the attributes of the periodic session PM-segment object contained in the periodic PM-store object managing the stored measurements or observations. The nomenclature code to identify the PM-segment class is MDC_MOC_PM_SEGMENT.

Table 29—Common PM-segment object attributes

Attribute name	Extended configuration			
Attribute name	Value	Qualifier		
Instance-Number	See IEEE Std 11073-20601-2014.	M		
PM-Segment-Entry-Map	See IEEE Std 11073-20601-2014.	M		
PM-Seg-Person-Id	See IEEE Std 11073-20601-2014.	С		
Operational-State	See IEEE Std 11073-20601-2014.	M		
Sample-Period	See IEEE Std 11073-20601-2014.	С		
Segment-Label	See IEEE Std 11073-20601-2014.	0		
Segment-Start-Abs-Time	See IEEE Std 11073-20601-2014.	С		
Segment-End-Abs-Time	See IEEE Std 11073-20601-2014.	С		
Date-and-Time-Adjustment	See IEEE Std 11073-20601-2014.	С		
Segment-Start-BO-Time	See IEEE Std 11073-20601-2014.	С		
Segment-End-BO-Time	See IEEE Std 11073-20601-2014.	С		
Segment-Usage-Count	See IEEE Std 11073-20601-2014.	M		
Segment-Statistics	See IEEE Std 11073-20601-2014.	0		
Fixed-Segment-Data	Segment data transferred as an array of entries in a format as specified in the PM-Segment-Entry-Map attribute.	M		
Confirm-Timeout	See IEEE Std 11073-20601-2014.	0		
Transfer-Timeout	See IEEE Std 11073-20601-2014.	M		

NOTE 1—See IEEE Std 11073-20601-2014 for information on whether an attribute is static or dynamic.

NOTE 2—See 6.3 for a description of the qualifiers.

The Fixed-Segment-Data attribute serves as the container of the stored measurements or observations. When the Fixed-Segment-Data attribute is transmitted, all entries in the event report are formatted according to the PM-Segment-Entry-Map. Each entry contains an optional header and one or more elements. Each element holds data from one or more metric measurements.

6.11 Schedule-store objects

6.11.1 General

In the context of personal health devices, insulin pump agent devices can store schedules at a time when agent/manager are associated or when out of the network and agent/manager associations cannot be established. It is also common that a given set of schedules made by insulin pump agents could be uploaded to more than one manager, for example, in the home and at a medical facility.

The data held in schedule-store objects may be deleted by user actions via the manager or user interface on the device, and the capacity is limited only by the agent's data storage capabilities.

Health informatics—Personal health device communication—Part 10419: Device Specialization—Insulin Pump

6.11.2 Schedule-store model

The schedule-store model defined by this standard utilizes three optional schedule-store objects:

- the basal profile settings schedule-store,
- the insulin to carbohydrate profile settings schedule-store, and
- the insulin sensitivity factor profile settings schedule-store for the storage of basal profile settings, insulin to carbohydrate ratio profile settings, and insulin sensitivity factor profile settings, respectively (see Figure 7).

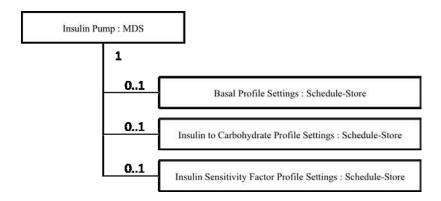


Figure 7—Insulin pump—schedule-store model

Note that the schedule-store objects are not part of standard configurations defined in this standard.

Following the guides provided in this standard should enable an implementer to store and retrieve the data within this model, but the specifics for determining the specific nature of the data layout and the subsequent visualization, mining, or other managing of the retrieved data is outside the scope of this standard.

6.11.3 Basal profile settings schedule-store object attributes

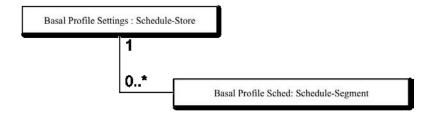


Figure 8—Insulin pump—basal profile settings schedule-store model

For an insulin pump, the basal profile settings schedule-store model (see Figure 8) may consist of zero or more basal profile schedule-segments to account for the collection of multiple basal profile settings.

If the insulin pump device offers the option of programming multiple basal profiles, a unique basal profile identifier is captured in the basal profile sched schedule-segment and can be used to identify the different

Health informatics—Personal health device communication—Part 10419: Device Specialization—Insulin Pump

collections of basal profile settings. The active basal profile sched is captured in the basal profile settings schedule-store.

Within a basal profile sched schedule-segment, at least one entry shall be used to account for a basal rate schedule setting.

Figure 9 is an example of a basal rate profile setting schedule-store.

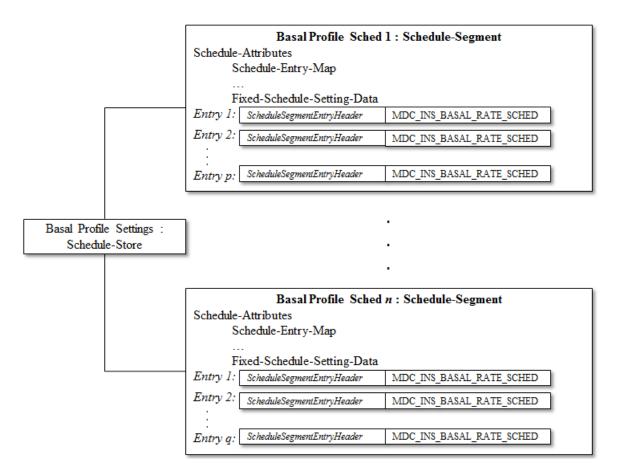


Figure 9—Basal profile settings schedule-store

Note that the basal profile settings schedule-store object is not part of standard configurations defined in this standard.

Table 30 summarizes the attributes of the basal profile settings schedule-store object. The nomenclature code to identify the basal profile settings schedule-store object is MDC_MOC_VMO_SCHEDSTORE.

Health informatics—Personal health device communication—Part 10419: Device Specialization—Insulin Pump

Table 30—Basal profile settings schedule-store object attributes

Attribute name	Extended configuration		
Attribute name	Value	Qualifier	
Schedule-Handle	See Annex E.	M	
Active-Schedule-Segment-Instance-Number	See Annex E.	M	
Updated-Schedule-Segment-Instance-Number-List	See Annex E.	M	
Schedule-Store-Capab	See Annex E.	M	
Schedule-Store-Capacity-Count	See Annex E.	M	
Schedule-Store-Usage-Count	See Annex E.	M	
Schedule-Store-Operational-Status	See Annex E.	M	
Schedule-Store-Label	See Annex E.	0	
Number-Of-Schedule-Segments	See Annex E.	M	

The basal profile settings Schedule-Store-Capab attribute shall set the following bits as indicated:

— schedsc-var-no-of-segm:

If the agent creates new schedule segments either due to storing data of multiple schedules or due to time changes as described in "Comparable time" subclause 8.12.2.2 of IEEE Std 11073-20601-2014, then schedsc-var-no-of-segm shall be set.

— schedsc-epi-seg-entries:

The schedsc-epi-seg-entries bit shall be set.

— schedsc-peri-seg-entries:

The schedsc-peri-seg-entries bit shall not be set.

The remaining bits of the Schedule-Store-Capab attribute are agent specific and shall be set appropriately.

NOTE 1—See Annex E for information on whether an attribute is static or dynamic.

NOTE 2—See 6.3 for a description of the qualifiers.

6.11.4 Insulin to carbohydrate profile settings schedule-store object attributes

For an insulin pump, the insulin to carbohydrate (I:CHO) profile settings schedule-store model (see Figure 10) may consist of zero or more I:CHO profile sched schedule segments to account for the collection of multiple I:CHO profile settings.

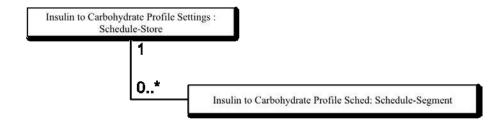


Figure 10—Insulin to carbohydrate profile settings schedule-store model

Health informatics—Personal health device communication—Part 10419: Device Specialization—Insulin Pump

If the insulin pump device offers the option of programming multiple I:CHO profiles, the unique I:CHO profile identifier is captured in the I:CHO profile sched schedule segment and can be used to identify the different I:CHO profile settings. The active I:CHO profile sched is captured in the I:CHO profile sched schedule-store.

Within an I:CHO profile sched schedule segment, at least one entry shall be used to account for an I:CHO schedule setting.

Figure 11 is an example of an I:CHO profile setting schedule-store.

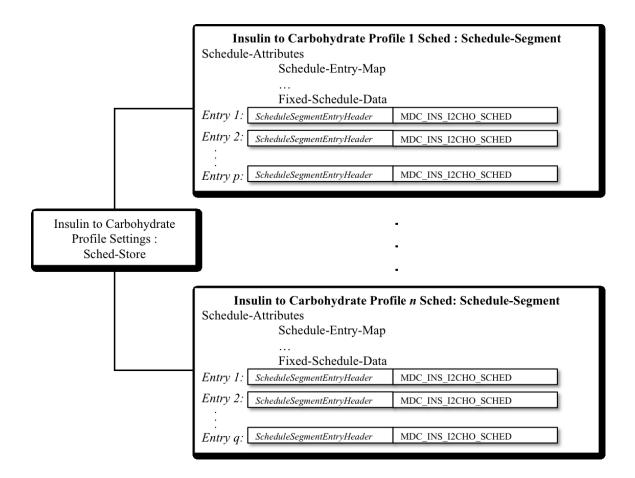


Figure 11—Insulin to carbohydrate profile settings schedule-store

Note that the I:CHO profile settings schedule-store object is not part of standard configurations defined in this standard.

Health informatics—Personal health device communication—Part 10419: Device Specialization—Insulin Pump

Table 31 summarizes the attributes of the I:CHO profiles schedule-store object. The nomenclature code to identify the I:CHO profile settings schedule-store object is MDC_MOC_VMO_SCHEDSTORE.

Table 31 —Insulin to carbohydrate profile settings schedule-store object attributes

Attribute name	Extended configuration		
Attribute name	Value	Qualifier	
Schedule-Handle	See Annex E.	M	
Active-Schedule-Segment-Instance-Number	See Annex E.	M	
Updated-Schedule-Segment-Instance-Number-List	See Annex E.	M	
Schedule-Store-Capab	See Annex E.	M	
Schedule-Store-Capacity-Count	See Annex E.	M	
Schedule-Store-Usage-Count	See Annex E.	M	
Schedule-Store-Operational-Status	See Annex E.	M	
Schedule-Store-Label	See Annex E.	0	
Number-Of-Schedule-Segments	See Annex E.	M	

The I:CHO profile settings Schedule-Store-Capab attribute shall set the following bits as indicated:

— schedsc-var-no-of-segm:

If the agent creates new schedule-segments either due to storing data of multiple schedules or due to time changes as described in the "Comparable time" subclause 8.12.2.2 of IEEE Std 11073-20601-2014, then schedsc-var-no-of-segm shall be set.

— schedsc-epi-seg-entries:

The schedsc-epi-seg-entries bit shall be set.

— schedsc-peri-seg-entries:

The schedsc-peri-seg-entries bit shall not be set.

The remaining bits of the Schedule-Store-Capab attribute are agent specific and shall be set appropriately.

NOTE 1—See Annex E for information on whether an attribute is static or dynamic.

NOTE 2—See 6.3 for a description of the qualifiers.

6.11.5 Insulin sensitivity factor profiles schedule-store object attributes

For an insulin pump, the insulin sensitivity factor (ISF) profile settings schedule-store model (see Figure 12) may consist of zero or more ISF profile sched schedule-segments to account for the data collection of multiple ISF profile settings.

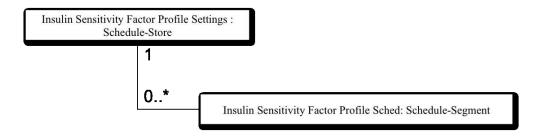


Figure 12—Insulin sensitivity factor profile settings schedule-store model

Health informatics—Personal health device communication—Part 10419: Device Specialization—Insulin Pump

If the insulin pump device offers the option of programming multiple ISF profiles, the unique ISF profile identifier is captured in the ISF profile sched schedule-segment and can be used to identify the different ISF profile settings. The active ISF profile sched is captured in the ISF profile sched schedule-store.

Within an ISF profile sched schedule-segment, at least one entry shall be used to account for an ISF schedule setting.

Figure 13 is an example of an ISF profile setting schedule-store.

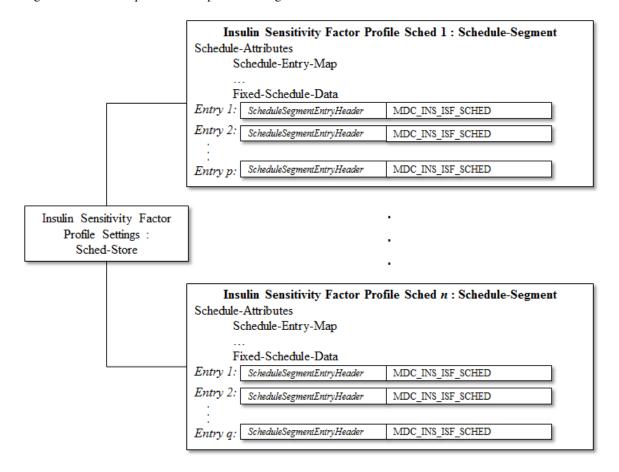


Figure 13—ISF profile settings schedule-store

Note that the ISF profile settings schedule-store object is not part of standard configurations defined in this standard.

Table 32 summarizes the attributes of the ISF profile settings schedule-store object. The nomenclature code to identify the ISF profile settings schedule-store object is MDC_MOC_VMO_SCHEDSTORE.

Health informatics—Personal health device communication—Part 10419: Device Specialization—Insulin Pump

Table 32—Insulin sensitivity factor profile settings schedule-store object attributes

Attribute name	Extended configuration		
Attribute name	Value	Qualifiers	
Schedule-Handle	See Annex E.	M	
Active-Schedule-Segment-Instance-Number	See Annex E.	M	
Updated-Schedule-Segment-Instance-Number-List	See Annex E.	M	
Schedule-Store-Capab	See Annex E.	M	
Schedule-Store-Capacity-Count	See Annex E.	M	
Schedule-Store-Usage-Count	See Annex E.	M	
Schedule-Store-Operational-Status	See Annex E.	M	
Schedule-Store-Label	See Annex E.	0	
Number-Of-Schedule-Segments	See Annex E.	M	

The ISF profile settings Schedule-Store-Capab attribute shall set the following bits as indicated:

— schedsc-var-no-of-segm:

If the agent creates new schedule-segments either due to storing data of multiple schedules or due to time changes as described in the "Comparable time" subclause 8.12.2.2 of IEEE Std 11073-20601-2014, then schedsc-var-no-of-segm shall be set.

— schedsc-epi-seg-entries:

The schedsc-epi-seg-entries bit shall be set.

— schedsc-peri-seg-entries:

The schedsc-peri-seg-entries bit shall not be set.

The remaining bits of the Schedule-Store-Capab attribute are agent specific and shall be set appropriately.

NOTE 1—See Annex E for information on whether an attribute is static or dynamic.

NOTE 2—See 6.3 for a description of the qualifiers.

6.11.6 Schedule-store object methods

This subclause applies to basal profile settings, insulin to carbohydrate profile settings, and insulin sensitivity factor profile settings schedule-store objects.

Table 33 defines the methods of a schedule-store object.

Table 33 — Common schedule-store object methods

Service	Subservice type name	Mode	Subservice type (action-type)	Parameters (action-info-args)	Results (action-info-args)
ACTION	Get- Schedule- Segment- Info	Confirmed	MDC_ACT_SCHED_SEG _GET_INFO	SchedSegmSelection	SchedSegmentInfoList
	Get- Schedule- Segment-Id- List	Confirmed	MDC_ACT_SCHED_SEG_GET _ID_LIST	(empty)	SchedSegmIdList
	Trig- Schedule- Segment- Data-Xfer	Confirmed	MDC_ACT_SCHED_SEG _TRIG_XFER	TrigSchedSegmData- XferReq	TrigSchedSegmDataXferRsp

Health informatics—Personal health device communication—Part 10419: Device Specialization—Insulin Pump

— Get-Schedule-Segment-Info

This method allows the manager to retrieve schedule-segment attributes.

Get-Schedule-Segment-Id-List

This method allows the manager to retrieve a list of the instance numbers of all the schedule-segments of a schedule-store.

— Trig-Schedule-Segment-Data-Xfer

This method allows the manager to initiate the transfer of the data entries stored in a schedule-segment object.

Refer to Annex E for details.

6.11.7 Schedule-store object events

This subclause applies to basal profile settings, insulin to carbohydrate profile settings, and insulin sensitivity factor profile settings schedule-store objects.

Table 34 defines the events sent by schedule-store objects.

Table 34 —Insulin pump schedule-store object events

Service	Subservice type	Mode	Subservice type	Parameters	Results
	name		(event-type)	(event-info)	(event-reply-info)
EVENT	Schedule-Segment-	Confirmed	MDC_NOTI_SCH	ScheduleSegmentD	ScheduleSegmentD
REPORT	Data-Event		ED_SEGMENT_D	ataEvent	ataResult
			ATA		

— Schedule-Segment-Data-Event

This event allows the agent to send the data entries stored in a schedule-segment object. This event is triggered by the manager using the Trig-Schedule-Segment-Data-Xfer action. Refer to E.1.5.

6.11.8 Schedule-store object services

This subclause applies to basal profile settings, insulin to carbohydrate profile settings, and insulin sensitivity factor profile settings schedule-store objects.

6.11.8.1 GET service

The GET service shall be provided by an agent implementing schedule-store objects. This service shall be available only while the agent is in the Operating state. Refer to E.1.4 for details.

6.11.8.2 SET service

There are currently no SET services defined for schedule-store objects in this standard.

6.11.9 Schedule-segment objects

This subclause applies to basal profile sched, insulin to carbohydrate profile sched, and insulin sensitivity factor profile sched schedule-segment objects.

Health informatics—Personal health device communication—Part 10419: Device Specialization—Insulin Pump

Each of the schedule-store objects contains at least one or more corresponding schedule-segment object.

Table 35 defines the attributes schedule-segment object contained in the schedule-store object managing the stored schedules. The nomenclature code to identify the schedule-segment class is MDC_MOC_SCHEDULE_SEGMENT.

Table 35—Schedule-segment object attributes

A 44-214	Extended configuration		
Attribute name	Value	Qualifier	
Schedule-Segment-Instance-Number	See Annex E.	M	
Schedule-Segment-Entry-Map	See Annex E.	M	
Schedule-Segment-Period	See Annex E.	M	
Schedule-Segment-Entry-Interval	See Annex E.	О	
Schedule-Segment-Person-Id	See Annex E.	О	
Schedule-Segment-Entry-Count	See Annex E.	О	
Schedule-Segment-Label	See Annex E.	О	
Schedule-Segment-LastUpdated-Abs-Time	See Annex E.	С	
Schedule-Segment-LastUpdated-HiRes-Time	See Annex E.	С	
Schedule-Segment-LastUpdated-BO-Time	See Annex E.	R	
Schedule-Segment-Reference-Abs-Time	See Annex E.	C	
Schedule-Segment-Reference-BO-Time	See Annex E.	R	
Schedule-Segment-Start-Abs-Time	See Annex E.	C	
Schedule-Segment-End-Abs-Time	See Annex E.	С	
Schedule-Segment-Start-BO-Time	See Annex E.	R	
Schedule-Segment-End- BO-Time	See Annex E.	R	
Fixed-Schedule-Segment-Data	Schedule-segment data transferred as an array of entries in a format as specified in the Schedule-Segment- Entry-Map attribute.	M	
Schedule-Segment-Confirm-Timeout	See Annex E.	О	
Schedule_Transfer-Timeout	See Annex E.	M	

NOTE 1—See Annex E for information on whether an attribute is static or dynamic.

NOTE 2—See 6.3 for a description of the qualifiers.

The Fixed-Schedule-Segment-Data attribute serves as the container of the schedules. When the Fixed-Schedule-Segment-Data attribute is transmitted, all entries in the event report are formatted according to the Schedule-Segment-Entry-Map. Each entry stores a single sample point, which may consist of a set of attributes.

6.12 Scanner objects

Scanner objects are not required by this standard.

6.13 Class extension objects

In this standard, no class extension objects are defined with respect to IEEE Std 11073-20601-2014.

Health informatics—Personal health device communication—Part 10419: Device Specialization—Insulin Pump

6.14 Insulin pump information model extensibility rules

The insulin pump domain information model of this standard may be extended by including elements defined in IEEE Std 11073-20601-2014 as well as vendor-specific elements. Any object or attribute extensions implemented should follow the guidelines of this standard as closely as possible.

An insulin pump agent having a configuration with extensions beyond the standard configuration, as specified in this standard, shall use a configuration ID in the range of IDs reserved for extended configurations (see IEEE Std 11073-20601-2014).

7. Insulin pump service model

7.1 General

The service model defines the conceptual mechanisms for data exchange services. These services are mapped to messages that are exchanged between the agent and the manager. Protocol messages within the ISO/IEEE 11073 series of standards are defined in ASN.1. See IEEE Std 11073-20601-2014 for a detailed description of the personal health device service model. Subclauses 7.2 and 7.3 define the specifics of object access and event reporting services for an insulin pump agent according to this standard.

7.2 Object access services

The object access services of IEEE Std 11073-20601-2014 are used to access the objects defined in the domain information model of the insulin pump.

The following generic object access services are supported by an insulin pump agent according to this standard:

- GET service: used by the manager to retrieve the values of the agent MDS object or PM-store attributes or schedule-store attributes. The list of insulin pump attributes is given in 6.6.4.1 for the MDS object, and 6.10.3 for the PM-store objects, and 6.11.3, 6.11.4, and 6.11.5 for the schedule-store objects.
- SET service: used by the manager to set the values of the agent object attributes. There are no settable attributes defined for an insulin pump agent according to this standard.
- Event report service: used by the agent to send configuration reports and measurement data to the manager. The list of event reports for the insulin pump device specialization is given for the MDS object in 6.6.3, for the PM-store object in 6.10.5, and for the schedule-store object in 6.11.7.
- Action service: used by the manager to invoke actions (or methods) supported by the agent. The list of MDS objects actions for setting the time is given in 6.6.2. The list of PM-store object actions is given in 6.10.4. The list of schedule-store object actions is given in 6.11.6.

Table 36 summarizes the object access services described in this standard.

Table 36—Insulin pump object access services

Service	Subservice type name	Mode	Subservice type	Parameters	Result	Remarks
	<na></na>	<implied Confirmed></implied 	<na></na>	GetArgumentSi mple = (obj-handle = 0), attribute-id- list <optional></optional>	GetResultSimple = (obj-handle = 0), attribute-list	Allows the manager to retrieve the value of attributes of the MDS object in the agent.
GET	<na></na>	<implied confirmed=""></implied>	<na></na>	GetArgumentSi mple = (obj-handle = handle of PM- store object), attribute-id-list <optional></optional>	GetResultSimple = (obj-handle = handle of PM- store object), attribute-list	Allows the manager to retrieve the values of attributes of a PM-store object in the agent.
	<na></na>	<implied confirmed=""></implied>	<na></na>	GetArgumentSi mple = (obj-handle = handle of Schedule-store object), attribute-id-list <optional></optional>	GetResultSimple = (obj-handle = handle of Schedule-store object), attribute- list	Allows the manager to retrieve the values of attributes of a Schedule-store object in the agent.
	MDS- Configuratio n-Event	Confirmed	MDC_NOTI _CONFIG	ConfigReport	ConfigReportRsp	Configuration Report to inform manager of the configuration of the agent.
	MDS- Dynamic- Data- Update-Var	Confirmed	MDC_NOTI _SCAN_REP ORT_VAR	ScanReportInfo Var	_	Data Report to provide dynamic data to manager for some or all of the agent's objects in variable format.
	MDS- Dynamic- Data- Update- Fixed	Confirmed	MDC_NOTI _SCAN_REP ORT_FIXED	ScanReportInfo Fixed	_	Data Report to provide dynamic data to manager for some or all of the agent's objects in fixed format.
EVENT	MDS- Dynamic- Data- Update-MP- Var	Confirmed	MDC_NOTI _SCAN_REP ORT_MP_V AR	ScanReportInfo MPVar	_	This is the same as MDS- Dynamic-Data-Update-Var but allows inclusion of data from multiple people.
REPORT	MDS- Dynamic- Data- Update-MP- Fixed	Confirmed	MDC_NOTI _SCAN_REP ORT_MP_FI XED	ScanReportInfo MPFixed	_	This is the same as MDS- Dynamic-Data-Update-Fixed but allows inclusion of data from multiple people.
	Segment- Data-Event	Confirmed	MDC_NOTI _SEGMENT _DATA	SegmentDataEv ent	SegmentDataRes ult	PM-store object event to provide data stored in the Fixed-Segment-Data of a PM-segment from the agent to the manager.
	Schedule- Segment- Data-Event	Confirmed	MDC_NOTI SCHED_SE GMENT_DA TA	ScheduleSegme ntDataEvent	ScheduleSegment DataResult	Schedule-store object event to provide data stored in the Fixed-Schedule-Segment-Data of a Schedule-segment from the agent to the manager.

Health informatics—Personal health device communication—Part 10419: Device Specialization—Insulin Pump

Table 36—Insulin pump object access services (continued)

Service	Subservice type name	Mode	Subservice type	Parameters	Result	Remarks
	Set-Time	Confirmed	MDC_ACT_ SET_TIME	SetTimeInvoke	_	Manager method to invoke the agent to set time in absolute time format to requested value.
	Set-Base- Offset-Time	Confirmed	MDC_ACT_ SET_BO_TI ME	SetBOTime Invoke	_	Manager method to invoke the agent to set time in base offset time format to requested value.
	Clear -Segments	Confirmed	MDC_ACT_ SEG_CLR	SegmSelection	_	Allows the manager to delete data stored in selected PM-segments in the agent.
	Get- Segment- Info	Confirmed	MDC_ACT_ SEG_GET_I NFO	SegmSelection	SegmentInfoList	Allows the manager to retrieve the value of PM-segment attributes of one or more PM-segments in the agent.
ACTION	Get- Segment-Id- List	Confirmed	MDC_ACT_ SEG_GET_I D_LIST	(empty)	SegmIdList	Allows the manager to retrieve a list of the instance numbers of all the PM-segments of a PM-store.
	Get- Schedule- Segment- Info	Confirmed	MDC_ACT_ SCHED_SE G_GET_INF O	SehedSegmSele ction	SchedSegmentInf oList	Allows the manager to retrieve the value of Schedule-segment attributes of one or more Schedule-segments in the agent.
	Get- Schedule- Segment-Id- List	Confirmed	MDC_ACTS CHED_SEG_ GET_ID_LIS T	(empty)	SchedSegmIdList	Allows the manager to retrieve a list of the instance numbers of all the Schedulesegments of a Schedule-store.
	Trig- Segment- Data-Xfer	Confirmed	MDC_ACT_ SEG_TRIG_ XFER	TrigSegmDataX ferReq	TrigSegmData XferRsp	Allows the manager to start the transfer of the Fixed- Segment-Data attribute of a PM-segment in the agent.
	Trig- Schedul- Segment- Data-Xfer	Confirmed	MDC_ACT_ SCHED_SE G_TRIG_XF ER	TrigSchedSegm DataXferReq	TrigSchedSegmD ata XferRsp	Allows the manager to start the transfer of the Fixed- Schedule-Segment-Data attribute of a Schedule- segment in the agent.

7.3 Object access event report services

The event report service (see Table 36) is used by the agent to report its information (e.g., measurements). Event reports in this standard are a property of the MDS object only. The event reports used in this standard are defined in IEEE Std 11073-20601-2014.

The following conditions apply for an insulin pump agent according to this standard:

- MDS event reports shall be used in confirmed mode.
- Agent initiated mode shall be supported for measurement data transmission.
- Persistently stored metric mode may be supported for measurement data transmission.
- Manager initiated mode may be support for measurement data transmission.

An insulin pump agent, which is designed to operate in an environment where data may be collected from multiple people, may use one of the multiple-person event report styles to transmit all the data from each

Health informatics—Personal health device communication—Part 10419: Device Specialization—Insulin Pump

person in a single event. If this functionality is not required, the agent may use the single-person event report styles, which have reduced overhead.

A manager shall support both single-person and multiple-person event reports. An insulin pump agent may support either one or both single-person and multiple-person event reports. The formats for single-person and multiple-person reports are described in IEEE Std 11073-20601-2014.

8. Insulin pump communication model

8.1 Overview

This clause describes the general communication model and procedures of the insulin pump agent as defined in IEEE Std 11073-20601-2014. Therefore, the respective parts of IEEE Std 11073-20601-2014 are not reproduced; rather the specific choices and restrictions with respect to optional elements (e.g., objects, attributes, and actions) and specific extensions (e.g., nomenclature terms) are specified.

For an illustrative overview of the various message transactions during a typical measurement session, see the sequence diagram for the example use case in Annex D and the corresponding protocol data unit (PDU) examples in Annex I.

8.2 Communications characteristics

In this subclause, limits on the size of an application protocol data unit (APDU) transmitted or to be received by an insulin pump agent are defined. Small limits allow for simple implementations in terms of low cost and complexity.

An insulin pump agent implementing only this device specialization shall not transmit any APDU larger than N_{tx} and shall be capable of receiving any APDU up to a size of N_{rx} . For this standard, N_{tx} shall be 5120 octets for implementations supporting persistent metric storage. In the absence of the persistent metric storage capability, an insulin pump agent implementing only this device specialization shall not transmit any APDU larger than N_{tx} and shall be capable of receiving any APDU up to a size of N_{rx} . For this standard, N_{tx} shall be 7168 octets and N_{rx} shall be 224 octets.

For an insulin pump agent implementing functions from other device specializations, an upper bound estimation of the APDU sizes brings the following: An agent shall not transmit any APDU larger than the sum of N_{tx} of all the device specializations implemented and shall be capable of receiving any APDU up to the sum of N_{rx} of all the device specializations implemented. If these numbers are higher than the maximum size determined in IEEE Std 11073-20601-2014, the latter shall be applied.

In case the APDU size limit does not allow for the inclusion of a certain amount of multiple pending measurements at the agent, they shall be sent using multiple event reports. See 8.5.3 for the maximum number of measurements allowed for inclusion in a single event report.

Health informatics—Personal health device communication—Part 10419: Device Specialization—Insulin Pump

8.3 Association procedure

8.3.1 General

Unless otherwise stated, the association procedure for an insulin pump agent and manager according to this standard shall be pursued as specified in IEEE Std 11073-20601-2014.

8.3.2 Agent procedure—association request

In the association request sent by the agent to the manager:

- The version of the association procedure used by the agent shall be set to assoc-version1 (i.e., assoc-version = 0x80000000).
- The DataProtoList structure element of the data protocol identifier shall be set to data-proto-id-20601 (i.e., *data-proto-id* = 0x5079).
- The *data-proto-info* field shall contain a PhdAssociationInformation structure that shall contain the following parameter values:
 - 1) The agent shall support protocol-version3. Support for any other version may be indicated by setting additional bits. When protocols higher than protocol-version3 are used, the agent shall continue to use only features as specified in this standard. When protocols lower than protocol-version3 are used, the agent shall use only features in that protocol.
 - 2) At least the MDERs shall be supported (i.e., encoding-rules = 0x8000).
 - 3) The version of the nomenclature used shall be set to nom-version1 (i.e., *nomenclature-version* = 0x80000000).
 - 4) The field *functional-units* may have the test association bits set but shall not have any other bits set.
 - 5) The field *system-type* shall be set to sys-type-agent (i.e., *system-type* = 0x00800000).
 - 6) The *system-id* field shall be set to the value of the System-Id attribute of the MDS object of the agent. The manager may use this field to determine the identity of the insulin pump with which it is associating and, optionally, to implement a simple access restriction policy.
 - 7) The *dev-config-id* field shall be set to the value of the Dev-Configuration-Id attribute of the MDS object of the agent.
 - 8) If the agent supports only the insulin pump specialization, then the field indicating the data request modes (*data-req-mode-capab*) supported by the insulin pump agent shall be set to *data-req-supp-init-agent*.
 - 9) If the agent supports only the insulin pump specialization, then *data-req-init-manager-count* shall be set to zero, and *data-req-init-agent-count* shall be set to 1.

8.3.3 Manager procedure—association response

In the association response message sent by the manager:

— The *result* field shall be set to an appropriate response from those defined in IEEE Std 11073-20601-2014. For example, if all other conditions of the association protocol are satisfied, *accepted* is returned when the manager recognizes the *dev-config-id* of the agent and *accepted-unknown-config* otherwise.

Health informatics—Personal health device communication—Part 10419: Device Specialization—Insulin Pump

- In the DataProtoList structure element, the data protocol identifier shall be set to data-proto-id-20601 (i.e., data-proto-id = 0x5079).
- The *data-proto-info* field shall be filled in with a PhdAssociationInformation structure that shall contain the following parameter values:
 - 1) The manager following this specialization shall support protocol-version2. The manager may support additional protocol versions and select them if the agent offers them.
 - 2) The manager shall respond with a single selected encoding rule that is supported by both agent and manager. The manager shall support at least the MDERs.
 - 3) The version of the nomenclature used shall be set to nom-version1 (i.e., *nomenclature-version* = 0x80000000).
 - 4) The field *functional-units* shall have all bits reset except for those relating to a test association.
 - 5) The field *system-type* shall be set to sys-type-manager (i.e., *system-type* = 0x80000000).
 - 6) The *system-id* field shall contain the unique system ID of the manager device, which shall be a valid EUI-64 type identifier.
 - 7) The field *dev-config-id* shall be manager-config-response (0).
 - 8) The field *data-req-mode-capab* shall be 0.
 - 9) If the agent supports only the insulin pump specialization, *data-req-initagent-count* shall be 1 and *data-req-init-manager-count* shall be 0.

8.4 Configuring procedure

8.4.1 General

The agent enters the Configuring state if it receives an association response of accepted-unknown-config. In this case, the configuration procedure as specified in IEEE Std 11073-20601-2014 shall be followed. Subclause 8.4.2 specifies the configuration notification and response messages for insulin pump agent with standard configuration ID 1900 (0x076C). Normally, a manager would already know the standard configuration. However, for the purposes of this example, it does not.

8.4.2 Insulin pump—standard configuration

8.4.2.1 Agent procedure

The agent performs the configuration procedure using a "Remote Operation Invoke | Confirmed Event Report" message with an MDC_NOTI_CONFIG event to send its configuration to the manager (see IEEE Std 11073-20601-2014). The ConfigReport structure is used for the *event-info* field (see Table 3). For an insulin pump agent with standard configuration ID 1900 (0x076C), the format and contents of the configuration notification message are as follows:

$0 \times E7$	0x00	APDU CHOICE Type (PrstApdu)
0x00	0×70	CHOICE.length = 112
0x00	0x6e	OCTET STRING.length = 110
0x00	0x02	invoke-id = 2 (start of DataApdu. MDER encoded.)
0x01	0x01	CHOICE(Remote Operation Invoke Confirmed Event Report)

0x00 0x68		CHOICE.length = 104
$0x00 \ 0x00$		obj-handle = 0 (MDS object)
$0x00 \ 0x00 \ 0x00$	0x00	event-time = 0
0x0D 0x1C		event-type = MDC_NOTI_CONFIG
0x00 0x2e		event-info.length = 46 (start of ConfigReport)
0x07 0x6C		config-report-id (Dev-Configuration-Id value)
$0x00 \ 0x01$		config-obj-list.count = 2 Measurement objects will be "announced"
$0x00 \ 0x58$		config-obj-list.length = 88
$0x00 \ 0x06$		obj-class = MDC_MOC_VMO_METRIC_NU
$0x00 \ 0x01$		obj-handle = 1 (\rightarrow 1 st object is bolus delivered)
$0x00 \ 0x04$		attributes.count = 4
$0x00 \ 0x24$		attributes.length = 36
0x09 0x2F		attribute-id = MDC_ATTR_ID_TYPE
$0x00 \ 0x04$		attribute-value.length = 4
0x00 0x80 0x74	0x1C	MDC_PART_PHD DM MDC_INS_BOLUS
0x0A 0x46		attribute-id = MDC_ATTR_METRIC_SPEC_SMALL
$0x00 \ 0x02$		attribute-value.length = 2
$0x40 \ 0x40$		Metric-Spec-Small (mss-avail-stored-data, mss-acc-agent-initiated)
0x09 0x96		attribute-id = MDC_ATTR_UNIT_CODE
$0x00 \ 0x02$		attribute-value.length = 2
$0x15 \ 0x60$		MDC_DIM_INTL_UNIT
0x0A 0x55		attribute-id = MDC_ATTR_ATTRIBUTE_VAL_MAP
$0x00 \ 0x0C$		attribute-value.length = 12
$0x00 \ 0x02$		AttrValMap.count = 2
$0x00 \ 0x08$		AttrValMap.length = 8
0x0A 0x4C 0x00	0x02	MDC_ATTR_NU_VAL_OBS_BASIC, 2
0x0A 0x82 0x00		MDC_ATTR_TIME_STAMP_BO, 8
$0x00 \ 0x06$		obj-class = MDC_MOC_VMO_METRIC_NU
$0x00 \ 0x02$		obj-handle = $2 \leftrightarrow 2^{\text{nd}}$ measurement is current basal rate setting)
$0x00 \ 0x04$		attributes.count = 4
$0x00 \ 0x24$		attributes.length = 36
0x09 0x2F		attribute-id = MDC_ATTR_ID_TYPE
$0x00 \ 0x04$		attribute-value.length = 4
$0x00 \ 0x02 \ 0x74$	0×0.1	MDC_PART_SCADA MDC_INS_BASAL_RATE_SETTING
0x0A 0x46	0.1.0 1	attribute-id = MDC_ATTR_METRIC_SPEC_SMALL
$0x00 \ 0x02$		attribute-value.length = 2
$0x40 \ 0x40$		Metric-Spec-Small (mss-avail-stored-data, mss-acc-agent-initiated)
$0x09 \ 0x96$		attribute-id = MDC_ATTR_UNIT_CODE
0x00 0x02		attribute-value.length = 2
$0x16 \ 0x40$		MDC_DIM_INTL_UNIT_PER_HR
0x0A 0x55		attribute-id = MDC_ATTR_ATTRIBUTE_VAL_MAP
0x0110x95 $0x00 0x0C$		attribute-value.length = 12
$0x00 \ 0x02$		AttrValMap.count = 2
$0x00 \ 0x02 \ 0x00 \ 0x08$		AttrValMap.length = 8
0x00 0x00 0x00 0x0A 0x4C 0x00	0×0.2	MDC_ATTR_NU_VAL_OBS_BASIC, 2
0x0A 0x4C 0x00 0x0A 0x82 0x00		MDC_ATTR_TIME_STAMP_BO, 8
0A0A 0A02 0A00	0.000	11DC_1111K_11111L_01/11111 _DO, 0

8.4.2.2 Manager procedure

The manager shall respond to a configuration notification message using a "Remote Operation Response | Confirmed Event Report" data message with an MDC_NOTI_CONFIG event using the ConfigReportRsp

Health informatics—Personal health device communication—Part 10419: Device Specialization—Insulin Pump

structure for the *event-info* field (see Table 3). As a response to the standard configuration notification message in 8.4.2.1 the format and contents of the manager's configuration notification response message are as follows:

0xE7	0x0			APDU CHOICE Type (PrstApdu)
0x00	0x16			CHOICE.length = 22
0x00	0x14			OCTET STRING.length = 20
0x00	0x02			invoke-id = $0x0002$ (mirrored from invocation)
0x02	0x01			CHOICE (Remote Operation Response Confirmed Event Report)
0x00	0x0E			CHOICE.length = 14
0x00	0x0			obj-handle = 0 (MDS object)
0x00	0x0	0x00	0x00	currentTime = 0
0x0D	0x1C			event-type = MDC_NOTI_CONFIG
0x00	0×04			event-reply-info.length = 4
0x07	0x6C			ConfigReportRsp.config-report-id = 1900
0x00	0x0			ConfigReportRsp.config-result = accepted-config

8.5 Operating procedure

8.5.1 General

Measurement data and status information are communicated from the insulin pump agent during the Operating state. If not stated otherwise, the operating procedure for an insulin pump agent of this standard shall be as specified in IEEE Std 11073-20601-2014.

8.5.2 GET insulin pump MDS attributes

See Table 4 for a summary of the GET service.

If the *attribute-id-list* field in the roiv-cmip-get service message is empty, the insulin pump agent shall respond with a rors-cmip-get service message in which the attribute-list contains a list of all implemented attributes of the MDS object.

If the manager requests specific MDS object attributes, indicated by the elements in *attribute-id-list*, and the agent supports this capability, the insulin pump agent shall respond with a rors-cmip-get service message in which the attribute-list contains a list of the requested attributes of the MDS object that are implemented. It is not required for an insulin pump agent to support this capability. If this capability is not implemented, the insulin pump agent shall respond as specified in the MDS object attributes clause in IEEE Std 11073-20601-2014.

8.5.3 Measurement data transmission

See Table 3 and Table 28 for a summary of the event report services available for measurement data transfer.

To limit the amount of data being transported within an APDU, the insulin pump agent shall not include more than 25 temporarily stored measurements in a single event report. If more than 25 pending measurements are available for transmission, they shall be sent using multiple event reports or using the PM-store mechanism. If multiple insulin pump measurements are available, up to 25 measurements should be transmitted within a single event report. Alternatively, they may be transmitted using a single event

Health informatics—Personal health device communication—Part 10419: Device Specialization—Insulin Pump

report for each insulin pump measurement. However, the former strategy is recommended to reduce overall message size and power consumption.

8.6 Time synchronization

Time synchronization between an insulin pump agent and a manager may be used to coordinate the clocks used when reporting physiological events. Note that the mechanism for synchronizing an agent to a manager is outside the scope of this standard. If time synchronization is used, then this shall be reported in the Mds-Time-Info attribute of the MDS object.

9. Test associations

The Test association provides a manufacturer the mechanism to test or demonstrate features of a product in a comprehensive manner. This clause defines the behavior of the standard insulin pump agent during a test association. Support for test association is optional.

9.1 Behavior with standard configuration

An agent or manager entering a test association using the configuration ID for the standard insulin pump device of this standard shall enter the Operating state in test mode. When in test mode, where possible, this should be indicated visually to any user. Normal functionality shall be suspended, and any test data generated shall not be processed by the device as physiological data.

The insulin pump agent shall send a single simulated current basal rate set value of 199.99 IU/h and a single simulated bolus delivery value of 299.99 IU (values never seen in normal usage and outside normal range) within 30 seconds of entering the Operating state. If the measurement-status attribute of the numeric object is implemented, then the test-data bit shall be set.

The test association is terminated in a manner consistent with the agent's normal behavior for terminating an association.

9.2 Behavior with extended configurations

This specification does not define a test association that uses an extended configuration.

10. Conformance

10.1 Applicability

This standard shall be used in conjunction with IEEE Std 11073-20601-2014.

An implementation or a system can conform to the following elements of this standard:

 Domain information model class hierarchy and object definitions (object attributes, notifications, methods, and data type definitions)

Health informatics—Personal health device communication—Part 10419: Device Specialization—Insulin Pump

- Nomenclature code values
- Protocol and service models
- Communication service model (association and configuration)

10.2 Conformance specification

This standard offers levels of conformance with respect to strict adherence to the standard device and the use of extensions for:

- Information model of a specific device
- Use of attributes, value ranges, and access methods

A vendor shall specify the level of conformance for an implementation based on this standard and provide details of the way in which the definitions of this standard and any extensions are applied.

Specifications shall be provided in the form of a set of implementation conformance statements (ICS) as detailed in 10.4.

This standard is used in conjunction with IEEE Std 11073-20601-2014. It is recommended that the ICS for this standard be created first so that the ICS created for IEEE Std 11073-20601-2014 may refer to the ICS for this standard where applicable.

10.3 Levels of conformance

10.3.1 General

This standard defines the following levels of conformance.

10.3.2 Conformance level 1: base conformance

The application uses elements of the information, service, and communication models (object hierarchy, actions, event reports, and data type definitions) and the nomenclature scheme defined in IEEE Std 11073-20601-2014 and the ISO/IEEE 11073-104zz documents. All mandatory features defined in the object definition tables and in the ICS tables are implemented. Furthermore, any conditional, recommended, or optional features that are implemented shall follow the requirements in the IEEE Std 11073-20601-2014 and ISO/IEEE 11073-104zz documents.

10.3.3 Conformance level 2: extended nomenclature (ASN.1 and/or ISO/IEEE 11073-10101:2004 [B6])

Conformance level 2 meets conformance level 1 but also uses or adds extensions in at least one of the information, service, or nomenclature models. Extensions to nomenclature codes shall conform to the ISO/IEEE 11073-10101:2004 [B6] framework and lie within the private nomenclature extension range (0xF000-0xFFFF).

Extensions to the information or service models shall be fully defined using ASN.1 where appropriate and have their behavior fully described following the framework of the IEEE Std 11073-20601-2014 and/or

Health informatics—Personal health device communication—Part 10419: Device Specialization—Insulin Pump

ISO/IEEE 11073-20101:2004 [B8]. All extensions shall be specified and include reference to the definition for the extension, or where no publicly available reference is available, the definition of the extension should be appended to the conformance statement.

10.4 Implementation conformance statements

10.4.1 General format

The ICSs are provided as an overall conformance statement document that comprises a set of tables in the form given by the templates in the following clauses.

Each ICS table has the following columns:

The table column headings have the following meaning:

- Index: an identifier (e.g., a tag) of a specific feature.
- Feature: briefly describes the characteristic for which a conformance statement is being made.
- Reference: to the clause/paragraph within this document or an external source for the definition of the feature (may be empty).
- Req/Status: specifies the conformance requirement (e.g., mandatory, recommended, etc.)—in some
 cases, this standard does not specify conformance requirements but requests the status of a
 particular feature be provided.
- Support: specifies the presence or absence of a feature and any description of the characteristics of the feature in the implementation. This column is to be filled out by the implementer.
- Comment: contains any additional information on the feature. This column is to be filled out by the implementer.

Subclauses 10.4.2 through 10.4.6 specify the format of the specific ICS tables.

10.4.2 General implementation conformance statement

The general ICS specifies the versions/revisions that are supported by the implementation and high-level system behavior.

Table 37 shows general ICSs.

Table 37—IEEE Std 11073-10419 general ICS table

Indexa	Feature	Reference	Req./Status	Support	Comment
GEN 11073- 10419-1	Implementation description	_	Identification of the device/ application. Description of functionality.		
GEN 11073- 10419-2	Standards followed and their revisions	(standard documents)	(set of existing revisions)	(set of supported revision)	
GEN 11073- 10419-3	Nomenclature document used and revision	(standard documents)	(set of existing revisions)	(set of supported revisions)	
GEN 11073- 10419-4	Conformance adherence Level 1	See 10.3.2	Base conformance declaration that device meets the following IEEE Std 11073-10419 conformance requirements: a) All mandatory requirements shall be implemented. b) If implemented, conditional, recommended, and optional requirements shall conform to this standard.	Yes/No (No is not expected as no implies that the implementation is non-conformant)	
GEN 11073- 10419-5	Conformance adherence Level 2	See 10.3.3	In addition to GEN 11073-10419-4, if the device implements extensions and/or additions, they shall conform to nomenclature codes from ASN.1 and/or 10101 framework. These extensions should also be defined in ICS tables pointing toward their reference.	Yes/No	
GEN 11073- 10419-6	Object containment tree	See 6.3	Provide object containment diagram showing relations between object instances used by the application. A conforming implementation uses only object relations as defined in the DIM.		
GEN 11073- 10419-7	Nomenclature document used and revision	(standard documents)	(set of existing revisions)	(set of supported revision)	
GEN 11073- 10419-8	Data structure encoding	_	_	Description of encoding method(s) for ASN.1 data structures	
GEN 11073- 10419-9	Use of private objects		Does the implementation use objects that are not defined in the DIM?	Yes/No (If yes: explain in Table 38)	

Health informatics—Personal health device communication—Part 10419: Device Specialization—Insulin Pump

Table 37—IEEE Std 11073-10419 general ICS table (continued)

Indexa	Feature	Reference	Req./Status	Support	Comment
GEN 11073- 10419-10	Use of private nomenclature extensions		Does the implementation use private extensions to the nomenclature (i.e., 0xF000–0xFFFF codes from ISO/IEEE 11073-10101:2004 [B6])? Private Nomenclature extensions are <i>only</i> allowed if the standard nomenclature does not include the specific terms required by the application.	Yes/No (If yes: explain in the Table 41)	
GEN 11073- 10419-11	11073-20601 Conformance		Provide the conformance report required by the IEEE Std 11073-20601-2014.		

^aThe prefix GEN11073-10419 is used for the index in the general ICSs table.

10.4.3 DIM MOC implementation conformance statement

The DIM MOC ICS defines which objects are implemented. Information on each object shall be provided as a separate row in the template of Table 38.

Table 38—Template for DIM MOC ICS table

Index	Feature	Reference	Req./Status	Support	Comment
MOC-n	Object	Reference to the	Implemented	Specify	
	description	clause in the		restrictions, e.g.,	
		standard or other		max. number of	
		location where the		supported	
		object is defined.		instances	

The n in the Index column should be the object handle for implementations that have predefined objects. Otherwise the Index column shall simply be a unique number (1..m).

All private objects should be specified and include either a reference to the definition for the object, or where no publicly available reference is available, the definition of the object should be appended to the conformance statement.

The Support column should indicate any restrictions for the object implementation.

An object containment diagram (class instance diagram) should be provided as part of the DIM MOC ICS.

10.4.4 MOC attribute ICS

The MOC attribute ICS defines which attributes, including any inherited attributes, are used/supported in each object of an implementation. Information on each attribute of an object shall be provided as a separate row in the template of Table 39. A separate MOC attribute ICS shall be provided for each object.

Health informatics—Personal health device communication—Part 10419: Device Specialization—Insulin Pump

Table 39—Template for MOC attribute ICS table

Index	Feature	Reference	Req./Status	Support	Comment
ATTR-n-x	Attribute name. Extended attributes shall include the attribute ID also.	Fill in the reference to the ASN.1 structure if the attribute is not defined in this standard.	M = Mandatory C = Conditional R = Recommended O = Optional (as per definition in attribute definition tables)	Implemented? Yes/No Static/Dynamic Specify restrictions, (e.g., value ranges). Describe how attribute is accessed (e.g., Get, Set, sent in config event report, sent in a data event report). Describe any specific restrictions.	

The Support column shall specify whether the attribute is implemented; for extension attributes, whether the attribute value is static or dynamic; any value ranges; restrictions on attribute access or availability; and any other information.

The n in the Index column refers to the ID of the managed object for which the table is supplied (i.e., the index of the managed object as specified in the MOC ICS). There is one separate table for each supported managed object.

The x in the Index column is a unique serial number (1..m).

10.4.5 MOC notification implementation conformance statement

The MOC notification ICS specifies all implemented notifications (typically in the form of the event report service) that are emitted by the agent. Table 40 provides a template for use. One table has to be provided for each object that supports special object notifications. One row of the table shall be used for each notification.

Table 40 — Template for MOC notification ICS table

Index	Feature	Reference	Req./Status	Support	Comment
NOTI-n-x	Notification Name	Reference to the clause		The Support column	
	and Notification ID	in the standard or other		shall specify how the	
		location where the		notification is sent and	
		event is defined.		any restrictions.	

The n in the Index column refers to the ID of the managed object for which the table is supplied (i.e., the index of the managed object as specified in the MOC ICS). There is one separate table for each managed object that supports specific object notifications (i.e., events).

The x in the Index column is a unique serial number (1..m).

Health informatics—Personal health device communication—Part 10419: Device Specialization—Insulin Pump

All private notifications should be specified and include reference to the definition for the notification. Where no publicly available reference is available, the definition of the notification should be appended to the conformance statement.

10.4.6 MOC nomenclature conformance statement

The MOC nomenclature ICS specifies all nonstandard nomenclature codes that are utilized by the agent. Table 41 provides a template for use. One row of the table is to be used for each nomenclature element.

Table 41 —Template for MOC nomenclature ICS table

Index	Feature	Reference	Req./Status	Support	Comment
NOME-n	Nomenclature	Reference to the		Describe how the	
	Name and	clause in the		nomenclature is used.	
	Nomenclature value	standard or other		Describe any specific	
		location where the		restrictions.	
		nomenclature is			
		defined or used			

The n in the Index column is a unique serial number (1..m).

Health informatics—Personal health device communication—Part 10419: Device Specialization—Insulin Pump

Annex A

(informative)

Bibliography

Bibliographical references are resources that provide additional or helpful material but do not need to be understood or used to implement this standard. Reference to these resources is made for informational use only.

- [B1] IEC 60601-1:2005, Ed. 3, Medical electrical equipment—Part 1: General requirements for basic safety and essential performance.⁷
- [B2] IEC 60601-2, Medical electrical equipment—Part 2: Particular requirements for the basic safety and essential performance for specific device. [See the entire series of standards, Part 2-1 through Part 2-51.]
- [B3] IEC 62304:2006/EN 62304:2006, Medical device software—Software life-cycle processes.8
- [B4] IEC 80001-1:2010, Application of risk management for IT-networks incorporating medical devices—Part 1: Roles, responsibilities, and activities.
- [B5] ISO 14971:2007, Medical devices—Application of risk management to medical devices.
- [B6] ISO/IEEE 11073-10101:2004, Health informatics—Point-of-care medical device communication— Part 10101: Nomenclature. 10
- [B7] ISO/IEEE 11073-10201:2004, Health informatics—Point-of-care medical device communication— Part 10201: Domain information model.
- [B8] ISO/IEEE 11073-20101:2004, Health informatics—Point-of-care medical device communication— Part 20101: Application profile—Base standard.
- [B9] ITU-T Rec. X.680-2002, Information technology—Abstract Syntax Notation One (ASN.1): Specification of basic notation.¹¹
- [B10] Miles A. A., M. V. Musset, and W. L. M. Perry, "Third international standard for insulin." of the World Health Organization, vol. pp. 445–459, http://www.ncbi.nlm.nih.gov/pmc/articles/PMC2554140/pdf/bullwho00630-0071.pdf.

⁷ IEC publications are available from the Sales Department of the International Electrotechnical Commission, Case Postale 131, 3 rue de Varembé, CH-1211, Genève 20, Switzerland/Suisse (http://www.iec.ch/). IEC publications are also available in the United States from the Sales Department, American National Standards Institute, 11 West 42nd Street, 13th Floor, New York, NY 10036, USA. ⁸ EN publications are available from the European Committee for Standardization (CEN), 36, rue de Stassart, B-1050 Brussels,

Belgium (http://www.cenorm.be). ISO publications are available from the ISO Central Secretariat, Case Postale 56, 1 rue de Varembé, CH-1211, Genève 20, Switzerland/ Suisse (http://www.iso.ch/). ISO publications are also available in the United States from the Sales Department,

American National Standards Institute, 11 West 42nd Street, 13th Floor, New York, NY 10036, USA (http://www.ansi.org/) ¹⁰ ISO/IEEE publications are available from the ISO Central Secretariat, 1, ch. de la Voie-Creuse, Case postale 56, CH-1211, Geneva 20, Switzerland (http://www.iso.ch/). ISO/IEEE publications are also available in the United States from the Institute of Electrical and Electronics Engineers, 445 Hoes Lane, Piscataway, NJ 08854-4141, USA (http://standards.ieee.org/).

¹¹ ITU publications are available from the International Telecommunications Union, Place des Nations, 1211 Geneva 20, Switzerland (http://www.itu.in/).

Health informatics—Personal health device communication—Part 10419: Device Specialization—Insulin Pump

Annex B

(normative)

Any additional ASN.1 definitions

B.1 Device status and insulin pump status bit mapping

The extension to the enumeration class for operational status requires the following ASN.1 structure definition:

```
OpStat ::= BITS-16 {
         insulin-device-op-undetermined (0),
         insulin-device-op-off (1),
         insulin-device-op-standby (2),
         insulin-device-op-preparing (3),
         insulin-device-op-priming (4),
         insulin-device-op-waiting (5),
         insulin-device-op-ready (6),
         -- reserved for future extension (7),
         -- reserved for future extension (8),
         insulin-device-therapy-stop (9),
         insulin-device-therapy-pause (10),
         insulin-device-therapy-run (11)
         -- reserved for future extension (12),
         -- reserved for future extension (13),
         -- reserved for future extension (14),
         -- reserved for future extension (15),
}
```

The extension to the enumeration class for PHD DM status requires the following ASN.1 structure definition:

```
PHDDMStat ::= BITS-32 {
         device-status-undetermined (0),
         device-status-reset (1),
         -- reserved for future extension (2),
         -- reserved for future extension (3),
         -- reserved for future extension (4),
         device-status-error (5),
         device-status-error-mechanical (6),
         device-status-error-electronic (7),
         device-status-error-software (8),
         device-status-error-battery (9),
         -- reserved for future extension (10),
         -- reserved for future extension (11),
         -- reserved for future extension (12),
         -- reserved for future extension (13),
         -- reserved for future extension (14),
         device-status-service (15),
         device-status-service-time-sync-required (16)
```

Health informatics—Personal health device communication—Part 10419: Device Specialization—Insulin Pump

```
device-status-service-calibration-required (17),
        device-status-service-replenishment-required (18),
        -- reserved for future extension (19),
        -- reserved for future extension (20),
        -- reserved for future extension (21),
        -- reserved for future extension (22),
        -- reserved for future extension (23),
        -- reserved for future extension (24),
        device-status-battery-low (25),
        device-status-battery-depleted (26),
        device-status-battery-replaced (27),
        device-status-battery-interrupted (28)
        -- reserved for future extension (29),
        -- reserved for future extension (30),
        -- reserved for future extension (31),
}
```

The extension to the enumeration class for insulin pump status requires the following ASN.1 structure definition.

```
InsPumpStat ::= BITS-32 {
                 air-pressure-out-of-range (0)
                 bolus-canceled (1)
                 delivery-max (2)
                 infusion-set-detached (3)
                 infusion-set-incomplete (4)
                 occlusion-detected (5)
                 power-insufficient (6)
                 priming-issue (7)
                 reservoir-empty (8)
                 reservoir-issue (9)
                 reservoir-low (10)
                 reservoir-attached (11)
                 temp-basal-canceled (12)
                 temp-basal-expired (13)
                 temperature-out-of-range (14)
```

Annex C

(normative)

Allocation of identifiers

C.1 General

This annex contains the nomenclature codes used in this document and not found in IEEE Std 11073-20601-2014. For those not contained in this annex, the normative definition is found in IEEE Std 11073-20601-2014.

C.2 Definitions of terms and codes

The format used here follows that of ISO/IEEE 11073-10101:2004 [B6].

/*************************************	*****	**********	**
* From Communication Infrastructure (MDC_PART_IN			
************			,
#define MDC_DEV_SPEC_PROFILE_INSULIN_PUMI	24115	/* Insulin pump	*/
/*************************************	*****	**********	k*
* From Object Infrastructure (MDC_PART_OBJ)			
**************************************	****	**********	**/
#define MDC MOC VMO SCHEDSTORE	81	/*	*/
#define MDC_MOC_VMO_SCHEDSTORE #define MDC_MOC_SCHEDULE_SEG	82	/ /*	*/
#define MDC_ATTR_SCHED_STORE_HANDLE	2800	<i>'</i>	*/
#define MDC_ATTR_SCHED_STORE_ACTIVE_INST		,	,
#defile MDC_ATTR_SCHED_STORE_ACTIVE_INST	2801	/*	*/
#define MDC_ATTR_SCHED_STORE_UPDATED_IN			,
#define MDC_ATTR_SCHED_STORE_OFDATED_IN	2802	/*	*/
#define MDC_ATTR_SCHED_STORE_CAPAB	2803	•	*/
#define MDC_ATTR_SCHED_STORE_CAPAC_CNT	2804	•	*/
#define MDC_ATTR_SCHED_STORE_USAGE_CNT		·	*/
#define MDC_ATTR_SCHED_STORE_OP_STAT	2806		*/
#define MDC_ATTR_SCHED_STORE_LABEL_STRIN		,	,
WIGHING WIDE_ATTR_SCHED_STORE_EADEL_STRIP	2807	/*	*/
#define MDC ATTR SCHED SEG NUM	2816	•	*/
#define MDC_ATTR_SCHED_SEG_INSTNO	2817		*/
#define MDC_ATTR_SCHED_SEG_MAP	2818	•	*/
#define MDC ATTR SCHED SEG PERIOD	2819	•	*/
#define MDC_ATTR_SCHED_SEG_ENTRY_INTERV		,	,
weemie wie e_ntrik_seniee_see_entriki_ntrekty	2820	/*	*/
#define MDC_ATTR_SCHED_SEG_PERSON_ID	2821	·	*/
# #define MDC ATTR SCHED SEG ENTRY CNT	2822	<i>'</i>	*/
#define MDC_ATTR_SCHED_SEG_LABEL_STRING	2823	•	*/
#define MDC_ATTR_SCHED_SEG_LAST_UPDATED			
	2828		*/
#define MDC_ATTR_SCHED_SEG_LAST_UPDATED			•
	2829	_	*/
	/		

#define MDC_ATTR_SCHED_SEG_LAST_UPDATED	во т	TME	
	2830		k/
#define MDC_ATTR_SCHED_SEG_REF_ABS_TIME	2831	/*	* /
#define MDC_ATTR_SCHED_SEG_REF_BO_TIME	2832	/*	* /
#define MDC_ATTR_SCHED_SEG_START_ABS_TIME	ΙE		
	2833	/* *	¢/
#define MDC_ATTR_SCHED_SEG_START_BO_TIME	3		
	2834	•	*/
#define MDC_ATTR_SCHED_SEG_END_ABS_TIME	2835		k/
#define MDC_ATTR_SCHED_SEG_END_BO_TIME	2836	•	k/ k/
#define MDC_ATTR_SCHED_SEG_FIXED_DATA #define MDC_ATTR_SCHED_SEG_CONFIRM_TIME(2840	<i></i>	`/
#defille MDC_ATTK_SCHED_SEG_CONTRM_TIMES	2841	/*	* /
#define MDC_ATTR_SCHED_SEG_TRANSFER_TIMI		,	′
"define NIDC_NTTR_SCHED_SEG_TRANSFER_TIME	2842	/*	* /
#define MDC_ACT_SCHED_SEG_GET_INFO	3108	·	*/
#define MDC_ACT_SCHED_SEG_GET_ID_LIST	3109	/*	k/
#define MDC_ACT_SCHED_SEG_TRIG_XFER	3110	/*	* /
#define MDC_ACT_SCHED_SEG_DATA	3111	/*	* /
/ ************************************			
* From Medical supervisory control and data acquisition			
*****************	******	***********	۴/
/*************************************	*****	**********	
* From Personal Health Device Disease Management (M			
**************************************	****	******************************	k /
#define MDC_INS_BASAL			, k/
#define MDC_INS_BASAL_RATE_SETTING			*/
#define MDC_INS_BASAL_PRGM		•	k/
#define MDC_INS_BASAL_TEMP_ABS	29694	1/* Temporary basal rate, absolute *	k/
#define MDC_INS_BASAL_TEMP_REL	29695	5/* Temporary basal rate, relative *	¢/
#define MDC_INS_BASAL_UNDETERMINED			^k /
#define MDC_INS_BASAL_DEVICE			*/
#define MDC_INS_BASAL_REMOTE		,	^k /
#define MDC_INS_BASAL_AP_CTRL			k/
#define MDC_INS_BASAL_OTHER	29700) /* Rate set by an other source *	k/
#define MDC_INS_BASAL_RATE_SCHED	20712) /* Decel wate cohodule cotting *	
#define MDC_INS_BASAL_KATE_SCHED	29/12	2 /* Basal rate schedule setting *	k /
#define MDC_INS_BOLUS_SETTING			k/
"define MDC_INS_BOECS_SETTING	29724	1 /* Rolus amount set *	
#define MDC_INS_BOLUS			*/
#define MDC_INS_BOLUS #define MDC_INS_BOLUS_FAST	29736	6/* Delivered bolus insulin *	*/ */
#define MDC_INS_BOLUS_FAST	29736 29737	6/* Delivered bolus insulin */ 7/* Fast bolus *	*/
	29736 29737 29738	5/* Delivered bolus insulin *7/* Fast bolus *8/* Extended bolus *	*/ */ */
#define MDC_INS_BOLUS_FAST #define MDC_INS_BOLUS_EXT	29736 29737 29738 29739	5/* Delivered bolus insulin 7/* Fast bolus 8/* Extended bolus 9/* Correction bolus **	*/ */ */
#define MDC_INS_BOLUS_FAST #define MDC_INS_BOLUS_EXT #define MDC_INS_BOLUS_CORR	29736 29737 29738 29739 29740	5/* Delivered bolus insulin 7/* Fast bolus 8/* Extended bolus 9/* Correction bolus 9/* Meal bolus **	*/ */ */ */
#define MDC_INS_BOLUS_FAST #define MDC_INS_BOLUS_EXT #define MDC_INS_BOLUS_CORR #define MDC_INS_BOLUS_MEAL #define MDC_INS_BOLUS_UNDETERMINED	29736 29737 29738 29739 29740 29741	5/* Delivered bolus insulin 7/* Fast bolus 8/* Extended bolus 9/* Correction bolus 9/* Meal bolus 1/* Undetermined bolus **	*/ */ */ */
#define MDC_INS_BOLUS_FAST #define MDC_INS_BOLUS_EXT #define MDC_INS_BOLUS_CORR #define MDC_INS_BOLUS_MEAL #define MDC_INS_BOLUS_UNDETERMINED #define MDC_INS_BOLUS_MANUAL	29736 29737 29738 29739 29740 29741	5 /* Delivered bolus insulin 7 /* Fast bolus 8 /* Extended bolus 9 /* Correction bolus 9 /* Meal bolus 1 /* Undetermined bolus 2 /* Manual, user defined bolus	*/ */ */ */ */
#define MDC_INS_BOLUS_FAST #define MDC_INS_BOLUS_EXT #define MDC_INS_BOLUS_CORR #define MDC_INS_BOLUS_MEAL #define MDC_INS_BOLUS_UNDETERMINED #define MDC_INS_BOLUS_MANUAL #define MDC_INS_BOLUS_RECOMMENDED	29736 29737 29738 29739 29740 29741 29742 29743	5 /* Delivered bolus insulin 7 /* Fast bolus 8 /* Extended bolus 9 /* Correction bolus 0 /* Meal bolus 1 /* Undetermined bolus 2 /* Manual, user defined bolus 8 /* Recommended bolus **	*/ */ */ */ */
#define MDC_INS_BOLUS_FAST #define MDC_INS_BOLUS_EXT #define MDC_INS_BOLUS_CORR #define MDC_INS_BOLUS_MEAL #define MDC_INS_BOLUS_UNDETERMINED #define MDC_INS_BOLUS_MANUAL	29736 29737 29738 29739 29740 29741 29742 29743	5 /* Delivered bolus insulin 7 /* Fast bolus 8 /* Extended bolus 9 /* Correction bolus 9 /* Meal bolus 1 /* Undetermined bolus 2 /* Manual, user defined bolus 8 /* Recommended bolus 4 /* Recommended bolus changed by	*/ */ */ */ */
#define MDC_INS_BOLUS_FAST #define MDC_INS_BOLUS_EXT #define MDC_INS_BOLUS_CORR #define MDC_INS_BOLUS_MEAL #define MDC_INS_BOLUS_UNDETERMINED #define MDC_INS_BOLUS_MANUAL #define MDC_INS_BOLUS_RECOMMENDED #define MDC_INS_BOLUS_MANUAL_CHANGE	29736 29737 29738 29739 29740 29741 29742 29744	5 /* Delivered bolus insulin 7 /* Fast bolus 8 /* Extended bolus 9 /* Correction bolus 9 /* Meal bolus 1 /* Undetermined bolus 2 /* Manual, user defined bolus 8 /* Recommended bolus 4 /* Recommended bolus changed by user 8	*/ */ */ */ */ */
#define MDC_INS_BOLUS_FAST #define MDC_INS_BOLUS_EXT #define MDC_INS_BOLUS_CORR #define MDC_INS_BOLUS_MEAL #define MDC_INS_BOLUS_UNDETERMINED #define MDC_INS_BOLUS_MANUAL #define MDC_INS_BOLUS_RECOMMENDED #define MDC_INS_BOLUS_MANUAL_CHANGE #define MDC_INS_BOLUS_COMMANDED	29736 29737 29738 29739 29740 29741 29742 29743 29744	5 /* Delivered bolus insulin 7 /* Fast bolus 8 /* Extended bolus 9 /* Correction bolus 9 /* Meal bolus 1 /* Undetermined bolus 2 /* Manual, user defined bolus 8 /* Recommended bolus 4 /* Recommended bolus changed by user 5 /* Commanded bolus	*/ */ */ */ */ */ */
#define MDC_INS_BOLUS_FAST #define MDC_INS_BOLUS_EXT #define MDC_INS_BOLUS_CORR #define MDC_INS_BOLUS_MEAL #define MDC_INS_BOLUS_UNDETERMINED #define MDC_INS_BOLUS_MANUAL #define MDC_INS_BOLUS_RECOMMENDED #define MDC_INS_BOLUS_MANUAL_CHANGE	29736 29737 29738 29739 29740 29741 29742 29743 29745 29746	5 /* Delivered bolus insulin 7 /* Fast bolus 8 /* Extended bolus 9 /* Correction bolus 9 /* Meal bolus 1 /* Undetermined bolus 2 /* Manual, user defined bolus 8 /* Recommended bolus 4 /* Recommended bolus changed by user 5 /* Commanded bolus 6 /* Other bolus	*/ */ */ */ */ */

ISO/IEEE 11073-10419:2016(E)

IEEE Std 11073-10419-2015
Health informatics—Personal health device communication—Part 10419: Device Specialization—Insulin Pump

#define MDC_INS_I2CHO_SCHED #define MDC_INS_ISF_SCHED #define MDC_INS_RESERVOIR #define MDC_INS_CONC	29756 /* I:CHO Schedule Setting 29768 /* ISF Schedule Setting 29780 /* Insulin reservoir remaining 29792 /* Insulin concentration	*/ */ */
#define MDC_INS_PUMP_OP_STAT #define MDC_PHD_DM_DEV_STAT	29804 /* Operational status 20000 /* PHD DM Device status	*/ */
#define MDC_INS_PUMP_DEV_STAT	29836/* Insulin Pump Device status	*/
/*************************************		
#define MDC_DIM_DIMLESS	512 /* Dimensionless	*/
#define MDC_DIM_G	1728 /* Gram	*/
#define MDC_DIM_INTL_UNIT	5472 /* International unit	*/
#define MDC_DIM_INTL_UNIT_PER_HR	5696 /* International units per hour	*/
#define MDC_DIM_INTL_UNIT_PER_L	5568 /* International units per liter	*/
#define MDC_DIM_INTL_UNIT_PER_ML	5600 /* International units per milliter	*/
#define MDC_DIM_INTL_UNIT_PER_M_CUBE	5536 /* International units	
	per cubic meter	*/
#define MDC_DIM_INTL_UNIT_PER_CM_CUBE	5504 /* International units	
	per cubic centimeter	*/
#define MDC_DIM_MILLI_G_PER_DL	2130 /* Milligram per deciliter	*/
#define MDC_DIM_MILLI_MOLE_PER_L	4722 /* Millimole per liter	*/

C.3 Systematic derivations of terms and codes

Systematic derivations of terms and codes are outlined in Table C.1.

Table C.1—Systematic derivations of terms and codes

Systematic name	Common term	Acronym	Description/definition	Reference ID	Code
Disease Management Device Status	PHD DM Status		Object containing the general device status for PHD disease management	MDC_PHD_DM_DEV_STAT	20000
Personal Health Device			, and the second		
Insulin Basal	Insulin Basal		Object containing the basal delivered amount	MDC_INS_BASAL	29680
Insulin Basal Rate Setting	Insulin Basal Rate Setting		Object containing the current basal rate setting	MDC_INS_BASAL_RATE_SETTING	29692
Insulin Basal Programmed	Insulin Basal Rate Programmed		Basal rate set by the active basal profile schedule setting	MDC_INS_BASAL_PRGM	29693
Insulin Basal Temporary Absolute	Absolute Temporary Insulin Basal Rate		Basal rate set by an absolute temporary change	MDC_INS_BASAL_TEMP_ABS	29694
Insulin Basal Temporary Relative	Relative Temporary Insulin Basal Rate		Basal rate set by an relative temporary change	MDC_INS_BASAL_TEMP_REL	29695
Insulin Basal Undetermined	Undetermined Insulin Basal Rate Context		Undetermined basal rate context	MDC_INS_BASAL_UNDETERMINED	29696
Insulin Basal Device	Insulin Pump Device Set Basal Rate		Basal rate set by the insulin pump device	MDC_INS_BASAL_DEVICE	29697
Insulin Basal Remote	Insulin Pump Remote Control Set Basal Rate		Basal rate set by a remote control	MDC_INS_BASAL_REMOTE	29698
Insulin Basal Artificial Pancreas	Artificial Pancreas Controller Set Basal Rate		Basal rate set by an artificial pancreas controller	MDC_INS_BASAL_AP_CTRL	29699
Insulin Basal Other	Other Source Set Basal Rate		Basal rate set by another source. This option is used when the basal rate source doesn't match an available option.	MDC_INS_BASAL_OTHER	29670
Insulin Basal Rate Schedule	Insulin Basal Rate Segment		Object containing the basal rate segment for a basal rate schedule	MDC_INS_BASAL_RATE_SCHED	29712
Insulin Bolus Setting	Insulin Bolus Setting		Object containing the current bolus setting	MDC_INS_BOLUS_SETTING	29724
Insulin Bolus	Insulin Bolus		Object containing the bolus delivered amount	MDC_INS_BOLUS	29736
Insulin Bolus Fast	Insulin Fast Bolus		Insulin delivered as a fast bolus	MDC_INS_BOLUS_FAST	29737
Insulin Bolus Extended	Insulin Extended Bolus		Insulin delivered as an extended bolus	MDC_INS_BOLUS_EXT	29738
Insulin Bolus Correction	Insulin Correction Bolus		Insulin delivered as a bolus for a correction	MDC_INS_BOLUS_CORR	29739
Insulin Bolus Meal	Insulin Meal Bolus		Insulin delivered as a bolus for a meal	MDC_INS_BOLUS_MEAL	29740
Insulin Bolus Undetermined	Undetermined Insulin Bolus Context		Undetermined bolus context	MDC_INS_BOLUS_UNDETERMINED	29741
Insulin Bolus Manual	Manual Bolus		Bolus defined and set manually by a user	MDC_INS_BOLUS_MANUAL	29742
Insulin Bolus Recommended	Recommended Bolus		Bolus recommended to the user such a bolus calculator	MDC_INS_BOLUS_RECOMMENDED	29743
Insulin Bolus Manual Change	Manually Changed Bolus		Bolus recommended to the user such a bolus calculator but bolus amount changed by the user	MDC_INS_BOLUS_MANUAL_CHANGE	29744

Systematic name	Common term	Acronym	Description/definition	Reference ID	Code
Insulin Bolus Commanded	Commanded Bolus	, and the second	Bolus commanded by another device such as an artificial pancreas controller	MDC_INS_BOLUS_COMMANDED	29745
Insulin Bolus Other	Other Source Set Bolus		Bolus set by another source. This option is used when the bolus source doesn't match an available option.	MDC_INS_BOLUS_OTHER	29746
Insulin Bolus Pending Delay	Bolus Pending Delay		Object containing the bolus pending delay	MDC_INS_BOLUS_PENDING_DELAY	29747
Insulin Carbohydrate Schedule	Insulin to Carbohydrate Ratio		Object containing the insulin to carbohydrate ratio for an insulin to carbohydrate schedule	MDC_INS_I2CHO_SCHED	29756
Insulin Sensitivity Factor Schedule	Insulin Sensitivity Factor		Object containing the insulin sensitivity factor for an insulin sensitivity factor schedule	MDC_INS_ISF_SCHED	29768
Insulin Operational Status	Insulin Device Status		Object containing the insulin pump operational status	MDC_INS_PUMP_OP_STAT	29804
Insulin Device Status	Insulin Device Status		Object containing the insulin pump device status	MDC_INS_PUMP_DEV_STAT	29836
Attribute Schedule Store Handle	Reference Ideintifier		Schedule-store attribute representing a reference ID for this object. Each object shall have a unique ID assigned by the agent. The handle identifies the object in event reports sent to the manager and to address the object instance in messages invoking object methods.	MDC_ATTR_SCHED_STORE_HANDLE	2800
Attribute Schedule Store Instance Number Active Segment	Active Schedule Identifier		Schedule-store attribute representing the unique ID of the specific schedule-segment object that is active. If no schedule segment is currently active or there are currently no schedule segments, the value shall be 0.	MDC_ATTR_SCHED_STORE_ACTIVE_ INSTNO	2801
Attribute Schedule Store Instance Number Updated List	Updated Schedule Identifier(s)		Schedule-store attribute representing the list of schedule-segment object IDs reporting all updated schedule-segments (i.e., created, modified, or deleted)	MDC_ATTR_SCHED_STORE_UPDATED _INSTNO	2802
Attribute Schedule Store Capabilities	Schedule-Store Capabilities		Schedule-store attribute representing basic capabilities of the schedule-	MDC_ATTR_SCHED_STORE_CAPAB	
Attribute Schedule Store Count Capacity	Maximum Capacity Count		store object instance Schedule-store attribute representing the maximum number of stored schedule- segment entries (entries in all contained schedule- segments)	MDC_ATTR_SCHED_STORE_CAPAC_ CNT	2803

Systematic name	Common term	Acronym	Description/definition	Reference ID	Code
Attribute Schedule	Number of		Schedule-store attribute	MDC_ATTR_SCHED_STORE_USAGE_	
Store Count	Schedule		representing the actual	CNT	
Usage	Segments		number of currently stored		
			schedule-segment entries		
			(entries in all contained		
			schedule-segments)		2805
Attribute Schedule	Schedule-Store		Schedule-store attribute	MDC_ATTR_SCHED_STORE_OP_STAT	
Store Operational	Operational		indicating if new entries are		
Status	Status		currently being inserted in		
			any of the contained		
			schedule-segments. If any		
			schedule-segment contained		
			by this schedule-store is being modified, this attribute		
			shall be set to enabled.		
			Otherwise, it shall be set to		
			disabled.		2806
Attribute Schedule	Schedule-Store		Schedule-store attribute	MDC_ATTR_SCHED_STORE_LABEL_	2000
Store Label	Label		representing an application-	STRING	
1 Store Laber	Label		dependent label for the	STRING	
			schedule-store in printable		
			ASCII to indicate its		
			intended use and may be		
			used for display purposes		2807
Attribute Schedule	Number of		Schedule-store representing	MDC ATTR SCHED SEG NUM	2816
Store Number of	Schedule		the number of currently		
Schedule Segments	Segments		instantiated schedule-		
Ü			segments contained in the		
			schedule-store. Note that the		
			schedule-segment attribute		
			Instance-Number is NOT		
			related to this number (i.e.,		
			Instance-Number does not		
			need to be in the range from		
			0 to number-of-schedule-		
M-4-1101-11-1	C + C + L + L		segments)	MDC ACT SCHED SEC CET BIEG	2100
Method Schedule Store Get	Get Schedule		Schedule-store object method that allows the	MDC_ACT_SCHED_SEG_GET_INFO	3108
Schedule Segment	Segment Info		manager to retrieve		
Info			Schedule-segment		
mio			attributes of one or more		
			Schedule-segments, with		
			the exception of the Fixed-		
			Schedule-Segment-Data		
			attribute which contains the		
			actual schedule data and is		
			retrieved by using the Trig-		
			Schedule-Segment-Data-		
			Xfer method		
Method Schedule	Get Schedule		Schedule-store object	MDC_ACT_SCHED_SEG_GET_ID_LIST	3109
Store Get	Segment Data		method that allows the		
Schedule Segment	Transfer		manager to retrieve a list of		
Data Transfer			the instance numbers of all		
			the schedule-segments of a		
M d 1101 111	TD: C1 11	ļ	schedule-store	MDG ACT COHED ONG TRIC MOTE	2110
Method Schedule	Trigger Schedule		Schedule-store object	MDC_ACT_SCHED_SEG_TRIG_XFER	3110
Store Trigger Schedule Segment	Segment Info		method that allows the manager to start the transfer		
Info			of the Fixed-Schedule-		
IIIIO			Segment-Data attribute of a		
			specified schedule-segment		
	l .	1	specifica schedule-segnicit	1	1

Systematic name	Common term	Acronym	Description/definition	Reference ID	Code
Event Schedule	Schedule		Schedule-store object event	MDC_NOTI_SCHED_	3111
Store Schedule	Segment Data		that sends data stored in the	SEG_DATA	
Segment Data	Event		Fixed-Schedule-Segment-		
Event			Data of a schedule segment.		
Attribute Schedule	Schedule		Schedule-segment attribute	MDC_ATTR_SCHED_SEG_INSTNO	
Segment	Segment		representing the ID of a		
Identifier	Identifier		specific schedule-segment		
			object instance. Each		
			instance shall have a unique		
			number assigned by the		
			agent starting at 1. It is used		
			by the manager to address a		
			schedule-segment. Note, the		
			instance-number value of 0		
			shall be reserved for		
			reporting no active		
			schedule-segment in the		
			schedule-store.		2817
Attribute Schedule	Schedule		Schedule-segment attribute	MDC_ATTR_SCHED_SEG_MAP	
Segment Map	Segment Map		representing the format and	MBC_MMC_SCMBB_SEG_MM	
1 Segment 1 Wap	Segment Map		contents of one schedule		
			entry. An entry has a		
			conditional header		
			containing information		
			applicable to all elements in		
			the entry. If the schedule-		
			segment period, start time		
			and entry interval are used,		
			then the header information		
			shall not be used.		
			shan not be used.		
			The entry then contains one		
			or more elements, defined		
			by the handle and an		
			attribute value map defining		
			the object attributes for each		
			element in the schedule-		
			segment.		2818
Attribute Schedule	Schedule		Schedule-segment attribute	MDC_ATTR_SCHED_SEG_PERIOD	2010
Segment Period	Segment Period		representing the period of	MDC_ATTR_SCHED_SEG_TERIOD	
- Segment FEHOU	Segment Ferrod		the schedule-segments. If		
			the schedule-segment occurs		
			only once, then the value		
			shall be 0		2819
Attailanta Calaa J1-	Cahadula			MDC ATTR CCHED SEC ENTRY	2019
Attribute Schedule	Schedule		Schedule-segment attribute	MDC_ATTR_SCHED_SEG_ENTRY_ INTERVAL	
Segment Entry	Segment Entry Interval		representing the interval between each entry in	INTERVAL	
Interval	micrvai		1		2020
			schedule-segments		2820

Systematic name	Common term	Acronym	Description/definition	Reference ID	Code
Attribute Schedule	Person		Schedule-segment attribute	MDC_ATTR_SCHED_SEG_PERSON_ID	
Segment Person	Identification		representing the person ID.		
ID			This standard supports		
			devices that have simple		
			support for data from		
			multiple persons. A person		
			ID is used to differentiate		
			different persons. If the		
			schedule-segment is able to		
			have schedule data for		
			multiple persons, it shall set		
			the schedsc-multi-person bit		
			in the Sched-Store-Capab		
			attribute. If this bit is set, all		
			schedule-segment instances		
			shall support the Schedule-		
			Segment-Person-Id attribute.		
			Otherwise, this attribute is		
			not defined.		2821
Attribute Schedule	Schedule		Schedule-segment attribute	MDC_ATTR_SCHED_SEG_ENTRY_CNT	2822
Segment Count	Segment Entry		representing the actual		
Entry	Count		number of schedule-segment		
			entries.		
Attribute Schedule	Schedule		Schedule-segment attribute	MDC_ATTR_SCHED_SEG_LABEL_	2823
Segment Label	Segment Label		representing an application-	STRING	
	_		dependent label in printable		
			ASCII for the schedule		
			segment to indicate its		
			intended use and may be		
			used for display purposes.		
Attribute Schedule	Last Updated		Schedule-segment attribute	MDC_ATTR_SCHED_SEG_LAST_	2828
Segment Time	Time (Absolute)		representing the time the	UPDATED_ABS_TIME	
Last Update			schedule segment was last		
Absolute			updated. If this attribute is		
			used, neither the Schedule-		
			Segment-LastUpdated-		
			HiRes-Time nor Schedule-		
			Segment-LastUpdated-BO-		
			Time shall be used.		
Attribute Schedule	Last Updated		Schedule-segment attribute	MDC_ATTR_SCHED_SEG_LAST_	2829
Segment Time	Time (High		representing the time the	UPDATED _HIRES_TIME	
Last Update High	Resolution)		schedule segment was last		
Resolution	·		updated. If this attribute is		
			used, neither the Schedule-		
			Segment-LastUpdated-Abs-		
			Time nor Schedule-		
			Segment-LastUpdated-BO-		
	<u> </u>		Time shall be used.		
Attribute Schedule	Last Updated		Schedule-segment attribute	MDC_ATTR_SCHED_SEG_LAST_	2830
Segment Time	Time (Base		representing the time the	UPDATED_BO_TIME	
Last Update Base	Offset)		schedule segment was last		
Offset	1		updated. If this attribute is		
	1		used, neither the Schedule-		
			Segment-LastUpdated-Abs-		
	1		Time nor Schedule-		
	1		Segment- LastUpdated-		
1	1		HiRes-Time shall be used.		

Systematic name	Common term	Acronym	Description/definition	Reference ID	Code
Attribute Schedule Segment Time Reference Absolute	Reference Time (Absolute)		Schedule-segment attribute representing the time from which all entries within the schedule segment are referenced. This attribute shall be used, if the Schedule-Segment-Reference-BO-Time is not used.	MDC_ATTR_SCHED_SEG_REF_ABS_ TIME	2831
Attribute Schedule Segment Time Reference Base Offset	Reference Time (Base Offset)		Schedule-segment attribute representing the time from which all entries within the schedule segment are referenced. This attribute shall be used, if the Schedule-Segment-Reference-Abs-Time is not used.	MDC_ATTR_SCHED_SEG_REF_BO_ TIME	2832
Attribute Schedule Segment Time Start Absolute	Start Time (Absolute)		Schedule-segment attribute representing the time of the first execution of the schedule. If this attribute is used, the Schedule-Segment-Start-BO-Time shall not be used. If this attribute is used, the schedule-segment becomes the active-schedule-segment in the schedule-store at the start time. A schedule-segment shall not be active before the start time.	MDC_ATTR_SCHED_SEG_START_ABS _TIME	2833
Attribute Schedule Segment Time End Absolute	End Time (Absolute)		Schedule-segment attribute representing the expiration of the schedule. If this attribute is used, the Schedule-Segment-End-BO-Time shall not be used. If this attribute is used and the end time is reached, the schedule-segment shall no longer be the active-schedule-segment in the schedule-store.	MDC_ATTR_SCHED_SEG_END_ABS_ TIME	2835
Attribute Schedule Segment Time Start Base Offset	Start Time (Base Offset)		Schedule-segment attribute representing the time of the first execution of the schedule. If this attribute is used, the Schedule-Segment-Start-Abs-Time shall not be used. If this attribute is used, the schedule-segment becomes the active-schedule-segment in the schedule-store at the start time. A schedule-segment shall not be active before the start time.	MDC_ATTR_SCHED_SEG_START_BO_ TIME	2834

Systematic name	Common term	Acronym	Description/definition	Reference ID	Code
Attribute Schedule	End Time (Base		Schedule-segment attribute	MDC_ATTR_SCHED_SEG_END_BO_	2836
Segment Time	Offset)		representing the expiration	TIME	
End Base Offset			of the schedule. If this		
			attribute is used, the		
			Schedule-Segment-End-		
			Abs-Time shall not be used.		
			If this attribute is used and		
			the end time is reached, the		
			schedule-segment shall no		
			longer be the active-		
			schedule-segment in the		
			schedule-store.		
Attribute Schedule	Schedule		Schedule-segment attribute	MDC_ATTR_SCHED_SEG_FIXED_DATA	2840
Segment Fixed	Segment Fixed		representing the schedule	MDC_NTIK_SCHED_SEG_TIKED_DATA	2040
Data	Data		data transferred as an array		
Data	Data		of entries in a format as		
			specified in the Schedule-		
			Segment-Entry-Map		
			attribute. This is defined		
			here as an opaque data		
			structure without a defined		
			data type. Note that this		
			attribute is not directly		
			accessible; it is only		
			retrievable by the manager		
			using the Trig-Schedule-		
			Segment-Data-Xfer method.		
Attribute Schedule	Confirm Timeout		Schedule-segment attribute	MDC_ATTR_SCHED_SEG_CONFIRM_	2841
Segment Timeout			representing the minimum	TIMEOUT	
Confirm			time that the agent shall wait		
			for a response message from		
			the manager after issuing a		
			Confirmed Event Report		
			invoke message before		
			timing out and transitioning		
			to the Unassociated state.		
			This is an informational		
			attribute for the benefit of		
			the manager. If this attribute		
			is supplied, it shall match		
			the actual timeout value that		
			the agent uses for the		
			Confirmed Event Report		
			generated from the		
			Schedule-store object.		
			Senedule store object.		
			This attribute is		
			informational for the		
			manager in the sense that the		
			manager in the sense that the manager does not use this		
			attribute in an actual		
			implementation of the		
			protocol (i.e., the manager		
			does not timeout on an		
			agent-generated Confirmed		
			Event Report). However, the		
			manager might wish to use		
			this information to prioritize		
			its handling of a "short"		
			timeout agent over that of a		
]	"long" timeout agent.		1

Systematic name	Common term	Acronym	Description/definition	Reference ID	Code
Attribute Schedule	Transfer Timeout		Schedule-segment attribute	MDC_ATTR_SCHED_SEG_TRANSFER_	2842
Segment Timeout			representing the minimum	TIMEOUT	
Transfer			time that the manager shall		
			wait for the complete		
			transfer of schedule-segment		
			information.		
			If the timeout expires prior		
			to the reception of the		
			complete schedule-segment,		
			the manager shall transition		
			to the Unassociated state as		
			described in subclause		
			8.9.5.6 of IEEE Std 11073-		
			20601-2014.		

Health informatics—Personal health device communication—Part 10419: Device Specialization—Insulin Pump

Annex D

(informative)

Message sequence examples

Figure D.1 shows a sequence diagram of the messaging procedure corresponding to the following use case. The user of an insulin pump agent device intends to connect it to a manager device for the first time. The insulin pump is capable of performing basal rate set and bolus measurements.

- a) When the user connects the insulin pump, the manager does not recognize the agent's configuration and sends a response to the agent's association request with the result *accepted-unknown-config*. See I.2.2.2 and I.2.2.3 for the corresponding PDU examples.
- b) As a consequence of this, the agent negotiates its configuration information to the manager. After getting confirmation from the manager accepting the agent's configuration, the agent device is ready to send measurements. Both devices enter the Operating state. See I.3.2.2 and I.3.2.3 for the corresponding PDU examples.
- c) Subsequently, the manager may request the MDS object attributes of the agent by sending a data message with the "Remote Operation Invoke | Get" command. Note that the manager may request the MDS object attributes as soon as the agent enters the Associated state, including the Configuring and Operating substates. As a response, the agent reports its MDS object attributes to the manager using a Data message with the "Remote Operation Response | Get" command. See I.4.1.2 and I.4.1.3 for the corresponding PDU examples.
- d) As a next step, the user of the agent device takes several measurements over some period of time. The measurement data are transmitted to the manager using an unconfirmed event reports. See I.5.1 for the corresponding PDU example.
- e) The user ends the measurement session (e.g., by pushing a proper button on the device, or just by not using the device for a duration longer than a certain time period). As a consequence, the agent disassociates from the manager by sending an association release request. The manager responds with an association release response. See I.6.1 and I.6.2 for the corresponding PDU examples.
- f) When the agent requests to associate to the manager for the next measurement session (e.g., the next day), the result in the manager's response is *accepted*, as it already knows the agent's configuration from the previous measurement session. Both devices transition directly to the Operating state.
- g) Finally, the last two steps shown are similar as in item d) and item e). The user takes several unconfirmed measurements followed by releasing the association.

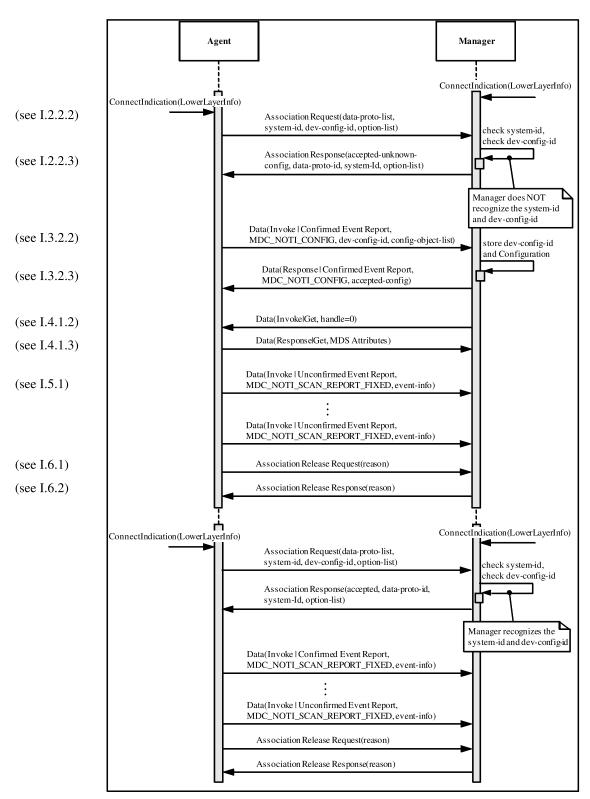


Figure D.1—Sequence diagram for insulin pump example use case

Health informatics—Personal health device communication—Part 10419: Device Specialization—Insulin Pump

Annex E

(normative)

Schedule-store class

E.1 Schedule-store class

E.1.1 General

An instance of the schedule-store class provides storage capabilities for a series of schedules. Data related to these schedules are stored in a variable number of schedule-segment objects (see E.2). The stored schedule data of the schedule-store object are requested from the agent by the manager using object access services (see 7.2). Anybody not familiar with the schedule-store concept may wish to read Annex G for a conceptual overview prior to reading the following subclauses.

Attribute values stored in a schedule-segment may require the use of additional attributes from that object to describe the schedule; a common example being the duration. If an attribute value in a schedule-segment depends on an attribute value not stored in the schedule-segment, then that dependent attribute shall not change value during the lifetime of the schedule-segment. Otherwise, the agent shall store the dependent attribute value in the schedule-segment.

E.1.2 Schedule-store class identification

The nomenclature code to identify the schedule-store class is MDC_MOC_VMO_SCHEDSTORE.

E.1.3 Schedule-store class attributes

Table E.1 defines the set of schedule-store attributes that are supported for personal health device communication:

Table E.1—Schedule-store attributes

Attribute name	Attribute ID	Attribute type	Remark	Qualifiers
Schedule- Handle	MDC_ATTR_SCHED_ STORE_HANDLE	HANDLE	The Schedule Handle attribute represents a reference ID for this object. Each object shall have a unique ID assigned by the agent. The handle identifies the object in event reports sent to the manager and to address the object instance in messages invoking object methods.	Mandatory Static
Active- Schedule- Segment- Instance- Number	MDC_ATTR_SCHED_ STORE_ACTIVE_ INSTNO	InstNumber	The Active-Schedule-Instance- Number is the unique ID of the specific schedule-segment object that is active. If no schedule segment is currently active or there are currently no schedule segments, the value shall be 0.	Mandatory Dynamic
Updated- Schedule- Segment- Instance- Number-List	MDC_ATTR_SCHED_ STORE_UPDATED_ INSTNO	InstNumberList	The Updated-Schedule-Instance- Number-list is the list of schedule-segment object IDs reporting all updated schedule segments (i.e., created, modified, or deleted).	Mandatory Observational
Schedule- Store-Capab	MDC_ATTR_SCHED_ STORE_CAPAB	SchedStoreCapab	This attribute defines basic capabilities of the schedule-store object instance.	Mandatory Static
Schedule- Store- Capacity- Count	MDC_ATTR_SCHED_ STORE_CAPAC_CNT	INT-U32	This attribute is the maximum number of stored schedule-segment entries (entries in all contained schedule segments).	Optional Static
Schedule- Store-Usage- Count	MDC_ATTR_SCHED_ STORE_USAGE_CNT	INT-U32	This attribute is the actual number of currently stored schedule-segment entries (entries in all contained Schedule segments).	Optional Dynamic
Schedule- Store- Operational- Status	MDC_ATTR_SCHED_ STORE_OP_STAT	OperationalState	The attribute indicates if new entries are currently being inserted in any of the contained schedule segments. If any schedule segment contained by this schedule-store is being modified, this attribute shall be set to enabled. Otherwise, it shall be set to disabled.	Mandatory Dynamic

Health informatics—Personal health device communication—Part 10419: Device Specialization—Insulin Pump

Table E.1—Schedule-store attributes (continued)

Attribute name	Attribute ID	Attribute type	Remark	Qualifiers
Schedule- Store-Label	MDC_ATTR_SCHED_ STORE_LABEL_STRING	OCTET STRING	This attribute is an application- dependent label for the schedule- store in printable ASCII to indicate its intended use and may be used for display purposes.	Optional Static
Number-Of- Schedule- Segments	MDC_ATTR_SCHED_ SEG _NUM	INT-U16	This attribute is the number of currently instantiated schedule segments contained in the schedule-store. Note that the schedule-segment attribute Instance-Number is NOT related to this number (i.e., Instance-Number does not need to be in the range from 0 to Number-Of-Schedule-Segments).	Mandatory Dynamic

E.1.4 Schedule-store object methods

Table E.2 defines the methods (actions) of a schedule-store object. These methods can be invoked using the ACTION service.

Table E.2—Schedule-store object methods

Method/Action	Mode	Action-type	action-info-args	Resulting action-info-args
Get-Schedule-	Confirmed	MDC_ACT_SCHED_SE	SchedSegmSelection	SchedSegmentInfoList
Segment-Info		G_GET_INFO		
Get-Schedule-	Confirmed	MDC_ACT_SCHED_SE	(empty)	SchedSegmIdList
Segment-Id-List		G_GET_ID_LIST		
Trig-Schedule-	Confirmed	MDC_ACT_SCHED_SE	TrigSchedSegmData	TrigSchedSegmDataXferRsp
Segment-Data-Xfer		G_TRIG_XFER	Xfer-Req	

If an agent supports the schedule-store class, the support of the Get-Schedule-Segment-Info or Get-Schedule-Segment-Id-List methods is mandatory, and support of the Trig-Segment-Data-Xfer method is mandatory.

If a manager supports the schedule-store class, the support of sending the Get-Schedule-Segment-Info, Get-Schedule-Segment-Id-List and Trig-Segment-Data-Xfer methods is mandatory.

— Get-Schedule-Segment-Info:

This method allows the manager to retrieve Schedule-segment attributes of one or more schedule segments, with the exception of the Fixed-Schedule-Segment-Data attribute which contains the actual schedule data and is retrieved by using the Trig-Schedule-Segment-Data-Xfer method. In particular, the Get-Schedule-Segment-Info method allows the manager to retrieve the attributes and their data contents from the Schedule-segment object instances identified by the SchedSegmSelection parameter.

The agent shall support the all-sched-segments choice in the SchedSegmSelection action-info-args of the Get-Schedule-Segment-Info method. The agent may support the sched-segm-id-list choice in the SchedSegmSelection action-info-args of the Get-Schedule-Segment-Info method. In this case the agent shall set the schedsc-segm-id-list-select flag in the Schedule-Store-Capab attribute. If the manager

Health informatics—Personal health device communication—Part 10419: Device Specialization—Insulin Pump

sends the Get-Schedule-Segment-Info method with the choice that the agent doesn't support, the agent shall reply with an unsupported-choice error (roer).

If the manager supports sending the Get-Schedule-Segment-Info method, the manager shall support at least the choice all-segments in the SchedSegmSelection action-info-args of the Get-Schedule-Segment-Info method. The manager may support additional choices.

If there is no schedule segment that matches the selection criteria in the SchedSegmSelection action-info-args such that no schedule segments are found by the action, then this is not an error and a normal response is sent, and the segment info list will just be empty.

If the choice of SchedSegmSelection in the Get-Schedule-Segment-Info method is sched-segm-id-list and the sched-segm-id-list is empty then the response shall be a sched-segment-info-list that is empty.

If the agent supports the Get-Schedule-Segment-Info method, the agent shall set the schedsc-get-segm-info-sup flag in the Schedule-Store-Capab attribute.

— Get-Schedule-Segment-Id-List:

This method allows the manager to retrieve a list of the instance numbers of all the schedule segments of a schedule-store. In particular, the Get-Schedule-Segm-Id-List method allows the manager to then retrieve the attributes of selected schedule segment object instances and their data contents without need to retrieve information of all schedule segments. This also allows the manager to retrieve multiple schedule segments as a series of requests.

If the agent supports the Get-Schedule-Segment-Id-List method, the agent shall set the schedsc-get-segm-id-list-sup flag in the Schedule-Store-Capab attribute.

If the agent supports the Get-Schedule-Segment-Id-List method, the agent shall support the action of Get-Schedule-Segment-Info by ID list.

— Trig-Schedule-Segment-Data-Xfer:

This method allows the manager to start the transfer of the Fixed-Schedule-Segment-Data attribute of a specified schedule segment. The agent indicates in the response if it accepts or denies this request. If the agent accepts the request, the agent sends Schedule-Segment-Data-Event messages as described in E.1.5. If this method is invoked on a schedule segment that has the Operational-State attribute set to enabled, the agent shall reply with a not-allowed-by-object error (roer) with a return code of MDC_RET_CODE_OBJ_BUSY.

E.1.5 Schedule-store object events

Table E.3 defines the potential events sent by a schedule-store object.

Table E.3—Schedule-store object events

Event	Mode	Event-type	Event-info parameter	Event-reply-info
Schedule-Segment-Data-	Confirmed	MDC_NOTI_SCHED_	ScheduleSegment	ScheduleSegment
Event		SEG_DATA	DataEvent	DataResult

— Schedule-Segment-Data-Event:

This event sends data stored in the Fixed-Schedule-Segment-Data of a schedule segment from the agent to the manager. The event is triggered by the manager by the Trig-Schedule-Segment-Data-Xfer method. Once the data transfer is triggered, the agent sends Schedule-Segment-Data-Event messages until the complete Fixed-Schedule-Segment-Data is transferred or the transfer is aborted by the manager or agent. See Transfer Schedule-segment content in H.1.1.1.2.2.1.1 for a full description.

Health informatics—Personal health device communication—Part 10419: Device Specialization—Insulin Pump

It is encouraged to place as many segment entries contained in a Schedule-Segment-Data-Event as possible to reduce the number of messages required for the transfer of the segment. The agent shall transfer all segment entries in order, first entry first (first in, first out).

Support for the event by the agent is mandatory if the agent supports schedule-store objects.

If confirmed event report is used by a schedule-store object, there shall be at most one unacknowledged confirmed event report outstanding from this object at any point in time.

E.1.6 Other schedule-store services

E.1.6.1 GET service

Support for the GET service shall be provided by any agent that supports one or more schedule-store objects only while in the Operating state. The manager uses the GET service to retrieve the values of all schedule-store object attributes. If a manager does not have the current value of a needed schedule-store attribute, then the GET service shall be used. An agent may also send scan event reports providing the manager with updates of the current attribute values, but this is not a mandated agent behavior.

The manager may request the schedule-store object attributes of the agent in which case the manager shall send the "Remote Operation Invoke | Get" command (see roiv-cmip-get in IEEE Std 11073-20601-2014 subclause A.10.2) with the handle value of the schedule-store object, as defined in the agent's configuration. The agent shall respond by reporting its schedule-store object attributes to the manager using the "Remote Operation Response | Get" response (see rors-cmip-get in IEEE Std 11073-20601-2014 subclause A.10.2). See schedule-store data transmission subclause)H.1.1.1.2.2.1.1) for a full explanation of the GET operation.

E.1.6.2 SET service

There are currently no uses of the schedule-store SET service defined in this standard.

E.2 Schedule-segment class

E.2.1 General

An instance of the schedule-segment represents stored settings that define a series of scheduled entries. A schedule-segment object is not part of the static agent configuration because the number of instantiated schedule-segment instances may dynamically change. The manager accesses schedule-segment objects indirectly by methods and events of the schedule-store object.

E.2.2 Schedule-segment identification

The nomenclature code to identify the schedule-segment class is MDC_MOC_SCHEDULE_SEG.

E.2.3 Schedule-segment attributes

Table E.4 defines the set of schedule-segment attributes that are supported for personal health device communication.

Table E.4—Schedule-segment attributes

Attribute name	Attribute ID	Attribute type	Remark	Qualifiers
Schedule- Segment- Instance- Number	MDC_ATTR_SCHED_ SEG_INSTNO	InstNumber	The Instance-Number is the ID of a specific schedule-segment object instance. Each instance shall have a unique number assigned by the agent starting at 1. It is used by the manager to address a schedule segment. Note, the Instance-Number value of 0 shall be reserved for reporting no active schedule-segment in the schedule-store.	Mandatory Static
Schedule- Segment- Entry-Map	MDC_ATTR_SCHED_ SEG_MAP	ScheduleSegment EntryMap	This attribute defines the format and contents of one schedule entry. An entry has a conditional header containing information applicable to all elements in the entry. If the schedule-segment period, start time, and entry interval are used, then the header information shall not be used. The entry then contains one	Mandatory Dynamic
			or more elements, defined by the handle and an attribute value map defining the object attributes for each element in the schedule segment.	
Schedule- Segment- Period	MDC_ATTR_SCHED_ SEG_PERIOD	HighResRelativeTime	This attribute specifies the period of the schedule segments. If the schedule segment occurs only once, then the value shall be 0.	Mandatory Dynamic
Schedule- Segment- Entry-Interval	MDC_ATTR_SCHED_ SEG_ENTRY_INTERVAL	HighResRelativeTime	This attribute specifies the interval between each entry in schedule segments.	Conditional Dynamic

Table E.4—Schedule-segment attributes (continued)

Attribute name	Attribute ID	Attribute type	Remark	Qualifiers
Schedule- Segment- Person-Id	MDC_ATTR_SCHED_ SEG_PERSON_ID	PersonId	This standard supports devices that have simple support for data from multiple persons. A person ID is used to differentiate different persons. If the schedule segment is able to have schedule data for multiple persons, it shall set the schedsc-multi-person bit in the Sched-Store-Capab attribute. If this bit is set, all schedule-segment instances shall support the Schedule-Segment-Person-Id attribute. Otherwise, this attribute is not defined.	Conditional Dynamic
Schedule- Segment- Entry-Count	MDC_ATTR_SCHED_ SEG_ENTRY_CNT	INT-U32	This attribute gives the actual number of schedule-segment entries.	Optional Dynamic
Schedule- Segment- Label	MDC_ATTR_SCHED_ SEG_LABEL_STRING	OCTET STRING	This attribute is an application-dependent label in printable ASCII for the schedule segment to indicate its intended use and may be used for display purposes.	Optional Static
Schedule- Segment- LastUpdated- Abs-Time	MDC_ATTR_SCHED_ SEG_LAST_UPDATED_ ABS_TIME	AbsoluteTime	This attribute defines the time the schedule segment was last updated. If this attribute is used, neither the Schedule-Segment-LastUpdated-HiRes-Time nor Schedule-Segment-LastUpdated-BO-Time shall be used.	Conditional Dynamic
Schedule- Segment- LastUpdated- HiRes-Time	MDC_ATTR_SCHED_ SEG_LAST_UPDATED_ HIRES_TIME	HiResRelativeTime	This attribute defines the time the schedule segment was last updated. If this attribute is used, neither the Schedule-Segment-LastUpdated-Abs-Time nor Schedule-Segment-LastUpdated-BO-Time shall be used.	Conditional Dynamic
Schedule- Segment- LastUpdated- BO-Time	MDC_ATTR_SCHED_ SEG_LAST_UPDATED_ BO_TIME	BaseOffsetTime	This attribute defines the time the schedule segment was last updated. If this attribute is used, neither the Schedule-Segment-LastUpdated-Abs-Time nor Schedule-Segment-LastUpdated-HiRes-Time shall be used.	Conditional Dynamic

Table E.4—Schedule-segment attributes (continued)

Attribute name	Attribute ID	Attribute type	Remark	Qualifiers
Schedule- Segment- Reference- Abs-Time	MDC_ATTR_SCHED_ SEG_REF_ABS_TIME	AbsoluteTime	This attribute defines the time from which all entries within the schedule segment are referenced. This attribute shall be used, if the Schedule-Segment-Reference-BO-Time is not used.	Conditional Dynamic
Schedule- Segment- Reference- BO-Time	MDC_ATTR_SCHED_ SEG_REF_BO_TIME	BaseOffsetTime	This attribute defines the time from which all entries within the schedule segment are referenced. This attribute shall be used, if the Schedule-Segment-Reference-Abs-Time is not used.	Conditional Dynamic
Schedule- Segment- Start-Abs- Time	MDC_ATTR_SCHED_ SEG_START_ABS_TIME	AbsoluteTime	This attribute defines the time of the first execution of the schedule. If this attribute is used, the Schedule-Segment-Start-BO-Time shall not be used. If this attribute is used, the schedule segment becomes the active-schedule-segment in the schedule-store at the start time. A schedule segment shall not be active before the start time.	Optional Dynamic
Schedule- Segment- End-Abs- Time	MDC_ATTR_SCHED_ SEG_END_ABS_TIME	AbsoluteTime	This attribute defines the expiration of the schedule. If this attribute is used, the Schedule-Segment-End-BO-Time shall not be used. If this attribute is used and the end time is reached, the schedule-segment shall no longer be the active-schedule-segment in the schedule-store.	Optional Dynamic
Schedule- Segment- Start-BO- Time	MDC_ATTR_SCHED_ SEG_START_BO_TIME	BaseOffsetTime	This attribute defines the time of the first execution of the schedule. If this attribute is used, the Schedule-Segment-Start-Abs-Time shall not be used. If this attribute is used, the schedule-segment becomes the active-schedule-segment in the schedule-store at the start time. A schedule segment shall not be active before the start time.	Optional Dynamic

Table E.4—Schedule-segment attributes (continued)

Attribute name	Attribute ID	Attribute type	Remark	Qualifiers
Schedule- Segment- End- BO- Time	MDC_ATTR_SCHED_ SEG_END_BO_TIME	BaseOffsetTime	This attribute defines the expiration of the schedule. If this attribute is used, the Schedule-Segment-End-Abs-Time shall not be used. If this attribute is used and the end time is reached, the schedule-segment shall no longer be the active-schedule-segment in the schedule-store.	Optional Dynamic
Fixed- Schedule- Segment- Data	MDC_ATTR_SCHED_ SEG_FIXED_DATA	N/A The data is stored internal to the device and so this data type never occurs in any protocol definition directly.	This attribute defines the schedule data transferred as an array of entries in a format as specified in the Schedule-Segment-Entry-Map attribute. This is defined here as an opaque data structure without a defined data type. Note that this attribute is not directly accessible; it is only retrievable by the manager using the Trig-Schedule-Segment-Data-Xfer method.	Mandatory Dynamic
Schedule- Segment- Confirm- Timeout	MDC_ATTR_SCHED_ SEG_CONFIRM_TIMEOUT	RelativeTime	This informational timeout attribute defines the minimum time that the agent shall wait for a response message from the manager after issuing a Confirmed Event Report invoke message before timing out and transitioning to the Unassociated state. This is an informational attribute for the benefit of the manager. If this attribute is supplied, it shall match the actual timeout value that the agent uses for the Confirmed Event Report generated from the schedule-store object. This attribute is informational for the manager in the sense that the manager does not use this attribute in an actual implementation of the protocol (i.e., the manager does not timeout on an agent-generated Confirmed Event Report). However, the manager might wish to use this information to prioritize its handling of a "short" timeout agent over that of a "long" timeout agent.	Optional Dynamic

Health informatics—Personal health device communication—Part 10419: Device Specialization—Insulin Pump

Table E.4—Schedule-segment attributes (continued)

Attribute name	Attribute ID	Attribute type	Remark	Qualifiers
Schedule- Segment- Transfer- Timeout	MDC_ATTR_SCHED_SEG_ TRANSFER_TIMEOUT	RelativeTime	This timeout attribute defines the minimum time that the manager shall wait for the complete transfer of schedule-segment information. If the timeout expires prior to the reception of the complete schedule segment, the manager shall transition to the Unassociated state as described in subclause 8.9.5.6 in IEEE Std 11073-20601-2014.	Mandatory Dynamic

NOTE—The attribute qualifiers "static," "dynamic," and "observational" are omitted from Table E.4 since schedule segments are dynamic (the object itself may change during the lifetime of an association).

The Fixed-Schedule-Segment-Data attribute may hold a sizable amount of data, depending on the agent capabilities and the application. An agent may choose to restrict the maximum size of the Fixed-Schedule-Segment-Data attribute in a way that is aligned with the maximum transmission unit of the transport system. In order to support this type of behavior, a manager that supports Schedule-segments shall be able to support the transfer of Fixed-Schedule-Segment-Data.

If two or more schedule segments have the same Schedule-Segment-Start-Abs-Time or Schedule-Segment-Start-BO-Time, then the agent shall choose one segment as the Active-Schedule-Segment.

If the Schedule-Segment-End-Abs-Time or Schedule-Segment-End-BO-Time of an active schedule has occurred, then the agent shall update the value of the Active-Schedule-Segment.

E.2.4 Schedule-segment object methods

There are currently no schedule-segment object methods defined in this standard.

E.2.5 Schedule-segment object events

There are currently no schedule-segment object events defined in this standard.

E.2.6 Other schedule-segment services

There are currently no uses of the schedule-segment SET or GET services defined in this standard.

Health informatics—Personal health device communication—Part 10419: Device Specialization—Insulin Pump

Annex F

(normative)

Schedule class ASN.1 definitions

F.1 ACTION-method-related data types

```
-- SchedSegmSelection selects the Schedule-segments that are subject to the method.
SchedSegmSelection::= CHOICE {
                 all-sched-segments [1] INT-U16, -- if this type is chosen to select all schedule
                                                    -- segments the actual contents of the field is "do not
                                                    -- care" and shall be zero
                 sched-segm-id-list [2] SchedSegmIdList,
                                                    -- using this list requires that the manager already
                                                    -- knows the Schedule-Instance-Number attributes
                                  of
                                                    -- the Schedule-segments, e.g., from a previous
                                                    -- Get-Schedule-Segment-Info or
                                                    -- Get-Schedul-Segment-Id-List method call.
}
-- SchedSegmIdList selects Schedule-segments by ID.
SchedSegmIdList ::= SEQUENCE OF InstNumber
-- ScheduleSegmentInfoList returns the object attributes (except the Fixed-Schedule-Segment-Data) of all
-- selected Schedule-segment object instances in response to the Get-Schedule-Segment-Info or
-- Get-Schedule-Segment-Id-List Schedule-store method.
-- This is required by the manager to retrieve the dynamic information about the schedule segments.
ScheduleSegmentInfoList ::= SEQUENCE OF ScheduleSegmentInfo
ScheduleSegmentInfo ::= SEQUENCE {
                 sched-seg-inst-no InstNumber,
                 sched-seg-info AttributeList
}
```

F.2 Data types for new object attributes and object services

F.2.1 Schedule store and schedule-segment related data types

- -- The instance number list is used to distinguish list of class or object instances of the same type or object
- -- instances that are not directly manageable (used, e.g., as a list of identification attribute for Schedule-store
- -- objects).

InstNumberList ::= SEQUENCE OF InstNumber

Health informatics—Personal health device communication—Part 10419: Device Specialization—Insulin Pump

```
-- The Schedule-Store-Capab attribute defines specific static capabilities and properties of the Schedule-
-- store object
-- instance. The default value of this attribute is 0 (no bits set).
-- All unassigned "SchedStoreCapab" bit values are reserved for future expansion and shall be set to zero.
SchedStoreCapab ::=BITS-32 {
                 schedsc-var-no-of-segm(0),
                                                     -- indicates that the number of Schedule-segments
                                                     -- contained in this Schedule-store is dynamic and
                                                     -- may change
                                                     -- Schedule-segments in the SchedSegmSelection
                 schedsc-segm-id-list-select(3),
                                                     -- data type can be selected by defining a list of
                                                     -- segment identifiers
                 schedsc-epi-seg-entries(4),
                                                     -- Some/ all Schedule-segments contain
                                                     -- episodic/aperiodic entries and therefore have
                                                     -- to contain explicit time stamp information
                 schedsc-peri-seg-entries(5),
                                                     -- Some/all Schedule-segments contain periodically
                                                     -- sampled entries and therefore the Schedule
                                                     -- segment or Schedule-store shall support the
                                                      -- Schedule-Sample-Period attribute
                                                      -- The Schedule-store supports Schedule-segment for
                 schedsc-multi-person(12)
                                                      -- more than one person
                 schedsc-get-segm-info-sup(13),
                                                      -- The Get-Schedule-Segment-Info method is
                                                      -- supported.
                 schedsc-get-segm-id-list-sup(14),
                                                     -- The Get-Schedule-Segment-Id-List method is
                                                      -- supported.
}
-- All entries in the segment shall follow the format defined by this attribute. First, the mandatory header
-- shall follow the description in sched-segm-entry-header. This allows each entry in the segment to be
-- preceded by a header (e.g., periodicity, start time and duration information) that is applicable to all
-- elements in an entry. Next, the elements shall follow the format and order described in
-- sched-segm-entry-elem-list. An element typically represents a setting. For each element, the stored
-- data is defined in the form of an attribute value map, in the same way as metric objects.
ScheduleSegmentEntryMap ::= SEQUENCE {
                 sched-segm-entry-header SchedSegmEntryHeader,
                                                                       -- defines start time type in front
                                                                        -- of each entry
                                                                        -- (ScheduleSegmentEntryHeader)
                 sched-segm-entry-elem-list SchedSegmEntryElemList
}
-- The following format defines data items that are in front of each schedule-segment entry. The
-- conditional header allows definition of data items that are common to all elements in the entry. If the
-- Schedule-Segment-Entry-Interval is used, then the header shall not be used. The schedule segment entry
-- event report shall begin with data from the first element.
SchedSegmEntryHeader ::= BITS-16 {
                 sched-seg-hdr-absolute-time(0),
                                                             -- entry preceded by absolute time
                                                             -- (data type AbsoluteTime)
                 sched-seg-hdr-relative-time(1),
                                                             -- entry preceded by relative time
                                                             -- (data type RelativeTime)
```

Health informatics—Personal health device communication—Part 10419: Device Specialization—Insulin Pump

```
sched-seg-hdr-hires-relative-time(2),
                                                            -- entry preceded by high resolution relative
                                                            -- time (data type HighResRelativeTime)
                 sched-seg-hdr-bo-time(3)
                                                            -- entry preceded by base offset time
                                                            -- (data type BaseOffsetTime)
                                                            -- option (0) and option (3) are mutually
                                                            -- exclusive
}
-- All unassigned "SchedSegmEntryStartTime" bit values are reserved for future expansion and shall be set
-- to zero. If any bits are set to one beyond the expected bits (e.g., a new bit was added in a later version),
-- the data shall not be retrieved since the offset to the first data element cannot be calculated.
-- If the Schedule-segment is periodic, the SchedSegmElemStartTime shall be relative time or
-- high-resolution relative time relative to the Schedule-Segment Reference Time attribute.
SchedSegmEntryElemList ::= SEQUENCE OF SchedSegmEntryElem
-- SchedSegmEntryElem shall reference a metric object instance in the agent configuration
-- using its handle value. This referenced object shall exist in the agent configuration.
SchedSegmEntryElem ::= SEQUENCE {
                 handle
                                   HANDLE,
                                                    -- handle of referenced object
                 attr-val-map
                                   AttrValMap
                                                    -- attribute value map describing the data
}
-- Request to start the transfer of the specified segment
TrigSchedSegmDataXferReq ::= SEQUENCE {
                 sched-seg-inst-no InstNumber
}
TrigSchedSegmDataXferRsp ::= SEQUENCE {
                 schedseg-inst-no
                                           InstNumber,
                 trig-sched-segm-xfer-rsp TrigSchedSegmXferRsp
}
-- All unassigned "TrigSchedSegmXferRsp" values are reserved for future expansion and shall not be used.
TrigSchedSegmXferRsp ::= INT-U16 {
                                                    -- Agent will start transfer of schedule segment
                 tschsxr-successful(0),
                 tschsxr-fail-no-such-segment(1), -- schedule-segment ID not found
                 tschsxr-fail-segm-empty(3),
                                                    -- the schedule segment being requested is empty
                 tschsxr-fail-not-otherwise-specified(512)
}
-- ScheduleSegmentDataEvent
-- Notes:
    the agent shall transfer all segment entries in order, first entry first (first in first out).
ScheduleSegmentDataEvent ::= SEQUENCE {
```

Health informatics—Personal health device communication—Part 10419: Device Specialization—Insulin Pump

```
sched-segm-data-event-descr SchedSegmDataEventDescr,
               sched-segm-data-event-entries OCTET STRING
                                                                    -- contains the specified segment
                                                                    -- entries in an opaque data structure.
                                                                    -- Only complete entries shall be
                                                                    -- included in this field.
}
ScheduleSegmentDataResult ::= SEQUENCE {
                 sched-segm-data-event-descr SchedSegmDataEventDescr
}
-- The Schedule Segment Data Event Descriptor defines which entries of the Schedule Segment Data are
-- communicated in the Event message.
SchedSegmDataEventDescr ::= SEQUENCE {
                 sched-segm-instance
                                         InstNumber, -- instance number of segment being transferred
                 sched-segm-evt-entry-index INT-U32, -- array index of the first entry in this event
                 sched-segm-evt-entry-count INT-U32, -- count of entries in this event
                 sched-segm-evt-status SchedSegmEvtStatus
-- All unassigned "SchedSegmEvtStatus" bit values are reserved for future expansion and shall be set to
-- zero.
SchedSegmEvtStatus ::= BITS-16 {
                 schsevtsta-first-entry(0),
                                                    -- this event contains the first segment entry
                 schsevtsta-last-entry(1),
                                                    -- this event contains the last segment entry (both first
                                                    -- and last bits can be set if all entries fit in one event)
                                                    -- transfer aborted by agent (manager shall reply
                 schsevtsta-agent-abort(4),
                                                    -- with the same status)
                 schsevtsta-manager-confirm(8),
                                                    -- set in reply if segment was received correctly (if
                                                    -- not set in reply, agent shall stop the segment
                                                    -- transfer and respond with an error (roer) code
                                                    -- of protocol-violation).
                                                    -- sent in reply by manager (agent shall stop sending
                 schsevtsta-manager-abort(12)
                                                     -- messages)
}
```

F.3 Data protocol definitions

F.3.1 GET service

```
-- For GET requests defined in this standard, obj-handle shall either be 0 to represent the MDS object or a
```

- -- handle representing a PM-store object or a handle representing a Schedule-store object.
- -- The attribute-id-list shall be left empty if the purpose is to query for all attributes of the MDS or PM-
- -- store object or Schedule-store object.
- -- Alternatively, specific attributes of an object may be queried by listing the desired
- -- Attribute IDs found in IEEE Std 11073-20601-2014 or Table E.1.

- -- For GET responses defined in this standard, obj-handle shall match the one in the corresponding request.
- -- The attribute-list contains all the requested attributes using the variable format.
- -- If a requested attribute ID does not exist within the MDS object, it shall not
- -- be returned in the attribute-list.

}

GetResultSimple ::= SEQUENCE {

obj-handle HANDLE, attribute-list AttributeList

Health informatics—Personal health device communication—Part 10419: Device Specialization—Insulin Pump

Annex G

(informative)

The schedule-store concept

G.1 General

The schedule-store concept provides storage capabilities for a series of scheduled segments stored in the agent. The information is organized in a hierarchical object model with capability to allow schedules to be stored with a structure appropriate to the nature of the data.

At a top level, the schedule-store object is the primary access point for all information about the schedules. An agent supporting schedules may instantiate one or more schedule-store objects. The schedule-store object is part of the device configuration and is directly accessed with the object access services defined in this standard.

Each schedule-store may contain 0, 1, or more schedule segments that are the actual data container objects. The number of schedule segments may change as a result of the operation of the agent. In other words, the agent may create new schedule segments based on time intervals, size of the stored data, or even manual controls of the user.

The schedule-store concept provides an information model with a two-level hierarchy with multiple schedule-segment objects within multiple schedule-store objects.

Typical-use cases for using multiple schedule-stores include the following:

— If the agent stores schedules with different characteristics (e.g., one-time schedule versus recurring schedule), separate schedule-store objects are used to define optimized data types for the stored data and thus conserve memory for the stored data.

Typical-use cases for using multiple schedule segments include the following:

— If the agent needs to structure the stored data in a more hierarchical form, it can use multiple instances of schedule-store objects with multiple instances of schedule-segment objects to model this hierarchy (e.g., use the schedule store to represent a particular type of schedule, and then use the schedule-segment to model individual schedules within this schedule type).

For the actual data storage, the attribute value map concept as used for the metric attributes is used here. A special mapping attribute allows defining the structure of the binary stored data and avoiding any overhead for identification, length fields, etc., in the actual stored and transmitted binary data. This assumes that stored data are essentially a large array of equally formatted data.

The transfer of the stored data is triggered by the manager after inspecting the information in the schedule-store objects. The manager can select the data episodes to transfer. The actual transfer is then done by the agent using acknowledged event report messages. The agent is expected to fill the ScheduleSegmentDataEvent data structure to the available maximum size.

Health informatics—Personal health device communication—Part 10419: Device Specialization—Insulin Pump

G.2 Schedule-store object hierarchy

G.2.1 General

A schedule-store consists of the following four key parts:

- Schedule-store—This object is at the top level, and it contains attributes about the storage object as well as zero or more schedule segments.
- Schedule segment—This object contains attributes that describe the segment as well as zero or more entries.
- **Entry**—Each entry holds an optional entry header and one or more elements.
- **Element**—Each element holds data from one or more metric settings.

Figure G.1 shows an example layout of these four parts, which are further described in the remainder of this annex.

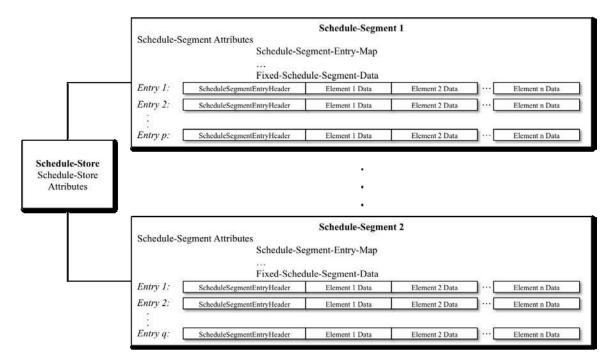


Figure G.1—Schedule-store with two schedule segments, fixed-segment-data within the segments

G.2.2 Schedule-store object

Support for the schedule-store object is optional. Only agents that wish to store schedules need to provide support for the schedule-store object, attributes, methods, and associated events. A manager becomes aware of all supported schedule-store objects as part of the agent configuration. Attributes of the schedule-store describe common characteristics of the schedule data (e.g., if the schedule is recurring or single instance).

An agent may provide more than one schedule-store. Multiple stores are used to represent schedule data with different formats or with different characteristics or to group schedule data into different logical

Health informatics—Personal health device communication—Part 10419: Device Specialization—Insulin Pump

buckets. The schedule-store object is also the access point for all methods related to the schedules (specifically, the manager retrieving the schedule data).

The agent controls the number of schedule-segments that exist in a schedule-store. A schedule-store may have zero elements when there are no data present. When data are present, the schedule-store has one or more schedule-segment objects. As the number of schedule segments is dynamic, the schedule-segment objects are not part of the agent configuration. Instead, the schedule-store object contains the information about available schedule-segments in the form of schedule-store attributes that are queried using the GET service.

G.2.3 Schedule-segment object

The basic format of the schedule-segment segment data is shown in Figure G.2.

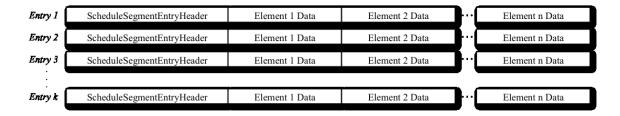


Figure G.2—Schedule-segment data format

The schedule segment contains k entries. The format of an entry is defined by the Schedule-Segment-Entry-Map attribute of the schedule segment. An entry represents the schedule at one particular point in time. Each entry is preceded by a conditional header (e.g., containing a start time) common to all elements of the entry. The entry then contains n elements; the format of each element is defined by an attribute value map. An entry typically contains a setting (e.g., numeric and enumeration). The resulting data structure does not contain any attribute ID or length fields and is, therefore, extremely compact.

The schedule segment typically represents one schedule episode. This episode has a time context (e.g., schedule data stored in this segment is from 12:00 to 15:00), some related attributes, and a storage array that contains the setting data for that episode contained in the Fixed-Schedule-Segment-Data attribute (see Figure G.3).

Health informatics—Personal health device communication—Part 10419: Device Specialization—Insulin Pump

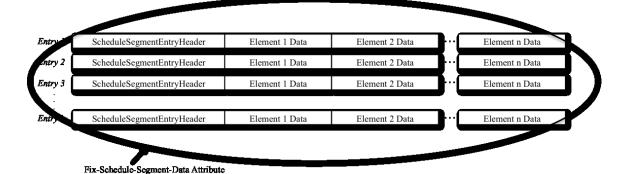



Figure G.3—Fixed-segment-data attribute containing actual stored data

A schedule-store may contain zero or more schedule segments (i.e., zero if data are not yet stored; one or more depending on the stored episodes and the capabilities of the agent).

For example, a running watch's schedule segment could contain the schedule data about one training cycle (e.g., a 5-mi run starting at 12:00). The device is able to store multiple schedule segments (i.e., multiple such training cycles).

G.2.4 Schedule-segment entry (within the fixed segment data)

The Fixed-Schedule-Segment-Data attribute contains both entries and elements. The entry items are depicted as rows in Figure G.4. All entries within a schedule segment have the same data structure that is defined by the Schedule-Segment-Entry-Map. This is very comparable to the Attribute-Value-Map that is defined for the metric objects. However, it allows grouping multiple settings in an entry item.

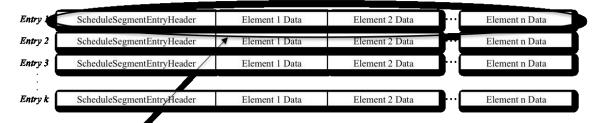

Fix-Schedule-Segment-Data Attribute

Figure G.4—Entry (array element in fixed-segment-data)

The Schedule-Segment-Entry-Map attribute defines the list of settings stored in one entry. For each setting, the list of attributes that are stored is also defined. Additionally, a common header (e.g., to include a common start time) that applies to the complete entry is optionally defined.

G.2.5 Schedule-segment entry element

An element contains the binary representation of the defined attributes of one metric object (see Figure G.5). The SchedSegmEntryElem (see F.1) within the Schedule-Segment-Entry-Map defines each entry element.

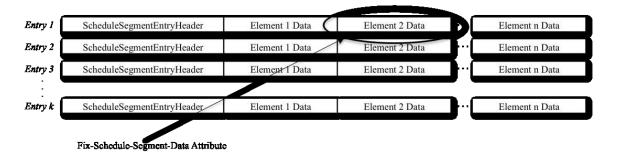


Figure G.5—Element: the set of attributes for one setting

Health informatics—Personal health device communication—Part 10419: Device Specialization—Insulin Pump

Annex H

(informative)

Scedule communication model

H.1 Operating procedure

H.1.1 General

The communication of health data and status information about the agent occurs during the Operating state.

H.1.1.1 Entry conditions

The agent and manager enter the Operating state either when the agent's configuration is already known by the manager or after the agent has communicated an acceptable configuration to the manager.

H.1.1.2 Normal procedures

H.1.1.2.1 General

Subclause H.1.1.2.2 describes procedures that can occur when in the Operating state.

H.1.1.2.2 Schedule-stored metric data transfer

H.1.1.2.2.1 General

When an agent implements one or more schedule-store objects, the agent reports about the existence of the schedule-store object during the configuration phase. The manager uses this information to query the schedule-store object(s) of the agent. The interactions between the manager and agent when retrieving the information in the schedule-store(s) is described in H.1.1.2.2.1.1.

Health informatics—Personal health device communication—Part 10419: Device Specialization—Insulin Pump

H.1.1.2.2.1.1 Schedule-stored metric data transmission

a) Retrieving the schedule-store attributes. When the agent and manager are in the Operating state, the manager can inspect the configuration negotiated with the agent to determine the number of schedule-store objects in the agent. The manager may query each schedule-store to determine the number of schedule-segments that exist within the schedule-store. Figure H.1 shows the sequence diagram of this operation. The manager sends a Get command to the agent requesting attribute information from a particular schedule-store. The manager uses the handle number to reference the desired schedule-store. The handle value is placed in the obj-handle field of the message (F.3) and is not present in the attribute-id list of the request or the attribute-list of the response. The attribute-id-list shall be left empty to query for all attributes of the schedule-store object. Alternatively, specific attributes of an object may be queried by listing the desired attribute IDs found in Table E.1. It is not required for an agent to support this capability. If this capability is not implemented, then the agent shall respond with an error (roer) message with an error-value of not-allowed-by-object.

The manager can inspect the attributes to learn about how the store is configured. For instance, the Schedule-Store-Capab describes the capabilities of the store, and Number-Of-Schedule-Segments defines how many segments are present in the Schedule-store. See Table E.1 for the full list of attributes and their definitions.

If the agent supports multiple schedule-store object instances, a Get request is required for each schedule-store.

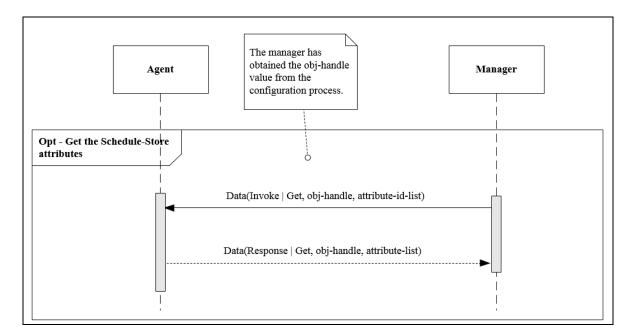


Figure H.1—Retrieving a schedule-store's attributes

Health informatics—Personal health device communication—Part 10419: Device Specialization—Insulin Pump

b) Retrieving the schedule-segment information. The manager retrieves information on the segments in a schedule-store by sending an ACTION.Get-Schedule-Segment-Info or ACTION.Get-Schedule-Segment-Id-List command to the specific schedule-store (see Figure H.2) with a request to return information from all segments, a particular list of segments, or any segments within a given time range. If there is no segment in any of these three cases, the agent responds with an empty list. The agent shall support the first selection criteria and may provide support for the second and third selection criteria. The manager is able to determine whether the agent provides support by inspecting schedsc-abs-time-select in the Schedule-Store-Capab attribute of the schedule-store information retrieved earlier.

The agent responds to the ACTION.Get-Schedule-Segment-Info command with a list of segment numbers followed by the full attribute list for each of the segments. The agent responds to the ACTION.Get-Schedule-Segment-Id-List command with a list of the instance numbers.

c) If the manager invokes one of the optional Get-Schedule-Segment-Info or Get-Schedule-Segment-Id-List methods but the agent does not support the particular optional action (list of segments), then the agent shall respond with a roer DataApdu with a RoerErrorValue of "not-allowed-by-object".

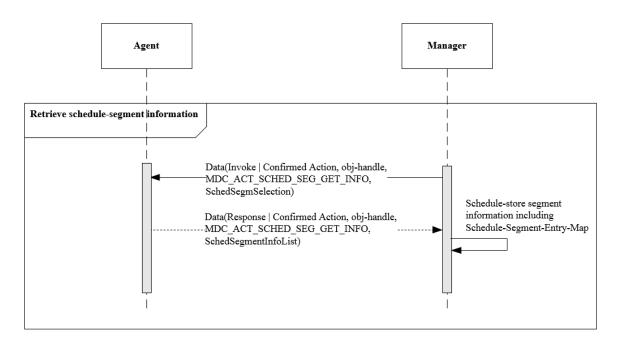


Figure H.2—Retrieve schedule-segment information using Get-Schedule-Segment-Info

d) **Transfer schedule-segment content.** The manager retrieves specific schedule-segments by using the Trig-Schedule-Segment-Data-Xfer ACTION method to initiate the data transfer (see Figure H.3). In the first step, the manager sends the ACTION method to the agent with the handle of the schedule-store object to access. The argument to this ACTION method is the instance number of the segment to transfer.

The agent shall decide whether the request can be honored. It checks for a valid segment number, available segment data (e.g., the segment could be in the update process), or any other error conditions. If there is an error, the agent shall return an appropriate error code in the response and ignore the transmit request. Otherwise, the agent shall send a tsxr-successful response code to indicate that it has received the request and it can be honored.

ISO/IEEE 11073-10419:2016(E)

IEEE Std 11073-10419-2015

Health informatics—Personal health device communication—Part 10419: Device Specialization—Insulin Pump

The manager may send the Trig-Schedule-Segment-Data-Xfer ACTION invoke message at any time

The agent shall send confirmed Schedule-Segment-Data-Event event reports until all entries in the schedule-segment are sent to the manager or the transfer is aborted by either the schsevtsta-agent-abort or schsevtsta-manager-abort bits described later. The agent fills in the Schedule-SegmentDataEvent structure with information about the segment being sent. The agent shall transfer all segment entries in order, first entry first (FIFO: first-in, first-out). The agent informs the manager of the schedule-store handle and uses the SchedSegmDataEventDescr to describe the segment number being transferred, the entry number of the first entry in the sched-segm-data-event-entries field, the number of entries in the message, and current status information. The agent shall always set any schsevtsta-manager-* bits to 0. If the message contains the first entry and/or the last entry of the data entries, then the agent shall set the schsevtsta-first-entry and/or schsevtsta-last-entry bits, respectively. If the agent wishes to abort the transfer, it shall set the schsevtsta-agent-abort bit to 1.

When transferring a segment, the agent uses the sched-segm-data-event-entries field to send all the entries. The agent shall start with the first entry collected, followed by the next entry, and so on. The agent should pack as many entries as possible into the event structure to optimize the transmissions. Each entry shall be formatted according to the structure defined in the schedule-segment Schedule-Segment-Entry-Map.

When the manager receives an event report, it shall reply with a ScheduleSegmentDataResult response that shall contain the same schedule-store-handle, sched-segm-instance number, sched-segm-evt-entry-index, and sched-segm-evt-entry-count. In the sched-segm-evt-status, the manager shall set the schsevtsta-manager-confirm bit.

If the agent sets the schsevtsta-agent-abort bit, then the manager shall confirm the agent abort by setting the same bit. If the manager wishes to abort the exchange, it shall set the schsevtsta-manager-abort bit.

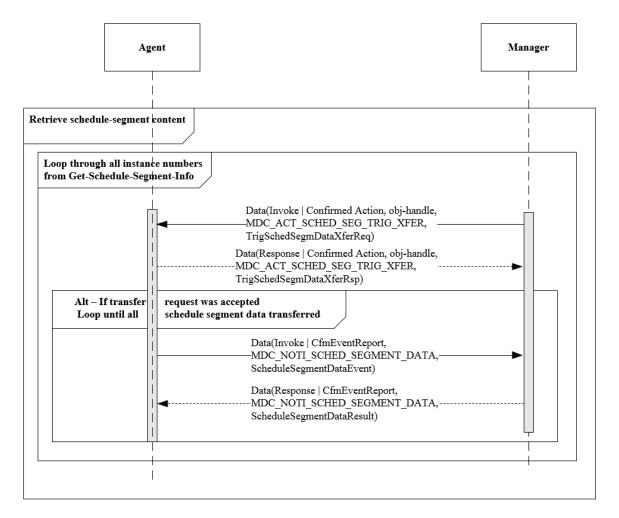


Figure H.3—Retrieve schedule-segment content

Health informatics—Personal health device communication—Part 10419: Device Specialization—Insulin Pump

Annex I

(informative)

Protocol data unit examples

I.1 General

This annex shows binary examples of messages exchanged between an insulin pump agent and a manager. Three different scenarios containing the association and configuration information exchanges are presented in I.2 and I.2.4. The first scenario illustrates the case when the agent intends to operate using an extended configuration. The manager does not have the configuration declared by the agent from a prior association. The second illustrates the agent presenting the same extended configuration to the manager, and the manager does have the configuration from the previously transferred configuration exchange. Finally, the agent presents a standard configuration to the manager, and the manager has the configuration because the manager has been preprogrammed with this configuration.

I.2 Association information exchange

I.2.1 General

When the transport connection is established between the manager and the agent, they both enter the Unassociated state. When the agent sends an association request, both manager and agent enter the Associating state.

I.2.2 Extended configuration

I.2.2.1 General

In this exchange, the agent sends an association request intending to use an extended configuration during measurement transfer. However, the manager does not have this configuration.

I.2.2.2 Association request

The insulin pump agent sends the following message to the manager. The agent intends to associate using an extended configuration.

0xE2	0x00			APDU CHOICE Type (AarqApdu)
0x00	0x32			CHOICE.length = 50
0x80	0x00	$0 \times 0 0$	0x00	assoc-version
0x00	0x01	$0 \times 0 0$	0x2A	data-proto-list.count = $1 \mid length = 42$
0x50	0x79			data-proto-id = 20601
0x00	0x26			data-proto-info length = 38
0x20	0x00	$0 \times 0 0$	0x00	protocolVersion
0x80	0x00			encoding rules = MDER
0x80	0x00	0x00	0x00	nomenclatureVersion

Health informatics—Personal health device communication—Part 10419: Device Specialization—Insulin Pump

0x00	0x0	$0 \times 0 0$	0x00	functionalUnits – no Test Association capabilities
0x00	08x0	0x00	0x00	systemType = sys-type-agent
0x00	0x08			system-id length = 8 and value (manufacturer-and device-specific)
0x31	0x32	0x33	0x34	
0x35	0x36	0x37	0x38	
0x40	0x00			dev-config-id – extended configuration
0x00	0x0			data-req-mode-flags
0x01	0x00			data-req-init-agent-count = 1 data-req-init-manager-count = 0
0x00	0x00	0x00	0x00	optionList.count = $0 \mid optionList.length = 0$

I.2.2.3 Association response

A manager responds to the agent that it can associate but does not have the insulin pump extended configuration (i.e., there is the need for the agent to send its configuration).

0xE3 0x	x00		APDU CHOICE Type (AareApdu)
0x00 0x	x2C		CHOICE.length = 44
0x00 0x	x03		result = accepted-unknown-config
0x50 0x	x79		data-proto-id = 20601
0x00 0x	x26		data-proto-info length = 38
0x20 0x	x00 0x	0x0	protocolVersion
0x80 0x	x00		encoding rules = MDER
0x80 0x	x00 0x	0x0	nomenclatureVersion
0x00 0x	x00 0x	0x0	functionalUnits – normal Association
0x80 0x	x00 0x	0x0	systemType = sys-type-manager
0x00 0x	80x		system-id length = 8 and value (manufacturer-and device-specific)
0x38 0x	x37 0x	36 0x3	
0x34 0x	x33 0x	32 0x3	
0x00 0x	x00		Manager's response to config-id is always 0
0x00 0x	x00		data-req-mode-flags
0x00 0x	x00		data-req-init-agent-count = 0 data-req-init-manager-count=0
0x00 0x	x00 0x	0x0	optionList.count = $0 \mid optionList.length = 0$

I.2.3 Previously known extended configuration

I.2.3.1 General

This exchange illustrates a transaction that takes place after a session beginning with an exchange like I.2.2 has occurred.

I.2.3.2 Association request

The insulin pump agent sends the following message to the manager. The agent intends to associate using an extended configuration.

0xE2	0x00			APDU CHOICE Type (AarqApdu)
0x00	0x32			CHOICE.length = 50
0x80	0x00	0x00	0x00	assoc-version

Health informatics—Personal health device communication—Part 10419: Device Specialization—Insulin Pump

0x00	0x01	0x00	0x2A	data-proto-list.count = 1 length = 42
0x50	0x79			data-proto-id = 20601
0x00	0x26			data-proto-info length = 38
0x20	0x00	0x00	0x00	protocolVersion
0x80	0x00			encoding rules = MDER
0x80	0x00	0x00	0x00	nomenclatureVersion
0x00	0x00	0x00	0x00	functionalUnits – no Test Association capabilities
0x00	0x80	0x00	0x00	systemType = sys-type-agent
0x00	0x08			system-id length = 8 and value (manufacturer-and device-specific)
0x31	0x32	0x33	0x34	
0x35	0x36	0x37	0x38	
0×40	0x00			dev-config-id – extended configuration
0x00	0x00			data-req-mode-flags
0x01	0x00			data-req-init-agent-count = 1 data-req-init-manager-count=0
0x00	0x00	0x00	0x00	optionList.count = 0 optionList.length = 0

I.2.3.3 Association response

A manager responds to the agent that it can associate with, recognizes, and accepts and has the insulin pump's extended configuration (i.e., there is no need for the agent to send its configuration).

0xE3	0x00			APDU CHOICE Type (AareApdu)
0x00	0x2C			CHOICE.length = 44
0x00	0x0			result = accepted
0x50	0x79			data-proto-id = 20601
0x00	0x26			data-proto-info length = 38
0x20	0x00	0x00	0x00	protocolVersion
0x80	0x00			encoding rules = MDER
0x80	0x00	0x00	0x00	nomenclatureVersion
0x00	0x00	0x00	0x00	functionalUnits – normal Association
0x80	0x00	0x00	0x00	systemType = sys-type-manager
0x00	0x08			system-id length = 8 and value (manufacturer-and device-specific)
0x38	0x37	0x36	0x35	
0x34	0x33	0x32	0x31	
0x00	0x00			Manager's response to config-id is always 0
0x00	0x00			data-req-mode-flags
0x00	$0 \times 0 0$			data-req-init-agent-count = 0 data-req-init-manager-count=0
0x00	0x00	0x00	0x00	optionList.count = 0 optionList.length = 0

I.2.4 Standard configuration

I.2.4.1 General

This transaction would occur if an agent presents an association request incorporating the dev-config-id corresponding to a standard configuration. The manager has the configuration because it has been programmed with this configuration according to the information presented in this standard.

I.2.4.2 Association request

The insuin pump agent sends the following message to the manager. The agent intends to associate using a standard configuration. The agent is willing to enter into a test association as defined in Clause 9.

0xE2	0x00			APDU CHOICE Type (AarqApdu)
0x00	0x32			CHOICE.length = 50
08x0	0x00	0x00	0x00	assoc-version
0x00	0x01	0x00	0x2A	data-proto-list.count = 1 length = 42
0x50	0x79			data-proto-id = 20601
0x00	0x26			data-proto-info length = 38
0x20	0x00	0x00	0x00	protocolVersion
0x80	0x00			encoding rules
0x80	$0 \times 0 0$	0x00	0×00	Nomenclature version
0x00	0x00	$0 \times 0 0$	0x00	Functional units
0x00	0x80	$0 \times 0 0$	0x00	systemType = sys-type-agent
0x00	0x08			system-id length = 8 and value (manufacturer-and device-specific)
	0x32			
	0x36	0x37	0x38	
	0x6C			dev-config-id: 1900
0x00				data-req-mode-flags
0x01				data-req-init-agent-count = 1 data-req-init-manager-count=0
$0 \times 0 0$	$0 \times 0 0$	$0 \times 0 0$	0x00	optionList.count = $0 \mid optionList.length = 0$

I.2.4.3 Association response

A manager responds to the agent that it can associate with, recognizes, and accepts and has the insulin pump standard configuration (i.e., there is no need for the agent to send its configuration). The manager does not start a test association.

0xE3	0x00			APDU CHOICE Type (AareApdu)
$0 \times 0 0$	0x2C			CHOICE.length = 44
0x00	$0 \times 0 0$			Result = accept
0x50	0x79			data-proto-id = 20601
0x00	0x26			data-proto-info length = 38
0x20	0x00	0x00	0x00	protocolVersion
0x80	0x00			encoding rules = MDER
0x80	0x00	0x00	0x00	nomenclatureVersion
0x00	0x00	0x00	0x00	functionalUnits
0x80	0x00	0x00	0x00	systemType = sys-type-manager
0x00	0x08			system-id length = 8 and value (manufacturer-and device-specific)
0x38	0x37	0x36	0x35	
0x34	0x33	0x32	0x31	
0×00	0x00			Manager's response to config-id is always 0
0×00	0x00			data-req-mode-flags
0x00	0x00			data-req-init-agent-count = 0 data-req-init-manager-count=0
0x00	0x00	0x00	0x00	optionList.count = 0 optionList.length = 0

Health informatics—Personal health device communication—Part 10419: Device Specialization—Insulin Pump

I.3 Configuration information exchange

I.3.1 General

If the association is not rejected or aborted, the agent and manager transition from the Associating state into one of two states. If the manager's AssociateResult code is accepted, the agent and manager enter the operating state. If the manager's AssociateResult code is accepted-unknown-config, the agent and manager enter the Configuring state.

I.3.2 Extended configuration

I.3.2.1 General

This exchange takes place when the manager returns the AssociateResult code of accepted-unknown-config. The agent presents a description of its configuration corresponding to the dev-config-id it presented in the association request.

I.3.2.2 Remote operation invoke event report configuration

The insulin pump agent sends the description of its extended configuration. It does this by sending a confirmed event report of type MDC_NOTI_CONFIG.

The insulin pump agent sends the description of its extended configuration. It does this by sending a confirmed event report of type MDC_NOTI_CONFIG.

0×01 0×32 CHOICE.length = 306 0×01 0×30 OCTET STRING.length = 304 0×12 0×34 invoke-id = 0×1234 (start of DataApdu. MDER encoded.) 0×01 0×01 CHOICE(Remote Operation Invoke Confirmed Event Report) 0×01 $0 \times 2A$ CHOICE.length = 298 0×00 0×00 obj-handle = 0 (MDS object) 0×01 0×01 event-time = $0 \times FFFFFFFF$ 0×00 0×10 event-time = 0×00 0×01 0×20 event-time = 0×00 0×01 0×00 config-report-id 0×00 0×00 config-obi-list.count = 0×00	0xE7	0x00			APDU CHOICE Type (PrstApdu)
$\begin{array}{cccccccccccccccccccccccccccccccccccc$	0x01	0x32			CHOICE.length = 306
0×01 0×01 CHOICE(Remote Operation Invoke Confirmed Event Report) 0×01 $0 \times 2A$ CHOICE.length = 298 0×00 0×00 obj-handle = 0 (MDS object) $0 \times FF$ $0 \times FF$ $0 \times FF$ event-time = $0 \times FFFFFFFFF$ 0×00 $0 \times 1C$ event-type = MDC_NOTI_CONFIG 0×01 0×20 event-info.length = 288 (start of ConfigReport) 0×40 0×00 config-report-id	0x01	0x30			OCTET STRING.length = 304
$ \begin{array}{cccccccccccccccccccccccccccccccccccc$	0x12	0x34			invoke-id = 0x1234 (start of DataApdu. MDER encoded.)
0×00 0×00 obj-handle = 0 (MDS object) $0 \times FF$ <td>0x01</td> <td>0x01</td> <td></td> <td></td> <td>CHOICE(Remote Operation Invoke Confirmed Event Report)</td>	0x01	0x01			CHOICE(Remote Operation Invoke Confirmed Event Report)
$0 \times FF$ $0 \times F$	0x01	0x2A			CHOICE.length = 298
$0 \times 0D$ $0 \times 1C$ event-type = MDC_NOTI_CONFIG 0×01 0×20 event-info.length = 288 (start of ConfigReport) 0×40 0×00 config-report-id	0x00	0x00			obj-handle = 0 (MDS object)
$0 \times 01 0 \times 20$ event-info.length = 288 (start of ConfigReport) config-report-id	0xFF	$0 \times FF$	$0 \times FF$	0×FF	event-time = 0xFFFFFFF
0x40 0x00 config-report-id	0x0D	0x1C			event-type = MDC_NOTI_CONFIG
÷ .	0x01	0x20			event-info.length = 288 (start of ConfigReport)
0v00 0v06 applied by list count - 6 Massurament objects will be "ennounced"?	0x40	0x00			config-report-id
config-obj-nst.count = 6 Measurement objects will be announced	0x00	0x06			config-obj-list.count = 6 Measurement objects will be "announced"
$0 \times 01 0 \times 1A$ config-obj-list.length = 282	0x01	$0 \times 1 A$			config-obj-list.length = 282
$0 \times 00 0 \times 06$ obj-class = MDC_MOC_VMO_METRIC_NU	0x00	0x06			obj-class = MDC_MOC_VMO_METRIC_NU
$0 \times 00 0 \times 01$ obj-handle = 1 ($\rightarrow 1^{st}$ object is current basal rate setting)	0x00	0x01			obj-handle = 1 (\rightarrow 1 st object is current basal rate setting)
$0 \times 00 0 \times 04$ attributes.count = 4	0x00	0×04			attributes.count = 4
$0 \times 00 0 \times 24$ attributes.length = 36	0x00	0x24			attributes.length = 36
$0 \times 09 0 \times 2F$ attribute-id = MDC_ATTR_ID_TYPE	0x09	0x2F			$attribute-id = MDC_ATTR_ID_TYPE$
$0 \times 00 0 \times 04$ attribute-value.length = 4	0x00	0×04			attribute-value.length = 4
0x00 0x80 0xE4 0x00 MDC_PART_PHD_DM MDC_INS_BASAL_RATE	0x00	0x80	$0 \times E4$	0x00	MDC_PART_PHD_DM MDC_INS_BASAL_RATE
$0 \times 0 = 0 \times 46$ attribute-id = MDC_ATTR_METRIC_SPEC_SMALL	0x0A	0x46			attribute-id = MDC_ATTR_METRIC_SPEC_SMALL
$0 \times 00 0 \times 02$ attribute-value.length = 2	0x00	0x02			attribute-value.length = 2
0xF0 0x40 intermittent, stored data, upd & msmt aperiodic, agent init	0xF0	0×40			intermittent, stored data, upd & msmt aperiodic, agent init
$0 \times 09 0 \times 96$ attribute-id = MDC_ATTR_UNIT_CODE	0x09	0x96			attribute-id = MDC_ATTR_UNIT_CODE

$0 \times 00 0 \times 02$ attribute-value.length = 2	
0x16 0x40 MDC_DIM_INTL_UNIT_PER_HR	
$0 \times 0 = 0 \times 55$ attribute-id = MDC_ATTR_ATTRIBUTE_V	AL_MAP
$0 \times 00 0 \times 0C$ attribute-value.length = 12	
$0 \times 00 0 \times 02$ AttrValMap.count = 2	
$0 \times 00 0 \times 08$ AttrValMap.length = 8	
0x0A 0x56 0x00 0x04 MDC_ATTR_NU_VAL_OBS_SIMP, 4	
0x0A 0x82 0x00 0x08 MDC_ATTR_TIME_STAMP_BO, 8	
0x00 0x06 obj-class = MDC_MOC_VMO_METRIC_NU	J
$0 \times 00 0 \times 02$ obj-handle = 2 (\Rightarrow 2 nd measurement is bolu	is delivered)
$0 \times 00 0 \times 04$ attributes.count = 4	ŕ
$0 \times 00 0 \times 24$ attributes.length = 36	
0x09 0x2F attribute-id = MDC_ATTR_ID_TYPE	
$0 \times 00 0 \times 04$ attribute-value.length = 4	
0x00 0x80 0xE4 0x0C MDC_PART_PHD_DM MDC_INS_BOLUS	S
$0 \times 0 = 0 \times 46$ attribute-id = MDC_ATTR_METRIC_SPEC_	
$0 \times 00 0 \times 02$ attribute-value.length = 2	_51/11/122
$0 \times 0 \times 0 \times 0$ intermittent, stored data, upd & msmt aperiod	ic agent init
$0 \times 0.9 0 \times 9.6$ miterimited, stored data, upd & hismit aperiod.	ie, agent init
$0 \times 00 0 \times 02$ attribute value.length = 2	
0x15 0x60	
$0 \times 0 \times 0 \times 55$ attribute-id = MDC_ATTR_ATTRIBUTE_V	ΔΙ ΜΔΡ
$0 \times 00 0 \times 0C$ attribute-value.length = 12	AL_MAI
$0 \times 00 0 \times 02$ AttrValMap.count = 2	
$0 \times 00 = 0 \times 02$ Attr ValMap.length = 8	
0x00 0x06 0x00 0x04 MDC_ATTR_NU_VAL_OBS_SIMP, 4	
0x0A 0x36 0x00 0x04 MDC_ATTR_NO_VAL_OBS_SIMF, 4 0x0A 0x82 0x00 0x08 MDC_ATTR_TIME_STAMP_BO, 8	
	т
$0 \times 00 0 \times 06$ obj-class = MDC_MOC_VMO_METRIC_NU obj-handle = 3 (\rightarrow 3 rd measurement is basa	
·	i delivered)
$0 \times 00 0 \times 04$ attributes.count = 4 $0 \times 00 0 \times 24$ attributes.length = 36	
$0 \times 09 0 \times 2F$ attribute-id = MDC_ATTR_ID_TYPE $0 \times 00 0 \times 04$ attribute-value.length = 4	
0x00 0x04 attribute-value.length = 4 0x00 0x80 0x73 0xF0 MDC_PART_PHD_DM MDC_INS_BASAl	т
	_SMALL
$0 \times 00 0 \times 02$ attribute-value.length = 2	ia acamtimit
0xF0 0x40 intermittent, stored data, upd & msmt aperiodi	ic, agent init
0x09 0x96 attribute-id = MDC_ATTR_UNIT_CODE	
$0 \times 00 0 \times 02$ attribute-value.length = 2	
0x15 0x60 MDC_DIM_INTL_UNIT	AT MAD
$0 \times 0 = 0 \times 55$ attribute-id = MDC_ATTR_ATTRIBUTE_V	AL_MAP
$0 \times 00 0 \times 0C$ attribute-value.length = 12	
$0 \times 00 0 \times 02$ AttrValMap.count = 2	
$0 \times 00 0 \times 08$ AttrValMap.length = 8	
0x0A 0x56 0x00 0x04 MDC_ATTR_NU_VAL_OBS_SIMP, 4	
0x0A 0x82 0x00 0x08 MDC_ATTR_TIME_STAMP_BO, 8	TT T
$0 \times 00 0 \times 05$ obj-class = MDC_MOC_VMO_METRIC_EN	
$0 \times 00 0 \times 04$ obj-handle = 4 ($\rightarrow 4^{th}$ object is operational	status)
$0 \times 00 0 \times 03$ attributes.count = 3	
$0 \times 00 0 \times 1E$ attributes.length = 30	
$0 \times 09 0 \times 2F$ attribute-id = MDC_ATTR_ID_TYPE	
$0 \times 00 0 \times 04$ attribute-value.length = 4	

 $\begin{array}{c} 117 \\ \text{Copyright @ 2015 IEEE. All rights reserved.} \end{array}$

0x00	0x80	0×74	0x30		MDC_PART_PHD_DM MDC_DM_INS_PUMP_OP_STAT
0x0A	0x46				attribute-id = MDC_ATTR_METRIC_SPEC_SMALL
0x00	0x02				attribute-value.length = 2
0xF0	0x40				intermittent, stored data, upd & msmt aperiodic, agent init
0x0A					attribute-id = MDC_ATTR_ATTRIBUTE_VAL_MAP
0x00					attribute-value.length = 12
0x00					AttrValMap.count = 2
0x00					AttrValMap.length = 8
	0x49	0~00	0202		MDC_ATTR_ENUM_OBS_VAL_SIMP_OID, 2
	0x43				MDC_ATTR_ENOM_OBS_VAL_SIMI_OBS, 2 MDC_ATTR_TIME_STAMP_BO, 8
		0.00	0.00		
0x00					obj-class = MDC_MOC_VMO_METRIC_ENUM
0x00					obj-handle = 5 (\rightarrow 5 th object is device status)
0x00					attributes.count = 3
	0x1E				attributes.length = 30
	0x2F				attribute-id = MDC_ATTR_ID_TYPE
0x00					attribute-value.length = 4
0×00	0x80	0×74	0x40		MDC_PART_PHD_DM MDC_DM_INS_PUMP_DEV_STAT
0x0A	0x46				attribute-id = MDC_ATTR_METRIC_SPEC_SMALL
0x00	0x02				attribute-value.length = 2
0xF0	0×40				intermittent, stored data, upd & msmt aperiodic, agent init
$0 \times 0 A$	0x55				attribute-id = MDC_ATTR_ATTRIBUTE_VAL_MAP
0x00	0x0C				attribute-value.length = 12
0x00	0x02				AttrValMap.count = 2
0x00					AttrValMap.length = 8
	0x49	0×00	0×02		MDC_ATTR_ENUM_OBS_VAL_SIMP_OID, 2
	0x82				MDC_ATTR_TIME_STAMP_BO, 8
0x00		01100	01100		obj-class = MDC_MOC_VMO_PMSTORE
0x00					obj-handle = 6 $(6^{th} \rightarrow PM\text{-store Observations})$
0x00					attributes.count = 8
0x00					attributes.length = 66
0x0A					attributes.iengur = 00 attribute-id = MDC_ATTR_PM_STORE_CAPAB
					attribute-value.length = 2
0x00					
0x09	UXUU				pm-store-capab = PMSC_EPI_SEG_ENTRIES
0 00	0 40				PMSC_CLEAR_SEGM_BY_LIST_SUP
0x09					attribute-id = MDC_ATTR_METRIC_STORE_SAMPLE_ALG
0x00					attribute-value.length = 2
0 x 0 4					store-sample-algorithm = ST_ALG_NO_DOWNSAMPLING
	0x41				attribute-id = MDC_ATTR_METRIC_STORE_CAPAC_CNT
0x00					attribute-value.length = 4
	$0 \times 0 0$	0x03	0xE8		store-capacity-count = 1000
0x09					attribute-id = MDC_ATTR_METRIC_STORE_USAGE_CNT
$0 \times 0 0$	0×4				attribute-value.length = 4
0×00	$0 \times 0 0$	0x00	0x35		storage-usage-count = 53
0x09	0x53				attribute-id = MDC_ATTR_OP_STAT
0x00	0x02				attribute-value.length = 2
0x00	0x01				operational-state = enabled
0x0A	0x57				attribute-id = MDC_ATTR_PM_STORE_LABEL_STRING
0x00	0×0E				attribute-value.length = 14
	0x0C				OCTET STRING.length = 12
		0x53	0x74	0x6E	0x72 0x65 0x4F 0x62 0x73 0x65 0x72
					OCTET STRING.data = "PMStoreObser"
0x09	0x51				attribute-id = MDC_ATTR_NUM_SEG

Health informatics—Personal health device communication—Part 10419: Device Specialization—Insulin Pump

0x0	0x02			attribute-value.length = 2
0x0	0x05			number-of-segments $= 5$
0x0A	0x63			attribute-id = MDC_ATTR_CLEAR_TIMEOUT
0x0	0x04			attribute-value.length = 4
0x00	0x00	0x5D	0xC0	clear-timeout = 24000 (3 seconds)

I.3.2.3 Remote operation response event report configuration

The manager responds that it can utilize the agent's configuration. The manager does by sending the confirmed event report response with a config-result of accepted-config.

0xE7	0x00			APDU CHOICE Type (PrstApdu)
0x00	0x16			CHOICE.length = 22
0x00	0×14			OCTET STRING.length = 20
0x12	0x34			invoke-id = 1234 (mirrored from invocation)
0x02	0x01			CHOICE (Remote Operation Response Confirmed Event Report)
0x00	0×0E			CHOICE.length = 14
0x00	0x00			obj-handle = 0 (MDS object)
0x00	0x00	0x00	0x00	currentTime = 0
0x0D	0x1C			event-type = MDC_NOTI_CONFIG
0x00	0×04			event-reply-info.length = 4
0x40	0x00			ConfigReportRsp.config-report-id = 0x4000
0x00	0x0			ConfigReportRsp.config-result = accepted-config

I.3.3 Known configuration

I.3.3.1 General

This exchange takes place when the manager returns the AssociateResult code of accepted because the manager had previously received and processed the configuration corresponding to the dev-config-id sent by the agent. In this case, there is no exchange of configuration information, and the manager and agent have moved into the Operating state.

I.3.3.2 Remote operation invoke event report configuration

Since the manager was already aware of the agent's configuration, the Configuring state is skipped, and no event report invocation is generated by the agent.

I.3.3.3 Remote operation response event report configuration

The Configuring state has been skipped. No event report invocation is generated by the agent, so the manager does not generate any response.

Health informatics—Personal health device communication—Part 10419: Device Specialization—Insulin Pump

I.3.4 Standard configuration

I.3.4.1 General

This exchange takes place when the manager returns the AssociateResult code of accepted because the manager had previously been programmed with the documented standard configuration corresponding to the dev-config-id sent by the agent. In this case, there is no exchange of configuration information, and the manager and agent have moved into the Operating state.

I.3.4.2 Remote operation invoke event report configuration

Since the manager had been programmed with the agent's configuration, the Configuring state is skipped, and no event report invocation is generated by the agent.

I.3.4.3 Remote operation response event report configuration

The Configuring state has been skipped. No event report invocation is generated by the agent, so the manager does not generate any response.

I.4 GET MDS attributes service

I.4.1.1 General

The GET MDS attributes is invoked at any time, when an agent is in Associated state.

I.4.1.2 Get all medical device system attributes request

The manager queries the agent for its MDS Object attributes.

$0 \times E7$	0x00	APDU CHOICE Type (PrstApdu)
0x00	0×0E	CHOICE.length = 14
0x00	0x0C	OCTET STRING.length = 12
0x34	0x56	invoke-id = $0x3456$ (differentiates this message from any other
		outstanding, choice is implementation specific)
0x01	0x03	CHOICE (Remote Operation Invoke Get)
0x00	0x06	CHOICE.length = 6
0x00	0x00	handle = 0 (MDS object)
0x00	0x00	attribute-id-list.count = 0 (all attributes)
0x0	0 x 0 0	attribute-id-list.length = 0

I.4.1.3 Get response with all MDS attributes

The insulin pump agent responds to the manager with its attributes. Furthermore, some optional fields are communicated as well.

0xE7 0x00	APDU CHOICE Type (PrstApdu)
0x00 0x69	CHOICE.length = 105

 $\frac{120}{\text{Copyright @ 2015 IEEE. All rights reserved.}}$

0x67			OCTET STRING.length = 103
0x56			invoke-id = $0x3456$ (mirrored from request)
0x03			CHOICE (Remote Operation Response Get)
0x64			CHOICE.length = 100
0x00			handle = 0 (MDS object)
0x06			attribute-list.count = 6
0x5E			attribute-list.length = 108
0x5A			attribute id = MDC_ATTR_SYS_TYPE_SPEC_LIST
0x08			attribute-value.length = 8
0x01			system-type-spec-list.count = 1
0x04			system-type-spec-list.length = 4
0x13			type = MDC_DEV_SPEC_PROFILE_INSULIN_PUMP
0x01			version = version 1 of the specialization
0x28			attribute id = MDC_ATTR_ID_MODEL
0x14			attribute-value.length = 20
	0x54	0x68	string length = 10 "TheCompany"
0x43	0×6F	0x6D	
0x61	0×6E	0x79	
0x06	0×54	0x68	string length = 6 "TheIPX"
0x49	0x50	0x58	
0x84			attribute-id = MDC_ATTR_SYS_ID
$0 \times 0 A$			attribute-value.length = 10
0x08	0x31	0x32	octet string length = 8 EUI-64
0x34	0x35	0x36	0x37 0x38
0×44			attribute-id = MDC_ATTR_DEV_CONFIG_ID
0x02			attribute-value.length = 2
$0 \times 0 0$			dev-config-id = $0x4000$ (extended-config-start)
0x2D			attribute id = MDC_ATTR_ID_PROD_SPECN
0x12			attribute-value.length = 18
0x01			ProductionSpec.count = 1
0x0E			ProductionSpec.length = 14
0x01			ProdSpecEntry.spec-type = 1 (serial-number)
0x00			ProdSpecEntry.component-id = 0
0x08	0×44	0x45	string length = 8 prodSpecEntry.prod-spec = "DE124567"
0x32	0x34	0x35	
0x37			
0x81			attribute id = MDC_ATTR_TIME_BO
0x08			attribute-value.length = 8
			Base-Offset-Time-Stamp = $2013-12-03T15:14:00.00$
0x00	$0 \times 0 0$	$0 \times 0 0$	
	0x43 0x61 0x06 0x49 0x84 0x0A 0x08 0x34 0x02 0x00 0x2D 0x12 0x01 0x0E 0x01 0x08 0x32 0x37 0x81 0x08 0x67	0x56 0x03 0x64 0x00 0x06 0x5E 0x5A 0x08 0x01 0x04 0x13 0x01 0x28 0x14 0x0A 0x43 0x6F 0x61 0x6E 0x06 0x54 0x49 0x50 0x84 0x0A 0x0A 0x0A 0x0A 0x35 0x44 0x0A 0x0A 0x35 0x44 0x0A 0x0A 0x34 0x35 0x44 0x0A 0x35 0x44 0x0A 0x35 0x44 0x0A 0x35 0x44 0x32 0x01 0x0E 0x01 0x0E 0x01 0x0B 0x34 0x34 0x35	0x56 0x03 0x64 0x00 0x06 0x5E 0x5A 0x08 0x01 0x04 0x13 0x01 0x28 0x14 0x0A 0x6E 0x5A 0x6E 0x6A 0x6E 0x5A 0x01 0x28 0x14 0x0A 0x13 0x01 0x28 0x14 0x0A 0x6F 0x6D 0x61 0x6E 0x79 0x06 0x54 0x68 0x49 0x50 0x58 0x84 0x0A 0x0A 0x0A 0x0A 0x0A 0x0A 0x0A 0x0

I.5 Data reporting

I.5.1 Unconfirmed measurement data transmission

The agent sends a spontaneous event report to the manager with measurement observations.

0xE7 0x00	APDU CHOICE Type (PrstApdu)
0x00 0x92	CHOICE.length = 146
0x00 0x90	OCTET STRING.length = 144
0x12 0x36	invoke-id = 0x1236

ISO/IEEE 11073-10419:2016(E)

IEEE Std 11073-10419-2015
Health informatics—Personal health device communication—Part 10419: Device Specialization—Insulin Pump

0x01	0x01			CHOICE(Remote Operation Invoke Confirmed Event Report)
0x00	0x8A			CHOICE.length = 138
0x00	0x0			obj-handle = 0 (MDS object)
0xFF	0×FF	0×FF	0xFF	event-time = 0xFFFFFFFF
0x0D	0x1D			event-type = MDC_NOTI_SCAN_REPORT_FIXED
0x00	0x80			event-info.length = 128
0xF0	0x00			ScanReportInfoFixed.data-req-id = $0xF000$
0x00	0x00			ScanReportInfoFixed.scan-report-no = 0
0x00	0×07			ScanReportInfoFixed.obs-scan-fixed.count = 7
0x00	0x78			ScanReportInfoFixed.obs-scan-fixed.length = 120
0x00	0×01			ScanReportInfoFixed.obs-scan-fixed.value[0].obj-handle = 1
0x00	0x0C			ScanReportInfoFixed.obs-scan-fixed.value[0].obs-val-data.length = 12
0xFE	0x00	0x00	0x00	Basic-Nu-Observed-Value = 0.0 IU
0xD4	0×67	0 x 4 0	0x38	Base-Offset-Time-Stamp = 2013-12-03T15:14:00.00
0x0	0x00	0x00	0x00	
0x00	0x02			ScanReportInfoFixed.obs-scan-fixed.value[1].obj-handle = 2
0x00	0x0C			ScanReportInfoFixed.obs-scan-fixed.value[1].obs-val-data.length = 12
0xFE	0x00	0x00	0x82	Basic-Nu-Observed-Value = 1.3 IU/hr
0xD4	0x67	0×40	0x38	Base-Offset-Time-Stamp = 2013-12-03T15:14:00.00
0x00	0x00	0x00	0x00	

I.6 Disassociation

I.6.1 Association release request


The insulin pump agent sends the following message to the manager.

0xE4 0x00	APDU CHOICE Type (RlrqApdu)
0x00 0x02	CHOICE.length = 2
0x00 0x00	reason = normal

I.6.2 Association release response

A manager responds to the agent that it can release association.

0xE5 0x00	APDU CHOICE Type (RlreApdu)
0x00 0x02	CHOICE.length = 2
0x00 0x00	reason = normal

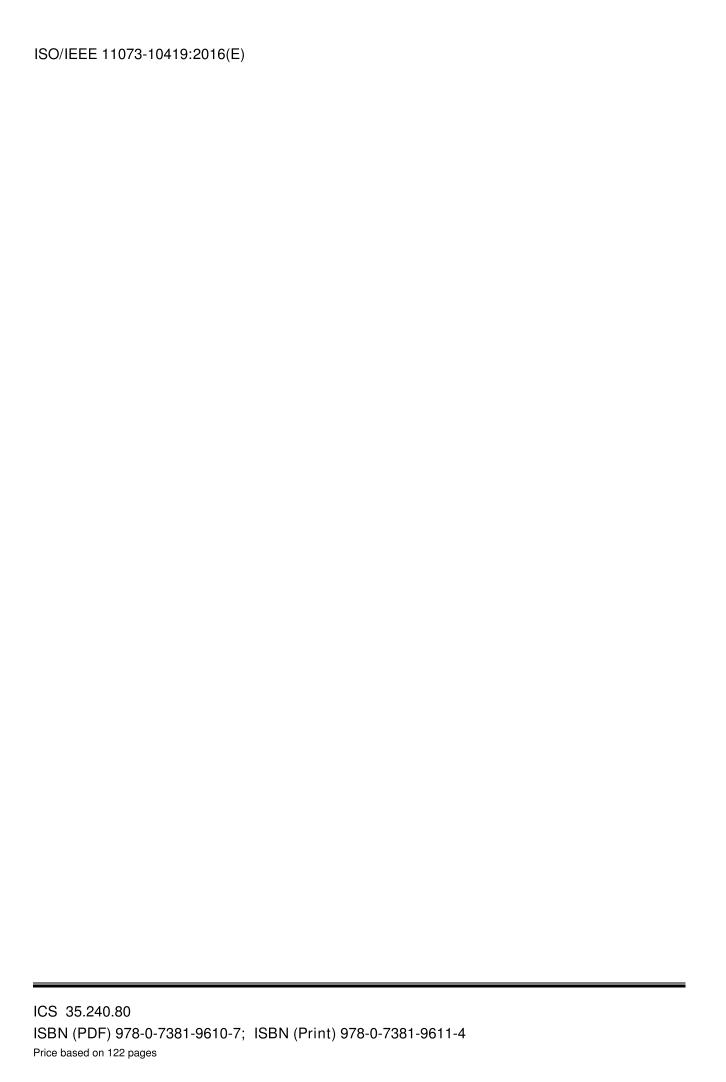

© ISO 2016 – All rights reserved © IEEE 2015 – All rights reserved