# TECHNICAL SPECIFICATION

ISO/TS 17182

First edition 2014-12-15

# Soil quality — Determination of some selected phenols and chlorophenols — Gas chromatographic method with mass spectrometric detection

Qualité du sol — Dosage de quelques phénols et chlorophénols sélectionnés — Méthode par chromatographie en phase gazeuse avec détection par spectrométrie de masse



Reference number ISO/TS 17182:2014(E)

ISO/TS 17182:2014(E)



#### COPYRIGHT PROTECTED DOCUMENT

All rights reserved. Unless otherwise specified, no part of this publication may be reproduced or utilized otherwise in any form or by any means, electronic or mechanical, including photocopying, or posting on the internet or an intranet, without prior written permission. Permission can be requested from either ISO at the address below or ISO's member body in the country of the requester.

ISO copyright office Case postale 56 • CH-1211 Geneva 20 Tel. + 41 22 749 01 11 Fax + 41 22 749 09 47 E-mail copyright@iso.org Web www.iso.org

Published in Switzerland

Cor	ntents	Page
Fore	eword	iv
1	Scope	1
2	Normative references	3
3	Principle	3
4	Interferences 4.1 Interference with sampling and extraction 4.2 Interference with GC-MS determination	3
5	Reagents 5.1 Products used in their commercially available form 5.2 Aqueous solutions 5.3 Standard solutions of phenols	4 4
6	Apparatus and equipment	5
7	Sampling	6
8	Procedure  8.1 Test portion  8.2 Dry matter content  8.3 Blank sample  8.4 Standard sample  8.5 Extraction  8.6 Acetylation  8.7 Gas chromatographic analysis  8.8 Calibration	6 6 6 7 7 7 7
9	<ul> <li>Calculation</li> <li>9.1 Calculation of the content of substance <i>i</i> in the extract</li> <li>9.2 Calculation of the content of selected phenol soil sample per dry matter (mg/kg dry matter)</li> </ul>	8
10	Test report	9
Anne	ex A (informative) Typical concentrations of the standard solutions	10
Anne	ex B (informative) Example of gas chromatographic conditions	11
Bibli	iography	15

#### **Foreword**

ISO (the International Organization for Standardization) is a worldwide federation of national standards bodies (ISO member bodies). The work of preparing International Standards is normally carried out through ISO technical committees. Each member body interested in a subject for which a technical committee has been established has the right to be represented on that committee. International organizations, governmental and non-governmental, in liaison with ISO, also take part in the work. ISO collaborates closely with the International Electrotechnical Commission (IEC) on all matters of electrotechnical standardization.

The procedures used to develop this document and those intended for its further maintenance are described in the ISO/IEC Directives, Part 1. In particular the different approval criteria needed for the different types of ISO documents should be noted. This document was drafted in accordance with the editorial rules of the ISO/IEC Directives, Part 2 (see <a href="https://www.iso.org/directives">www.iso.org/directives</a>).

Attention is drawn to the possibility that some of the elements of this document may be the subject of patent rights. ISO shall not be held responsible for identifying any or all such patent rights. Details of any patent rights identified during the development of the document will be in the Introduction and/or on the ISO list of patent declarations received (see <a href="www.iso.org/patents">www.iso.org/patents</a>).

Any trade name used in this document is information given for the convenience of users and does not constitute an endorsement.

For an explanation on the meaning of ISO specific terms and expressions related to conformity assessment, as well as information about ISO's adherence to the WTO principles in the Technical Barriers to Trade (TBT) see the following URL: Foreword - Supplementary information.

The committee responsible for this document is ISO/TC 190, *Soil quality*, Subcommittee SC 3, *Chemical methods and soil characteristics*.

# Soil quality — Determination of some selected phenols and chlorophenols — Gas chromatographic method with mass spectrometric detection

WARNING — Phenols and chlorophenols are both toxic and corrosive and should be handled with care. Methanol and acetonitrile are toxic and acetic acid is corrosive. Latex or nitrile gloves and eye protection should be worn at all times. Spills should immediately be wiped up with adsorbent tissue and placed in sealed containers used for the disposal of toxic chemicals. Samples should be treated as toxic and harmful. Extreme care shall be taken through all steps of the extraction procedure, which shall be performed in a fume hood. All solvent wastes shall be collected and treated as hazardous waste. Thus, there is a contamination risk in the laboratory.

#### 1 Scope

This International Standard describes the gas chromatographic determination of phenols, methyl phenols, di-methylphenols and tri-methylphenols (see <a href="Table 1">Table 1</a>) and selected chlorophenols (see <a href="Table 2">Table 2</a>) by using mass spectrometric detection in soil samples. This method can also be used for other solid samples, such as sediments and solid wastes. This International Standard describes an acidic liquid extraction of soil, followed by acetylation and then liquid/liquid extraction. Determination takes place by gas chromatography and mass spectrometric detection.

With this method, phenols and chlorophenols can be determined at the lowest of mass concentrations ranging from approximately 0.01 mg/kg to 0.1 mg/kg depending on the component sensitivity and the quantity of sample used. In some cases, complete separation of isomers cannot be achieved. Then the sum is reported.

NOTE With this method, other higher methylated phenols can also be analysed provided that the suitability and the validity of the method is proven.

Table 1 — Target phenolic compounds with relevant mass fractions of respective acetylated compounds for MS-detection

Commound	CAS-RNa	Chemical	Acetylated compounds				
Compound	CAS-RN <sup>a</sup>	formula	Fragmentation <sup>b</sup>				
			1st mass	Relative intensity %	2 <sup>nd</sup> mass	Relative intensity %	
phenol	108-95-2	C <sub>6</sub> H <sub>6</sub> O	94	100	66	26	
2-methylphenol (o-cresol)	95-48-7	C <sub>7</sub> H <sub>8</sub> O	108	100	107	68	
3-methylphenol (m-cresol)	108-39-4	C <sub>7</sub> H <sub>8</sub> O	108	100	107	85	
4-methylphenol (p-cresol)	106-44-5	C <sub>7</sub> H <sub>8</sub> O	108	100	107	92	
2,3-dimethylphenol	596-75-0	C <sub>8</sub> H <sub>10</sub> O	122	90	107	100	
2,4-dimethylphenol	105-67-9	C <sub>8</sub> H <sub>10</sub> O	122	100	107	85	
2,5-dimethylphenol	95-87-4	C <sub>8</sub> H <sub>10</sub> O	122	100	107	90	
2,6-dimethylphenol	576-261-1	C <sub>8</sub> H <sub>10</sub> O	122	100	107	92	

a Chemical abstract service registry number.

b Spectral database for organic compounds (SDBS).

**Table 1** (continued)

CAS-RN <sup>a</sup>	Chemical formula	Acetylated compounds  Fragmentation <sup>b</sup>					
		1st mass	Relative intensity %	2 <sup>nd</sup> mass	Relative intensity %		
95-65-8	C <sub>8</sub> H <sub>10</sub> O	122	99	107	100		
108-68-9	C <sub>8</sub> H <sub>10</sub> O	122	100	107	68		
697-82-5	C <sub>9</sub> H <sub>12</sub> O	136	87	121	100		
2416-94-6	C <sub>9</sub> H <sub>12</sub> O	136	90	121	100		
527-60-6	C <sub>9</sub> H <sub>12</sub> O	136	10	121	85		
527-54-8	C <sub>9</sub> H <sub>12</sub> O	136	70	121	100		
90-00-6	C <sub>8</sub> H <sub>10</sub> O	107	100	122	48		
620-17-7	C <sub>8</sub> H <sub>10</sub> O	122	48	107	100		
123-07-9	C <sub>8</sub> H <sub>10</sub> O	107	100	122	35		
645-56-7	C <sub>9</sub> H <sub>12</sub> O	107	100	136	25		
99-89-8	C <sub>9</sub> H <sub>12</sub> O	121	100	136	28		
90-43-7	C <sub>12</sub> H <sub>10</sub> O	170	100	141	23		
Internal standards							
367-12-4	C <sub>6</sub> H <sub>5</sub> FO	112	100	64	58		
371-41-5	C <sub>6</sub> H <sub>5</sub> FO	112	100	64	18		
452-72-2	C <sub>7</sub> H <sub>7</sub> FO	126	100	97	29		
367-32-8	C <sub>6</sub> H <sub>5</sub> FO <sub>2</sub>	128	100	75	28		
	95-65-8 108-68-9 697-82-5 2416-94-6 527-60-6 527-54-8 90-00-6 620-17-7 123-07-9 645-56-7 99-89-8 90-43-7 367-12-4 371-41-5 452-72-2 367-32-8	95-65-8	$\begin{array}{c ccccccccccccccccccccccccccccccccccc$	$\begin{array}{c ccccccccccccccccccccccccccccccccccc$	$\begin{array}{c ccccccccccccccccccccccccccccccccccc$		

Chemical abstract service registry number.

 $Table\ 2-Target\ chlorophenolic\ compounds\ with\ relevant\ mass\ fractions\ of\ respective$ acetylated compounds for MS-detection

Compound	CAS-RNa	Chemical formula	Acetylated compounds Fragmentation			
			1 <sup>st</sup> mass	Relative intensity %	2 <sup>nd</sup> mass	Relative intensity %
2-chlorophenol	95-57-8	C <sub>6</sub> H <sub>5</sub> OCl	128	100	130	34
3-chlorophenol	108-43-0	C <sub>6</sub> H <sub>5</sub> OCl	128	100	130	29
4-chlorophenol	106-48-9	C <sub>6</sub> H <sub>5</sub> OCl	128	100	130	35
2,6-dichlorophenol	87-65-0	C <sub>6</sub> H <sub>4</sub> OCl <sub>2</sub>	162	100	164	67
2,4-dichlorophenol	120-83-2	C <sub>6</sub> H <sub>4</sub> OCl <sub>2</sub>	162	100	164	66
2,5-dichlorophenol	583-78-8	C <sub>6</sub> H <sub>4</sub> OCl <sub>2</sub>	162		164	
3,5-dichlorophenol	591-35-5	C <sub>6</sub> H <sub>4</sub> OCl <sub>2</sub>	162	100	164	66
2,3-dichlorophenol	576-24-9	C <sub>6</sub> H <sub>4</sub> OCl <sub>2</sub>	162	100	164	66
3,4-dichlorophenol	95-77-2	C <sub>6</sub> H <sub>4</sub> OCl <sub>2</sub>	162	100	164	66
2,4,6-trichlorophenol	88-06-2	C <sub>6</sub> H <sub>3</sub> OCl <sub>3</sub>	196	100	198	98

Spectral database for organic compounds (SDBS).

#### 2 Normative references

The following documents, in whole or in part, are normatively referenced in this document and are indispensable for its application. For dated references, only the edition cited applies. For undated references, the latest edition of the referenced document (including any amendments) applies.

ISO 3696, Water for analytical laboratory use — Specification and test methods

ISO 10381-1, Soil quality — Sampling — Part 1: Guidance on the design of sampling programmes

 ${\sf ISO\,11465}$ , Soil quality — Determination of dry matter and water content on a mass basis — Gravimetric method

ISO 14507, Soil quality — Pretreatment of samples for determination of organic contaminants

ISO 18512, Soil quality — Guidance on long and short term storage of soil samples

ISO 22892, Soil quality — Guidelines for the identification of target compounds by gas chromatography and mass spectrometry

#### 3 Principle

The method described here is based on two steps. The first step includes a solid/liquid extraction: the phenols are extracted from field moist soil with methanol at low pH. The second step is the derivatization of phenols obtained in an aliquot of methanol in aqueous carbonate solution with acetic anhydride; the derivatives formed are extracted from this sample with cyclohexane. The cyclohexane fraction is analyzed by gas chromatography with mass spectrometric detection.

#### 4 Interferences

#### 4.1 Interference with sampling and extraction

Standard laboratory glassware, appropriately cleaned and free of interfering compounds. Do not use any kind of plastic containers, since phenols can be adsorbed to these. Plastic materials can also contribute their impurities to the sample material. Surfactants, present in the sample, might disturb the proper extraction of phenols.

#### Interference with GC-MS determination 4.2

Substances which have the same or nearly the same retention time of the target analytes and which give the same mass fragments (i.e. isomeric phenols other than stated in the table) might interfere with the determination.

#### 5 Reagents

During the analysis, unless otherwise stated, use only reagents of recognized analytical grade and water of at least grade 1 as defined in ISO 3696.

Only commercially available certified standards of high purity shall be used.

#### 5.1 Products used in their commercially available form

- Gas chromatography gases, of a purity as recommended by the gas chromatography manufacturer.
- 5.1.2 **Methanol**, CH<sub>3</sub>OH, with a mass fraction of 99,5 %.
- 5.1.3 **Cyclohexane**, C<sub>6</sub>H<sub>12</sub> or other hydrocarbon solvent like heptane.
- **Hydrochloric acid**, HCl, concentrated, with a mass fraction of 37 %. 5.1.4
- 5.1.5 Sodium hydroxide, NaOH.
- 5.1.6 **Potassium carbonate**, K<sub>2</sub>CO<sub>3</sub>.
- **Acetic anhydride**, (CH<sub>3</sub>CO)<sub>2</sub>O, anhydrous. 5.1.7
- 5.1.8 **Sodium sulfate**, Na<sub>2</sub>SO<sub>4</sub>, anhydrous.

Weigh portions of 2 g of sodium sulfate into test tubes with polytetrafluoroethylene (PTFE) lined caps. Dry the test tubes at  $(550 \pm 20)$  °C for 4 h ± 30 min without caps. Take them to a desiccator and let them cool. When cooled, put on the caps and store the tubes at room temperature. Sodium sulfate can also be dried in bigger portions and stored in a desiccator after cooling. Weigh portions of 2 g into tubes when needed.

- **5.1.9** Selected phenol and chlorophenol standards, listed in Table 1 and Table 2.
- **5.1.10 Selected fluorophenol standards**, listed in Table 1 as internal standard.

Other compounds such as isotopic and deuterated phenols which might not interfere with the GC-MS measurement can also be used as internal standards.

- 5.2 **Aqueous solutions**
- 5.2.1 **Sodium hydroxide solution**, NaOH, c(NaOH) = 0.1 mol/l.
- 5.2.2 **Sodium hydroxide solution**, NaOH, c(NaOH) = 0.5 mol/l.
- 5.2.3 **Potassium carbonate solution**,  $K_2CO_3$ ,  $c(K_2CO_3) = 0.1 \text{ mol/l}$ .
- 5.2.4 **Potassium carbonate solution**,  $K_2CO_3$ ,  $c(K_2CO_3) = 5.2 \text{ mol/l}$ .

#### 5.3 Standard solutions of phenols

#### 5.3.1 Internal standard solution

#### **5.3.1.1** General

Fluorophenols can be used as suitable internal standard.

#### 5.3.1.2 Stock solution

Prepare the stock solution by weighing fluorophenols (5.1.10) and dissolving it in methanol (5.1.2). A suitable concentration is 1 mg/ml. Divide the stock solution into 5 ml bottles with tight caps, 1,5 ml to each, and store at -18 °C.

NOTE Stock solutions are stable for at least half a year when stored in the dark at 4  $^{\circ}$ C. They are stable at least one year at –18  $^{\circ}$ C.

#### 5.3.1.3 Working solution

Prepare this working solution by diluting the stock solution (5.3.1.2) with distilled water. A suitable concentration is  $100 \,\mu\text{g/ml}$ .

#### 5.3.2 Individual standard phenol solutions

#### **5.3.2.1** General

Typical concentrations of standard solutions are given in Annex A.

#### 5.3.2.2 Stock solution

Weigh out each of the phenol standards (5.1.9) into the same or into separate measuring flasks and dissolve them in methanol (5.1.2). A suitable concentration is 0,5 mg/ml. Divide the stock solution(s) into 5 ml bottles with tight PTFE lined caps, 1,5 ml to each, and store at –18 °C.

NOTE Stock solutions are stable for at least half a year when stored in the dark at 4  $^{\circ}$ C. They are stable at least one year at -18  $^{\circ}$ C.

#### 5.3.2.3 Working solution

Dilute stock solution(s) (5.3.2.2) with distilled water to one working solution. A suitable concentration is  $50 \mu g/ml$ .

#### 6 Apparatus and equipment

#### 6.1 General requirements

Standard laboratory glassware appropriately cleaned and free of interfering compounds. Do not use any kind of plastic containers, since the phenols can adsorb to these or plastic materials can contribute their impurities to the sample material.

#### 6.2 Gas chromatography system

**6.2.1 Gas chromatograph**, equipped with a non-discriminating injection system, capillary column, and a mass spectrometric detector (GC-MS).

**6.2.2 Capillary columns**, each comprising a 5 % phenyl-methyl silicone stationary phase coated onto fused silica capillary column or an equivalent chemically bonded phase column. Their dimensions should be sufficient to separate the critical pairs mentioned below. In general, column length should be at least 30 m. Internal diameter 0,25 mm and film thickness 0,2  $\mu$ m.

Sufficient resolution (0,7) between the chromatographic peaks of critical pairs as m-cresol and p-cresol shall be set as quality criteria for the capillary column. Otherwise, the sum of both isomers is to be reported.

#### 6.3 Ultrasonic bath

**6.4 Shaking device**, with horizontal movement (200 strokes per min to 300 strokes per min).

#### 7 Sampling

Take samples according to ISO 10381-1 and pre-treat the samples according to ISO 14507.

Samples shall be taken in glass vessels with PTFE caps. It is recommended that the vessel be filled completely.

Store the samples in the dark in the laboratory according to ISO 18512. Phenols may be subject to microbial conversion under certain conditions. It is recommended that samples are frozen if they are stored for more than 2 days.

As phenols in soil are subject to biodegradation, some of the phenols are volatile; the sample containers should be filled on site with methanol to preserve the phenols. Bottles containing reagent should be preweighed in the laboratory before filled with sample material and re-weighed on receipt in the laboratory.

#### 8 Procedure

#### 8.1 Test portion

At least 10 g with an accuracy of 0,01 g of field moist soil is taken directly from the vessel. The sample is taken with a spoon or a spatula by mixing the soil as well as possible.

#### 8.2 Dry matter content

The dry matter content of the soil shall be determined in a subsample of the same vessel according to ISO 11465.

#### 8.3 Blank sample

Blank samples shall be treated exactly the same way as normal samples but the test portion (8.1) is replaced by 10 g of sample matrix known to be free of phenols.

#### 8.4 Standard sample

Standard samples shall be treated exactly the same way as normal samples but the test portion (8.1) is replaced by 50  $\mu$ l of standard working solution (5.3.2.3) and 10 g of sample matrix known to be free of phenols.

#### 8.5 Extraction

Take a test portion of 10 g soil and place it in a conical flask. Add 50  $\mu$ l of internal standard working solution (5.3.1.3), 75 ml of methanol (5.1.2), and 0,5 ml of hydrochloric acid (5.1.4).

In case of high concentrations, more of the internal standard should be added.

Sonicate the sample for 10 min and then shake for at least 30 min with 200 strokes per min to 300 strokes per min. Allow the sample to settle and transfer an aliquot of the methanol extract (e.g. 10 ml) to a separating funnel with 50 ml of potassium carbonate aqueous solution (5.2.3).

Some types of soil are acidic, contain carbonates, or have a high buffer capacity. In these cases, the amounts of acid and base added are not enough to reach sufficiently low pH. When phenols are to be analysed in these types of soil, check the pH-values to be pH < 3 in the extraction step.

#### 8.6 Acetylation

#### 8.6.1 Acetylation of extract

Acetylate the phenols and chlorophenols in 50 ml potassium carbonate solution (5.2.3) as follows: Add 2 ml of sodium hydroxide aqueous solution (5.2.2) and 1 ml of acetic anhydride (5.1.7) to the extract and shake the mixture vigorously for 2 min to release any carbon dioxide formed in the funnel. Let the mixture stand for 10 min while shaking occasionally, and then add 10 ml of cyclohexane (5.1.3). Shake the funnel and let the two phases separate. Transfer as large a portion as possible of the cyclohexane phase to a tube containing 2 g of Na<sub>2</sub>SO<sub>4</sub> (5.1.8) for drying. After shaking, transfer the cyclohexane solution to another vial with Na<sub>2</sub>SO<sub>4</sub> and store at 4 °C. Analysis of the acetylated phenols shall be done as soon as possible (within 48 h) since the acetates are rather labile towards hydrolysis.

#### 8.6.2 Acetylation of control samples

Take 50  $\mu$ l of standard working solution and add it to a separation funnel with 50 ml of potassium carbonate solution (5.2.3) and shake the funnel. Add 2 ml of sodium hydroxide aqueous solution (5.2.2) and 1 ml of acetic anhydride (5.1.7) to the separation funnel and shake the mixture vigorously for 2 min to release any carbon dioxide formed in the funnel. Let the mixture stand for 10 min while shaking occasionally, and then add 10 ml of cylohexane (5.1.3). Shake the funnel and let the two phases separate. Transfer as large a portion as possible of the cyclohexane phase to a tube containing 2 g of Na<sub>2</sub>SO<sub>4</sub> (5.1.8) for drying. After shaking, transfer the cyclohexane solution to another vial with Na<sub>2</sub>SO<sub>4</sub> and store at 4 °C. Analysis of control samples shall be done as soon as possible (within 48 h).

#### 8.7 Gas chromatographic analysis

Set up the GC-MS system equipped with appropriate columns (6.2.2). Optimize the gas flows to obtain sufficient separation. Ensure a stable condition. An example of gas chromatographic conditions and temperature programme is given in Annex B. The mass spectrometer is tuned in accordance with the manufacturer's instructions. Chromatograms are recorded in full scan or selected ion mode (SIM). Examples for characteristic masses of analytes are given in Table 1 and Table 2.

Use ISO 22892 for identification of the analytes.

#### 8.8 Calibration

Use extracted and acetylated standard samples (8.4) as calibration solutions. Calibration solutions should be prepared according to 8.4 by adding a constant definite volume of internal standard working solution (5.3.1.3) adding different volumes of working solution (5.3.2.3). For each compound, a separate calibration function shall be established, consisting of at least five measurement points. It is permissible to examine several compounds in one calibration.

The calibration function obtained for a particular compound is valid only for the established concentration range and the respective sample preparation used. It depends on the operational conditions of the entire analytical system.

Calibration shall be based on peak height or area and on the response of the internal standard.

#### ISO/TS 17182:2014(E)

Calculate the calibration function by regression using the ratio  $y_{ie}/y_{se}$  as a function of the ratio of  $\rho_{ie}/\rho_{se}$ according to Formula (1):

$$\frac{y_{ie}}{y_{se}} = a_{is} \cdot \frac{\rho_{ie}}{\rho_{se}} + b_{is} \tag{1}$$

where

- $y_{ie}$  is the (dependent variable) measured response of the compound i during calibration, depending on  $\rho_{ie}$ ;
- $y_{se}$  is the (dependent variable) measured response (signal area or signal height) of the internal standard s during calibrating, depending on  $\rho_{se}$ ;
- $\rho_{ie}$  is the (independent variable) mass concentration of the compound *i* in the calibration solution, in micrograms per litre (µg/l);
- $\rho_{se}$  is the (independent variable) mass concentration of the internal standard s in the calibration solution, in micrograms per litre (µg/l);
- $a_{is}$  is the slope of the calibration curve of the ratio of the measured values  $y_{ie}/y_{se}$  as a function of the corresponding ratio of the mass concentration  $\rho_{ie}/\rho_{se}$  (corresponds to the substance-specific response factor, often referred to as  $f_i$ ), dimensionless;
- $b_{is}$  is the axis intercept of the calibration curve on the ordinate, dimensionless.

#### **Calculation**

#### Calculation of the content of substance i in the extract

Calculate the concentration  $c_{i,j}$  of the substance i in the extract according to Formula (2):

$$c_{i,j} = \frac{(y_{i,j} / y_{s,j}) - b_{is}}{a_{is}} \cdot c_{s,j}$$
 (2)

where

- $c_{i,i}$  is the concentration of selected phenol i in the extract, in micrograms per millilitre (µg ml) =  $\rho_{iW}$ ;
- $y_{i,i}$  is the height or area of the selected phenol i;
- $y_{s,i}$  is the height or area of internal standard;
- $c_{s,i}$  is the amount of internal standard added to the extract, in micrograms per millilitre (µg ml);
- $b_{is}$  is the axis intercept of the calibration curve on the ordinate, dimensionless;
- $a_{is}$  is the slope of the calibration curve of the ratio of the measured values  $y_{ie}/y_{se}$  as a function of the corresponding ratio of the mass concentration  $\rho_{ie}/\rho_{se}$  (corresponds to the substance-specific response factor, often referred to as  $f_i$ ), dimensionless.

# 9.2 Calculation of the content of selected phenol soil sample per dry matter (mg/kg dry matter)

Calculate the content of a specific phenol the soil sample by using Formula (3):

$$w_{idm} = \frac{\rho_{iW} \cdot V_{E}}{m_{dm}} \tag{3}$$

where

 $w_{\text{idm}}$  is the content of the specific individual phenol "i" in the sample, in milligrams per kilogram, mg/kg, of dry matter;

 $\rho_{iW}$  is the concentration of the analyte "i" in the extract, in micrograms per millilitre (µg/ml);

 $m_{\rm dm}$  is the mass of the test sample of dry matter used for extraction, in grams (g);

 $V_{\rm E}$  is the total volume of the extract, i.e (volume of methanol added to the soil sample +volume of water present in the field moist sample obtained from the determination of dry matter content according to ISO 11465) in millilitres (ml).

Report the content of the individual phenols with two significant figures.

#### 10 Test report

At least the following information shall be included in the report:

- a) all information necessary for identification of the sample;
- b) a reference to this International Standard, including the year of publication (i.e. ISO 17182:2014);
- c) the concentration of each phenol in mg/kg dry matter according to Clause 9;
- d) the interval between the sampling and the extraction;
- e) any disturbances observed during determination, for instance other peaks;
- f) details of any operation not specified in this International Standard or regarded as optional, as well as any factor which might have affected the results.

## Annex A

(informative)

# Typical concentrations of the standard solutions

Typical concentrations for an internal standard solution and standard solutions are given in Table A.1 and Table A.2.

#### Table A.1 — Internal standard

Internal standard	Stock solution (5.3.1.2)	Working solution (5.3.1.3)
Fluorophenols	1 mg/ml	100 μg/ml

#### Table A.2 — Individual phenol standards

Phenols listed in <u>Table 1</u>	Stock solution (5.3.2.2)	Working solution (5.3.2.3)
Phenols	400 μg/ml	40 μg/ml

#### **Annex B**

(informative)

### **Example of gas chromatographic conditions**

#### **B.1 Phenols**

#### B.1.1 GC-MS conditions for acetylated phenols, Chromatogram B.1

Gas chromatograph: Agilent 6890<sup>1)</sup>

Mass spectrometer: Agilent 59731)

Autosampler: Agilent 76831)

Injection: split/splitless

Injection volume: 1 μl

Injection mode: splitless

Injection temperature: 250 °C

Capillary column: DB5ms<sup>2)</sup>, 30 m × 0,25 mm ID, film thickness 0,25  $\mu$ m

Carrier gas: Helium, 0,7 ml/min

Temperature program: 50 °C – 1 min

10 °C/min to 130 °C (5 min)

10 °C/min to 220 °C

10 °C/min to 280 °C (4 min)

Transfer-Line-Temperature: 290 °C

Mode: SIM (Single Ion Monitoring)

<sup>1)</sup> Agilent 6890, Agilent 5973 and Agilent 7683 are examples of suitable products available commercially. This information is given for the convenience of users of this International Standard and does not constitute an endorsement by ISO of these products.

<sup>2)</sup> DB5ms is an example of a suitable product available commercially. This information is given for the convenience of users of this International Standard and does not constitute an endorsement by ISO of this product.

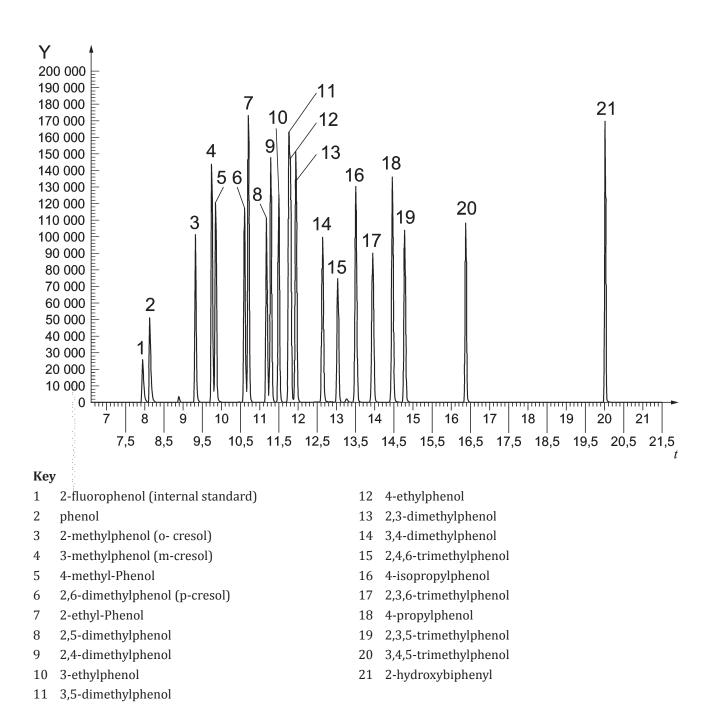


Figure B.1 — GC-MS chromatogram of acetylated alkylphenol standards in SIM-Mode  $(c = 30 \mu g/l per individual phenol)$ 

### **B.2 Chlorophenols**

#### B.2.1 GC-MS conditions for acetylated chlorophenols, Chromatogram B.2

Gas chromatograph: Agilent 68901) Mass spectrometer: Agilent 59731) Autosampler: Agilent 76831) Injection: split/splitless

Injection volume: 1 µl

Injection mode: splitless

Injection temperature: 250 °C

Capillary column: Phenomenex 7KG-G002-11 ZB5<sup>3</sup>)

Carrier gas: Helium, 1,9 ml/min

Temperature program: 60 °C – 0,20 min

30 °C/min to 150 °C (0 min)

8 °C/min to 179 °C (0 min)

20 °C/min to 240 °C (2 min)

50 °C/min to 320 °C (5 min)

Transfer-Line-Temperature: 310 °C

Mode: SIM (Single Ion Monitoring)

<sup>3)</sup> Phenomenex 7KG-G002-11 ZB5 is an example of a suitable product available commercially. This information is given for the convenience of users of this International Standard and does not constitute an endorsement by ISO of this product.

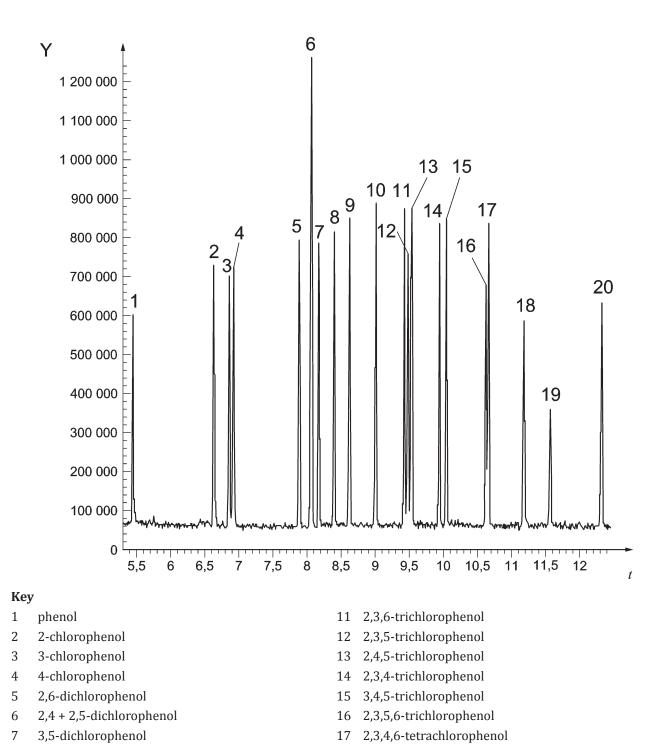


Figure B.2 — GC-MS chromatogram of acetylated chlorophenols in SIM-Mode

18

19

20

2,3-dichlorophenol

3,4-dichlorophenol

2,4,6-trichlorophenol

8

9

2,3,4,5-tetrachlorophenol

TBP (internal standard)

pentachlorophenol

# **Bibliography**

- $[1] \hspace{2em} \hbox{EN 12673, Water quality} \hspace{-0.5em} \hspace{-0.5em} \textit{Gas chromatographic determination of some selected chlorophenols in water}$
- [2] ISO 11709, Soil quality Determination of selected coal-tar-derived phenolic compounds using high performance liquid chromatography (HPLC)
- [3] ISO 14154, Soil quality Determination of some selected chlorophenols Gas-chromatographic method with electron-capture detection



ICS 13.080.10

Price based on 15 pages